CD317 (BST2, PDCA-1) Monoclonal Antibody (eBio927), PE-Cyanine7, eBioscience™

r roudet Details	
Size	100 µg
Species Reactivity	Mouse
Published Species	Mouse
Host/Isotype	Rat / IgG2b, kappa
Recommended Isotype Control	Rat IgG2b kappa Isotype Control (eB149/10H5), PE-Cyanine7, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	eBio927
Conjugate	PE-Cyanine7
Excitation/Emission Max	569/780 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2573440

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	0.25 μg/test	5 Publications

Product Specific Information

Description: The eBio927 monoclonal antibody reacts with PDCA-1 (BST2, CD317), a specific marker of plasmacytoid dendritic cells (pDC), also known as type I IFN-producing cells (IPC) in the naive mouse. mouse IPCs are typically CD11c+, CD11b-, B220+, Ly-6C+, and CD62L+. PDCA-1 is predominantly expressed by IPCs in the naive mouse which represents a very minor population (<0.5%) of splenocytes. It is upregulated on numerous cell types following stimulation which triggers an IFN response. PDCA-1 cycles between cell surface and intracellular compartments and may function to regulate trafficking of secreted cytokines. PDCA-1 (BST2) is the protein recognized by antibody 120G8.

The eBio927 monoclonal antibody has also been shown to have functional activity. The epitope recognized by eBio927 is distinct from eBio129c; thus, the antibodies can be used to costain, purify and identify pDCs.

Applications Reported: This eBio927 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This eBio927 antibody has been tested by flow cytometric analysis of mouse splenocytes. This can be used at less than or equal to 0.25 μ g per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Light sensitivity: This tandem dye is sensitive photo-induced oxidation. Please protect this vial and stained samples from light.


1

Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222) (100 μ L cell sample + 100 μ L IC Fixation Buffer) or 1step Fix/Lyse Solution (Product # 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 488-561 nm; Emission: 775 nm; Laser: Blue Laser, Green Laser, Yellow-Green Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For CD317 (BST2, PDCA-1) Monoclonal Antibody (eBio927), PE-Cyanine7, eBioscience™

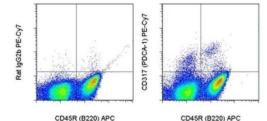

CD317 (BST2, PDCA-1) Antibody (25-3172-82) in Flow

Staining of SJL splenocytes with Anti-Human/Mouse CD45R (B220) APC (Product # 17-0452-82) and 0.125 µg of Rat IgG2b K Isotype Control PE-Cyanine7 (Product # 25-4031-82) (left) or 0.125 µg of Anti-Mouse CD317 (BST2, PDCA-1) PE-Cyanine7 (right). Cells in the lymphocyte gate were used for analysis.

5 References

Flow Cytometry (5)

Journal of immunology research pDC Activation by TLR7/8 Ligand CL097 Compared to TLR7 Ligand IMQ or TLR9 Ligand CpG. "25-3172-82 was used in Flow cytometry/Cell sorting to aid in the development of a new strategy for future tumor immunotherapy." Authors: Wu J,Li S,Li T,Lv X,Zhang M,Zang G,Qi C,Liu YJ,Xu L,Chen J	Year 2019 Species Mouse	
Stem cell research & therapy Circulating healing (CH) cells expressing BST2 are functionally activated by the injury-regulated systemic factor HGFA.	Year 2018 Species	
"25-3172 was used in Flow cytometry/Cell sorting to show that bone marrow stromal cell antigen 2 allows the isolation of a population of circulating progenitors, the circulating healing cells, characterized by a distinctive core signature." Authors: Lo Sicco C,Reverberi D,Villa F,Pfeffer U,Quarto R,Cancedda R,Tasso R	Mouse	

View more Flow references on thermofisher.com

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuate. Unless otherwise stated in the Documentation, this warranty is provided herein is valid only when used by properly trained individuates of the general type and quality of goods and does not represent that any Product will conform to such model or sample furnities. EXPRESS 0R IMPED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICURS DR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANT PERIOD IS LIMITED TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OPTION. THERE IS NO BILICATION TO REPARK, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OPTION. THERE IS NO BILICATION TO REPARK, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLICENCE OF OR BY BUYER', (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. UNless otherwise asystead on the Product in the documentation or animals.

2