TER-119 Monoclonal Antibody (TER-119), Brilliant Ultra Violet[™] 805, eBioscience[™]

Product Details	
Size	100 µg
Species Reactivity	Mouse
Host/Isotype	Rat / IgG2b, kappa
Recommended Isotype Control	Rat IgG2b kappa Isotype Control (eB149/10H5), Brilliant Ultra Violet™ 805, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	TER-119
Conjugate	Brilliant Ultra Violet™ 805
Excitation/Emission Max	349/804 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2896157

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	0.25 µg/test	-

Product Specific Information

Description: The TER-119 monoclonal antibody reacts with mouse erythroid cells from early proerythroblast to mature erythrocyte stages. The TER-119 antigen is present in yolk sac, fetal and newborn liver, but is not expressed by cells carrying BFU-E and CFU-E activities. Several erythroleukemia cell lines tested so far are negative for expression of TER-119 antigen even after dimethylsulfoxide stimulation. Biochemical and molecular analysis of the TER-119 antigen indicate that this molecule is associated with the surface glycophorin A, but is not a typical glycophorin.

Applications Reported: This TER-119 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This TER-119 antibody has been tested by flow cytometric analysis of mouse bone marrow cells. This may be used at less than or equal to 0.25 μ g per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Brilliant Ultra Violet[™] 805 (BUV805) is a tandem dye that emits at 797 nm and is intended for use on cytometers equipped with an ultraviolet (355 nm) laser. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright, Brilliant Violet[™], Brilliant Ultra Violet[™], or other polymer dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) or Brilliant Stain Buffer (Product # 00-4409-75) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer or Brilliant Stain Buffer for more information.

Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.


Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222) (100 μ L of cell sample + 100 μ L of IC Fixation Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Our internal testing suggests that Brilliant Ultra Violet[™] 805 (BUV805) is compatible with short-term methanol-based fixation, but should not be stored in buffers containing methanol for longer than one hour.

Excitation: 355 nm; Emission: 797 nm; Laser: Ultraviolet Laser.

BRILLIANT ULTRA VIOLET[™] is a trademark or registered trademark of Becton, Dickinson and Company or its affiliates, and is used under license. Powered by Sirigen[™].

Product Images For TER-119 Monoclonal Antibody (TER-119), Brilliant Ultra Violet™ 805, eBioscience™

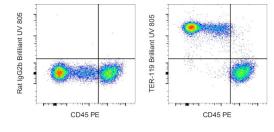

TER-119 Antibody (368-5921-82) in Flow

C57BL/6 mouse bone marrow cells were stained with CD45 Monoclonal Antibody, PE (Product # 12-0451-82) and 0.125 µg of Rat IgG2b kappa Isotype Control, Brilliant Ultra Violet 805 (BUV805) (Product # 368-4031-81) (left) or 0.125 µg of TER-119 Monoclonal Antibody, Brilliant Ultra Violet 805 (BUV805) (right). Total viable cells were used for analysis, as determined by 7-AAD (Product # 00-6993-50).

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of saids expression of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is performed to expression provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is minited to one year from date of shipment when the Product is subjected to normal, proper and linended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTES, EXPRESS OR INNPELD, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES OF MERCHANTABULTY, ITTHESS FOR ANY PARTICURCTS DURING THE WARRANT PERION IS LIMITED TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS DURING THE WARRANT PERION IS LIMITED TO REPAR, REPLACE COR REFUND FOR THE RON-CONFORMING PRODUCTS SOLE OPTION TO REPARIN, REPLACE OR REFUND FOR THE RON-CONFORMING PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (III) USE OF THE PRODUCTS IN A MANKER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS IN A MANKER FOR WHICH THEY WERE NOT IN INSUES, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANKER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product in the documentation accompanying the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro