Performance guarenteed

CD2 Monoclonal Antibody (RPA-2.10), eFluor™ 450, eBioscience™

Product	Details

Size	100 Tests
Species Reactivity	Baboon, Chimpanzee, Cynomolgus monkey, Human, Non-human primate, Pig, Rhesus monkey
Published Species	Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), eFluor™ 450, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	RPA-2.10
Conjugate	eFluor™ 450
Excitation/Emission Max	405/445 nm
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.2% BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2574006

Applications	Tested Dilution	Publications
Immunohistochemistry (Frozen) (IHC (F))	-	1 Publication
Flow Cytometry (Flow)	5 μL (0.06 μg)/test	5 Publications

Product Specific Information

Description: The RPA-2.10 monoclonal antibody reacts with human CD2, a 50 kDa cell surface receptor expressed by a majority of thymocytes, all mature T cells and subset of NK cells. CD2 is a ligand for CD58 in the human and is involved in adhesion and activation of T cells. RPA-2.10 blocks mixed lymphocyte reaction.

RPA-2.10 crossreacts to non-human primates and pigs.


Applications Reported: This RPA-2.10 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This RPA-2.10 antibody has been pre-titrated and tested by flow cytometric analysis of normal human peripheral blood cells. This can be used at 5 μ L (0.06 μ g) per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.

eFluor® 450 is an alternative to Pacific Blue®. eFluor® 450 emits at 445 nm and is excited with the Violet laser (405 nm). Please make sure that your instrument is capable of detecting this fluorochome.

Excitation: 405 nm; Emission: 445 nm; Laser: Violet Laser.

Product Images For CD2 Monoclonal Antibody (RPA-2.10), eFluor™ 450, eBioscience™

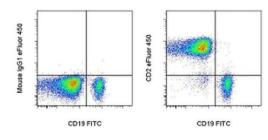

CD2 Antibody (48-0029-42) in Flow

Staining of normal human peripheral blood cells with Anti-Human CD19 FITC (Product # 11-0199-42) and Mouse IgG1 K Isotype Control eFluor® 450 (Product # 48-4714-82) (left) or Anti-Human/Non-Human Primate CD2 eFluor® 450 (right). Cells in the lymphocyte gate were used for analysis.

□ 6 References

Immunohistochemistry (Frozen) (1)

The journal of histochemistry and cytochemistry : official journal of the Histochemistry Society	Year 2009
A monoclonal antibody selection for immunohistochemical examination	
of lymphoid tissues from non-human primates.	
Authors: Kap YS,van Meurs M,van Driel N,Koopman G,Melief MJ,Brok HP,Laman JD,'t Hart BA	
Flow Cytometry (5)	
Biophysical journal	Year
Constraints on GPCR Heterodimerization Revealed by the Type-4	2019
Induced-Association BRET Assay.	
"Published figure using CD2 monoclonal antibody (Product # 48-0029-42) in Flow Cytometry"	
Authors: Felce JH,MacRae A,Davis SJ	
European journal of immunology	
Tissue factor expression by myeloid cells contributes to protective	2016
immune response against Mycobacterium tuberculosis infection.	Species
"48-0029 was used in Flow cytometry/Cell sorting to determine the role of tissue factor expressed by myeloid cells during Mycobacterium tuberculosis infection."	Human
Authors: Venkatasubramanian S,Tripathi D,Tucker T,Paidipally P,Cheekatla S,Welch E,Raghunath A,Jeffers A,	

Authors: Venkatasubramanian S, Tripathi D, Tucker T, Paidipally P, Cheekatla S, Welch E, Raghunath A, Jeffers A, Tvinnereim AR, Schechter ME, Andrade BB, Mackman N, Idell S, Vankayalapati R

View more Flow references on thermofisher.com

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is initiated to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample turnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTES, EXPRESS OR INMELD, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES OF MERCHANTABILITY, ITNESS FOR ANY PARTICURS OR NON INFRINGEMENT. BUYER'S SCLUSVE REMEUP FOR NON-CONFORMING PRODUCTS DURING THE WARRANT REPRIOD IS LIMITED TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS ON DURING THE WARRANTE, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OPTION. THERE IS NO BULGATION TO REPAR, REPLACE OR REFUND FOR THE RODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (III) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANKER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Na MANKER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. IN a MANKER FOR WHICH THEY were not press, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeuricus uses, or any type of consumption by or application to human or animals.

2