



## CD73 Monoclonal Antibody (AD2), eFluor™ 450, eBioscience™

100 Tests
Human
Human
Mouse / IgG1, kappa
Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), eFluor™ 450, eBioscience™
Monoclonal
Antibody
AD2
eFluor™ 450
405/445 nm
Liquid
5 μL/Test
Affinity chromatography
PBS, pH 7.2, with 0.2% BSA
0.09% sodium azide
4° C, store in dark, DO NOT FREEZE!
AB_11041811

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 μL (0.25 μg)/test	35 Publications

### **Product Specific Information**

Description: This AD2 monoclonal antibody reacts with human CD73, a 5'-ectonucleotidase that converts 5'-adenosine monophosphate to adenosine. CD73 is expressed on the surface of endothelial cells, as well as B and T cells, including some CD4+Foxp3+ regulatory T cells. Adenosine production by these cells has been linked to the inhibition of CD4 T cell effector functions such as proliferation and cytokine secretion.

Applications Reported: This AD2 antibody has been reported for use in flow cytometric analysis.


Applications Tested: This AD2 antibody has been pre-titrated and tested by flow cytometric analysis of normal human peripheral blood cells. This can be used at 5 µL (0.25 µg) per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.

eFluor® 450 is an alternative to Pacific Blue®. eFluor® 450 emits at 445 nm and is excited with the Violet laser (405 nm). Please make sure that your instrument is capable of detecting this fluorochome.

Excitation: 405 nm; Emission: 445 nm; Laser: Violet Laser.

Filtration: 0.2 µm post-manufacturing filtered.

## Product Images For CD73 Monoclonal Antibody (AD2), eFluor™ 450, eBioscience™

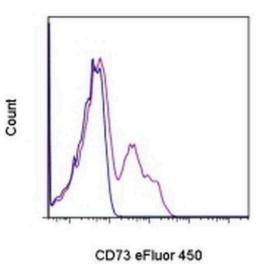

#### CD73 Antibody (48-0739-42) in Flow

Staining of normal human peripheral blood cells with Mouse IgG1 K Isotype Control eFluor® 450 (Product # 48-4714-82) (blue histogram) or Anti-Human CD73 eFluor® 450 (purple histogram). Cells in the lymphocyte gate were used for analysis.

## View more figures on thermofisher.com

#### **□** 35 References

## Flow Cytometry (35)

Frontiers in bioengineering and biotechnology

Impact of Fibronectin Knockout on Proliferation and Differentiation of Human Infrapatellar Fat Pad-Derived Stem Cells.

"Published figure using CD73 monoclonal antibody (Product # 48-0739-42) in Flow Cytometry" Authors: Wang Y,Fu Y,Yan Z,Zhang XB,Pei M

Molecular medicine (Cambridge, Mass.)

# hUCMSCs carrying exenatide prevent T1DM by improving intestinal microflora composition and islet tissue damage repair.

"Published figure using CD73 monoclonal antibody (Product # 48-0739-42) in Flow Cytometry"
Authors: Wang W,Wang Y,Chi J,Tan X,Hu J,Ma X,Sun X,Che K,Lv W,Wang Y

**Year** 2022

**Year** 2022

View more Flow references on thermofisher.com

#### More applications with references on thermofisher.com

For Research Use Only, Not for use in diagnostic procedures. Not for resale without express authorization, Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production obcumentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is imitted to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICUAR PURPOSE, OR NON INFRINGEMENT.

BUYER'S EXCLUSIVE REMEDY FOR NON-CORNORNOR PRODUCTS DURING THE WARRANTY PERIOD IS LIMITAGE. RECEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCTS, AT SELLER'S SOLE OFTION. THERE IS NO BUYER'S EXCLOSED THE NON-CONFORMING PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, or vivo or in vivo therapeutic uses