




# HA Tag Monoclonal Antibody (16B12), Alexa Fluor™ 488

<b>Product Details</b>		
Size	100 μL	
Species Reactivity	Tag	
Published Species	Tag	
Host/Isotype	Mouse / IgG1	
Class	Monoclonal	
Туре	Antibody	
Clone	16B12	
Conjugate	Alexa Fluor™ 488	
Excitation/Emission Max	499/520 nm	
Immunogen	twelve amino acid peptide CYPYDVPDYASL.	
Form	Liquid	
Concentration	1 mg/mL	
Purification	purified	
Storage buffer	PBS, pH 7.2, with 0.1% BSA	
Contains	5mM sodium azide	
Storage conditions	4° C, store in dark	
RRID	AB_2535829	

Applications	Tested Dilution	Publications
Western Blot (WB)	1-10 μg/mL	3 Publications
Immunohistochemistry (IHC)	-	3 Publications
Immunocytochemistry (ICC/IF)	5 μg/mL	13 Publications
Flow Cytometry (Flow)	-	7 Publications
Immunoprecipitation (IP)	-	3 Publications
Miscellaneous PubMed (Misc)	-	4 Publications

## Product Images For HA Tag Monoclonal Antibody (16B12), Alexa Fluor™ 488

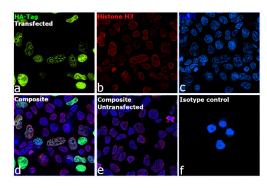

#### HA Tag Antibody (A-21287) in ICC/IF

Immunofluorescent analysis of HA Tag was performed using 70% confluent HA-H3 transfected HEK-293E cells. The cells were fixed with 4% paraformaldehyde for 10 minutes, permeabilized with 0.1% Triton™ X-100 for 15 minutes and blocked with 2% BSA for 1 hour at room temperature. The cells were labeled with HA Tag Monoclonal Antibody (16B12), Alexa Fluor 488 (Product # A-21287) at 5 μg/mL and Histone H3 Rabbit Polyclonal Antibody (Product # 711055) at 0.5 μg /mL in 0.1% BSA, incubated at 4 degree celsius overnight and then labeled with Goat anti-Rabbit IgG (H+L) Highly Cross-Adsorbed Secondary Antibody, Alexa Fluor Plus 647 (Product # A32733) respectively a dilution of 1:2,000 for 45 minutes at room temperature. Panel a (Nuclei: Green) represents HA Tag Monoclonal Antibody. Panel b (Nuclei: Red) represents Histone H3. Panel c (Nuclei: Blue) represents ProLong™ Diamond Antifade Mountant with DAPI (Product # P36962). Panel d represents the merged image showing the colocalization of nuclear signals in transfected cells. Panel e represents untransfected HEK-293E cells. Panel f represents control cells with isotype control antibody to assess background. The images were captured at 60X magnification.

View more figures on thermofisher.com

#### **□** 33 References

### Western Blot (3)

**Nature communications** 

Onsite GTP fuelling via DYNAMO1 drives division of mitochondria and peroxisomes.

"A-21287 was used in Western Blotting to study the role of dynamin-based ring motive-force organiser 1 in mitochondrial division and peroxisome-dividing machinery."

Authors: Imoto Y, Abe Y, Honsho M, Okumoto K, Ohnuma M, Kuroiwa H, Kuroiwa T, Fujiki Y

**Year** 2018

Species Tag

Cell host & microbe

Host Nitric Oxide Disrupts Microbial Cell-to-Cell Communication to Inhibit Staphylococcal Virulence.

"A-21287 was used in Western blotting following immunoprecipitation to demonstrate an anti-virulence mechanism for nitric oxide in innate immunity."

Authors: Urbano R, Karlinsey JE, Libby SJ, Doulias PT, Ischiropoulos H, Warheit-Niemi HI, Liggitt DH, Horswill AR, Fang FC

**Year** 2018

Species Tag

View more WB references on thermofisher.com

## Immunohistochemistry (3)

eLife

A toolbox of nanobodies developed and validated for use as intrabodies and nanoscale immunolabels in mammalian brain neurons.

"A-21287 was used in Immunohistochemistry to report the generation and validation of nAbs against a set of proteins prominently expressed at specific subcellular sites in mammalian brain neurons."

Authors: Dong JX,Lee Y,Kirmiz M,Palacio S,Dumitras C,Moreno CM,Sando R,Santana LF,Südhof TC,Gong B,Murray KD,Trimmer JS

**Year** 2019

Species Tag

The Journal of biological chemistry

Nicotinic acid receptor agonists differentially activate downstream effectors.

"A-21287 was used in immunohistochemistry to study the signaling pathways induced by GPR109A pyrazole agonists that inhibit lipolysis but fail to elicit a flushing response"

Authors: Richman JG,Kanemitsu-Parks M,Gaidarov I,Cameron JS,Griffin P,Zheng H,Guerra NC,Cham L,Maciejewski-Lenoir D,Behan DP,Boatman D,Chen R,Skinner P,Ornelas P,Waters MG,Wright SD,Semple G,Connolly DT

**Year** 2007

Species Tag

View more IHC references on thermofisher.com

## More applications with references on thermofisher.com

ICC/IF (13) Flow (7) IP (3) Misc (4)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production obcumentation specifications and/or accompanying package inserts ("Documentation"), No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is imitted to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any mode of sample is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample inches the product is used to sample inches the product of the product or any product of the documentation accompanying the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, or vito or in vivo therapeutic uses, or any type of consumption to human or animal diagnostic uses, or vivo or in vivo therapeutic uses, or any type of consumption to human or animal to the product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, or any type of consumption to human or animal to human or animal to the purpose.