# CD95 (APO-1/Fas) Monoclonal Antibody (SM1/23), FITC, eBioscience™

Product Details	
Size	100 Tests
Species Reactivity	Human
Published Species	Human
Host/Isotype	Mouse / IgG2b
Recommended Isotype Control	Mouse IgG2b kappa Isotype Control (eBMG2b), FITC, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	SM1/23
Conjugate	FITC
Excitation/Emission Max	498/517 nm
Form	Liquid
Concentration	10 µL/Test
Purification	Affinity chromatography
Storage buffer	TBS, pH 7.2-7.4, with 1% BSA
Contains	0.02% sodium azide
Storage conditions	4° C, store in dark
RRID	AB_10597907

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	Assay-Dependent	2 Publications

#### **Product Specific Information**

Description: BMS140 is suitable as primary antibody for staining in FACS analysis.

BMS140FI can be used for direct staining of cells. BMS140BT can be used in a two-step procedure using the biotin-(strept-) avidin system.

Functional studies:

BMS140 was found to block induction of apoptosis induced by clone SM1/1 or other commercial anti APO-1 antibodies.

Optimal dilutions should be determined by the individual laboratory for each application.


Applications Tested: Flow Cytometry.

Purity: >95%.

Excitation: 488 nm; Emission: 520 nm; Laser: Blue Laser


1

## Product Images For CD95 (APO-1/Fas) Monoclonal Antibody (SM1/23), FITC, eBioscience™

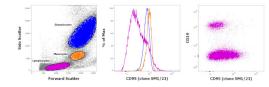

#### CD95 (APO-1/Fas) Antibody (BMS140FI) in Flow

Staining of human peripheral blood mononuclear cells with CD45 Pacific Blue, CD19 PE-Cy7 and CD95 FITC. As expected based on known relative expression patterns, CD95 clone SM1/23 stains a subset of lymphocytes (pink), monocytes (orange) and not granulocytes (blue).

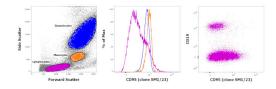

#### CD95 (APO-1/Fas) Antibody (BMS140FI)

Staining of human peripheral blood mononuclear cells with CD45 Pacific Blue, CD19 PE-Cy7 and CD95 FITC. As expected based on known relative expression patterns, CD95 clone SM1/23 stains a subset of lymphocytes (pink), monocytes (orange) and not granulocytes (blue). {RE}

### 2 References

#### Flow Cytometry (2)

Scientific reports The metabolism of cells regulates their sensitivity to NK cells depending on p53 status.	Year 2022 Species Human
"BMS140FI was used in Flow Cytometry to show that treatment of patients with DCA, or a biosimilar drug, could be a clinical option to increase the effectiveness of CAR T cell or allogeneic NK cell therapies."	numan
Authors: Belkahla S,Brualla JM,Fayd'herbe de Maudave A,Falvo P,Allende-Vega N,Constantinides M,Khan AUH, Coenon L,Alexia C,Mitola G,Massa P,Orecchioni S,Bertolini F,Mnif W,Hernandez J,Anel A,Villalba M	
Scientific reports	Year
Metformin sensitizes leukemic cells to cytotoxic lymphocytes by	2022
increasing expression of intercellular adhesion molecule-1 (ICAM-1).	Species
"BMS140FI was used in Flow Cytometry to explore the potential and underlying mechanisms of metaformin in the treatment of hematologic cancers."	Human
Authors: Allende-Vega N,Marco Brualla J,Falvo P,Alexia C,Constantinides M,de Maudave AF,Coenon L,Gitenay D, Mitola G,Massa P,Orecchioni S,Bertolini F,Marzo I,Anel A,Villalba M	

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will minited to incomod or sample. No OTHER WARRANTIES, EXERGIS OR IMPLED, ARERSS OR IMPLED, ARENTES FOR EXPLOY PARTICULAR PURPOSE, OR NON/ERRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR REDUD FOR THE RRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, [III] USE OF OR BY BUYER, [III] USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. USE for many other wine express, instand on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption