



Armenian Hamster IgG Isotype Control (eBio299Arm), Brilliant Ultra Violet™ 661, eBioscience™

Product Details	
Size	50 μg
Host/Isotype	Armenian hamster / IgG
Class	Control
Туре	Isotype Control
Clone	eBio299Arm
Conjugate	Brilliant Ultra Violet™ 661
Excitation/Emission Max	349/659 nm
Form	liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	Assay-dependent	-
Control (Ctrl)	Assay-dependent	-

Product Specific Information

Description: The eBio299Arm monoclonal antibody is useful as an isotype control immunoglobulin.

Applications Reported: This eBio299Arm antibody has been reported for use in flow cytometric analysis.

Applications Tested: This eBio299Arm antibody has been tested by flow cytometric analysis of normal human peripheral blood and mouse splenocytes. Use at the same concentration as the experimental antibody.

Brilliant Ultra Violet™ 661 (BUV661) is a tandem dye that emits at 660 nm and is intended for use on cytometers equipped with an ultraviolet (355 nm) laser. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright, Brilliant Violet[™], Brilliant Ultra Violet[™], or other polymer dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401-42) or Brilliant Stain Buffer[™] (Product # 00-4409-75) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer or Brilliant Stain Buffer for more information.


Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.

Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222-49) (100 μ L of cell sample + 100 μ L of IC Fixation Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333-54) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone-specific performance should be determined empirically.

Excitation: 350 nm; Emission: 660 nm; Laser: Ultraviolet Laser.

BRILLIANT ULTRA VIOLET™ is a trademark or registered trademark of Becton, Dickinson and Company or its affiliates, and is used under license. Powered by Sirigen.™

Product Images For Armenian Hamster IgG Isotype Control (eBio299Arm), Brilliant Ultra Violet™ 661, eBioscience™

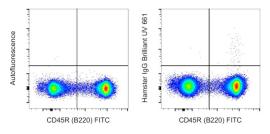

Armenian Hamster IgG Isotype Control (376-4888-81) in Flow C57BL/6 mouse splenocytes were stained with CD45R (B220) Monoclonal Antibody, FITC (Product # 11-0452-82) and staining buffer (autofluorescence) (left) or 0.25 µg of Armenian Hamster IgG Isotype Control, Brilliant Ultra Violet 661 (right). Viable cells in the lymphocyte gate were used for analysis, as determined by 7-AAD (Product # 00-6993-50).

For Research Use Only, Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. No OTHER WARRANTIES, EXPERESS OR IMPLECT LIMITATION, IMPLIED WARRANTIES, EXPERESS OR IMPLECT. IN IMPLIED WARRANTIES, EXPERESS OR IMPLECT. IN IMPLICATION, IMPLIED WARRANTIES, EXPERESS OR IMPLECT. IN IMPLICATION, IMPLIED WARRANTIES, EXPERESS OR IMPLECT. IN IMPLICATION, IMPL