Rat IgG2a kappa Isotype Control (eBR2a), Brilliant Violet[™] 421, eBioscience[™]

Product Details

Size	50 µg
Host/Isotype	Rat / IgG2a, kappa
Class	Monoclonal
Туре	Isotype Control
Clone	eBR2a
Conjugate	Brilliant Violet™ 421
Excitation/Emission Max	406/423 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2921020

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	Assay-Dependent	-
Control (Ctrl)	Assay-Dependent	-

Product Specific Information

Description: The monoclonal rat IgG2a kappa is useful as an isotype control immunoglobulin.

Applications Reported: This eBR2a antibody has been reported for use in flow cytometric analysis.

Applications Tested: This eBR2a antibody has been tested by flow cytometric analysis of mouse splenocytes and normal human peripheral blood cells and should be used at the same concentration as the experimental antibody.

Brilliant Violet[™] 421 (BV421) is a dye that emits at 423 nm and is intended for use on cytometers equipped with a violet (405 nm) laser. Please make sure that your instrument is capable of detecting this fluorochrome.


When using two or more Super Bright, Brilliant Violet[™], Brilliant Ultra Violet[™], or other polymer dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401-42) or Brilliant Stain Buffer[™] (Product # 00-4409-75) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer or Brilliant Stain Buffer for more information.

Excitation: 407 nm; Emission: 423 nm; Laser: Violet Laser.

BRILLIANT VIOLET[™] is a trademark or registered trademark of Becton, Dickinson and Company or its affiliates, and is used under license. Powered by Sirigen.[™].

1

Product Images For Rat IgG2a kappa Isotype Control (eBR2a), Brilliant Violet[™] 421, eBioscience[™]

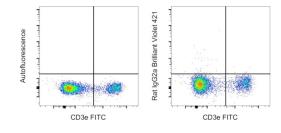

Rat IgG2a kappa Isotype Control (404-4321-81) in Flow

C57BL/6 mouse splenocytes were stained with CD3e Monoclonal Antibody, FITC (Product # 11-0031-82) and staining buffer (autofluorescence) (left) or 1.0 µg of Rat IgG2a kappa Isotype Control, Brilliant Violet 421 (right). Viable cells in the lymphocyte gate were used for analysis, as determined by 7-AAD (Product # 00-6993-50).

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and first accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is performed to effect at the ime of subject to homain, proper and intended usage. This warranty does not extend to anyone other than the buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO THER WARRANTES, EXPRESS OR INNELD, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES OF MERCHANTABILITY, FITNESS FOR ANY PARTICURST DAV PARTICURST DEVENDES, EXPRESS OR INNERINEED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES OF MERCHANTABILITY, FITNESS FOR ANY PARTICURST DAV PA