

IL-2 Monoclonal Antibody (MQ1-17H12), FITC, eBioscience™

Product Details	
Size	100 µg
Species Reactivity	Human
Published Species	Human
Host/Isotope	Rat / IgG2a, kappa
Recommended Isotype Control	Rat IgG2a kappa Isotype Control (eBR2a), FITC, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	MQ1-17H12
Conjugate	FITC
Form	Liquid
Concentration	0.5 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.1% gelatin
Contains	0.09% sodium azide
Storage Conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_465385

Applications	Tested	Dilution	Published
Flow Cytometry (Flow)	✓	1 µg/test	4 Publications

Product Specific Information


Description: The MQ1-17H12 antibody reacts with human interleukin-2 (IL-2), a 17 kDa T cell growth factor and a major immunoregulatory cytokine. The MQ1-17H12 antibody is a non-neutralizing antibody.

Applications Reported: MQ1-17H12 has been reported for use in intracellular flow cytometric analysis.

Applications Tested: This MQ1-17H12 antibody has been tested by intracellular staining and flow cytometric analysis of stimulated normal human peripheral blood cells. This can be used at less than or equal to 1 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Excitation: 488 nm; Emission: 520 nm; Laser: Blue Laser.

Filtration: 0.2 µm post-manufacturing filtered.

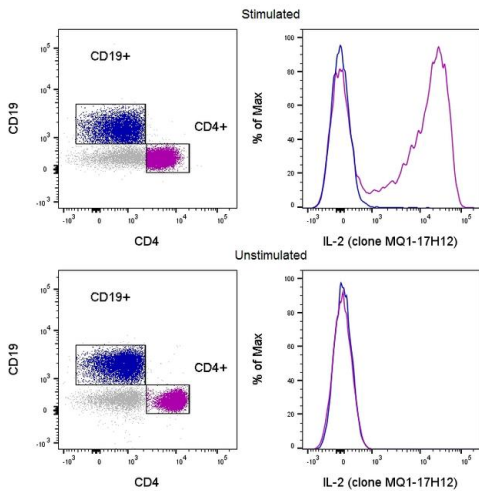

IL-2 Antibody (11-7029-82)

Intracellular staining of stimulated human peripheral blood cells. As expected based on known expression patterns, IL-2 clone MQ1-17H12 stains a subset of CD4+ T cells only after stimulation (top) and does not stain CD19+ B cells regardless of stimulation (top and bottom). Details: Normal human peripheral blood cells were cultured in the presence of Protein Transport Inhibitors (Unstimulated, bottom row) or Cell Stimulation Cocktail (plus protein transport inhibitors, 500X) for 5 hours (Stimulated, top row). Cells were fixed and permeabilized with the IC Fixation and Permeabilization Buffer Set and protocol followed by intracellular staining with CD19 (clone SJ25C1), CD4 (clone RPA-T4), and IL-2 (clone MQ1-17H12). Cells in the CD19+ (blue histogram) or CD4+ (purple histogram) gates were used for analysis. Cell Treatment validation info.

4 References

Flow Cytometry (4)

The Journal of experimental medicine

Anti-TNF drives regulatory T cell expansion by paradoxically promoting membrane TNF-TNF-RII binding in rheumatoid arthritis.

"11-7029 was used in Flow cytometry/Cell sorting to investigate whether TNF boosts or inhibits regulatory T cells (T reg cells), showing that anti-TNF drives regulatory T cell expansion by promoting membrane TNF-TNF-RII binding, in rheumatoid arthritis."

Authors: Nguyen DX, Ehrenstein MR

Species
Human

Dilution
Not Cited

Year
2016

Journal of immunology (Baltimore, Md. : 1950)

Transcriptional profile of tuberculosis antigen-specific T cells reveals novel multifunctional features.

"11-7029 was used in Flow cytometry/Cell sorting to study the role of CD4(+) T cells in latent tuberculosis infection, showing that tuberculosis antigen-specific T cells have novel multifunctional features."

Authors: Arlehamn CL, Seumois G, Gerasimova A, Huang C, Fu Z, Yue X, Sette A, Vijayanand P, Peters B

Species
Human

Dilution
Not Cited

Year
2014

[View more Flow references on thermofisher.com](https://thermofisher.com)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON-INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.