



# IFN gamma Monoclonal Antibody (XMG1.2), PE, eBioscience™

<b>Product Details</b>	
Size	50 μg
Species Reactivity	Mouse
Published Species	Mouse, Human
Host/Isotype	Rat / IgG1, kappa
Recommended Isotype Control	Rat IgG1 kappa Isotype Control (eBRG1), PE, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	XMG1.2
Conjugate	PE
Excitation/Emission Max	565/576 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_466192

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	1 Publication
Immunocytochemistry (ICC/IF)	-	1 Publication
Flow Cytometry (Flow)	0.25 µg/test	299 Publications
ELISA (ELISA)	-	1 Publication
Neutralization (Neu)	-	1 Publication
In Situ Hybridization (ISH) (ISH)	-	2 Publications

#### **Product Specific Information**

Description: The XMG1.2 antibody reacts with mouse interferon (IFN) gamma. The XMG1.2 antibody is a neutralizing antibody. Mouse IFN gamma is a 20 kDa factor produced by activated T, B and NK cells, and is an anti-viral and anti-parasitic cytokine. IFN gamma, in synergy with other cytokines such as TNF alpha, inhibits proliferation of normal and transformed cells. Immunomodulatory effects of IFN gamma are exerted on a wide range of cell types expressing the high affinity receptors for IFN gamma. Glycosylation of IFN gamma does not affect its biological activity.

Applications Reported: This XMG1.2 antibody has been reported for use in intracellular staining followed by flow cytometric analysis.


Applications Tested: This XMG1.2 antibody has been tested by intracellular staining and flow cytometric analysis of stimulated mouse splenocytes using the Intracellular Fixation & Permeabilization Buffer Set (cat. 88-8824) and protocol. Please refer to Best Protocols: Protocol A: Two step protocol for (cytoplasmic) intracellular proteins located under the Resources Tab online. This can be used at less than or equal to 0.25 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell

sample in a final volume of 100  $\mu$ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Excitation: 488-561 nm; Emission: 578 nm; Laser: Blue Laser, Green Laser, Yellow-Green Laser.


Filtration: 0.2 µm post-manufacturing filtered.

# Product Images For IFN gamma Monoclonal Antibody (XMG1.2), PE, eBioscience™

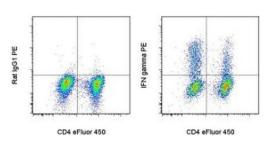

#### IFN gamma Antibody (12-7311-81) in Flow

Mouse splenocytes were stimulated with Con A for 2 days then restimulated with Cell Stimulation Cocktail (plus protein transport inhibitors) (Product # 00-4975-03) for 5 hours. Cells were then surface stained with Anti-Mouse CD4 eFluor® 450 (Product # 48-0042-82) followed by intracellular staining with 0.125 µg of Rat IgG2a K Isotype Control PE (Product # 12-4301-82) (left) or 0.125 µg of Anti-Mouse IFN gamma PE (right) using the Intracellular Fixation & Permeabilization Buffer Set (Product # 88-8824-00) and protocol. Total viable cells, as determined by Fixable Viability Dye eFluor® 506 (Product # 65-0866-14), were used for analysis.

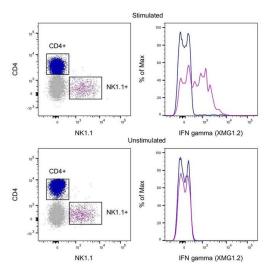

#### **IFN gamma Antibody (12-7311-81)**

Intracellular staining of stimulated mouse splenocytes. As expected based on known expression patterns, IFN gamma clone XMG1.2 stains a minor subset of CD4+ T cells and a larger subset of NK1.1+ NK cells with no staining observed without stimulation. Details: Mouse splenocytes were cultured in the presence of Protein Transport Inhibitors (500X) (Unstimulated, bottom row) or Cell Stimulation Cocktail (plus protein transport inhibitors, 500X) for 5 hours (Stimulated, top row). Cells were fixed and permeabilized with the IC Fixation & Permeabilization Buffer Set and protocol followed by intracellular staining with CD4 (clone RM4-5), NK1.1 (clone PK136) and IFN gamma (clone XMG1.2). Cells in the CD4+ (blue histogram) or NK1.1+ (purple histogram) gates were used for analysis. {TM}

View more figures on thermofisher.com

#### **□ 305 References**

## Immunohistochemistry (1)

**PLoS** pathogens

IL-22 and IDO1 affect immunity and tolerance to murine and human vaginal candidiasis.

"12-7311 was used in Immunohistochemistry to study host resistance and tolerance to the commensal fungus, Candida albicans, during vulvovaginal candidiasis."

Authors: De Luca A,Carvalho A,Cunha C,Iannitti RG,Pitzurra L,Giovannini G,Mencacci A,Bartolommei L,Moretti S, Massi-Benedetti C,Fuchs D,De Bernardis F,Puccetti P,Romani L

**Year** 2014

Species Human Mouse

## Immunocytochemistry (1)

Journal of immunology (Baltimore, Md.: 1950)

Transplantation survival is maintained by granzyme B+ regulatory cells and adaptive regulatory T cells.

"12-7311 was used in Immunocytochemistry to implicate granzyme B expression by regulatory T cells in sustaining long-lived graft survival."

Authors: Gondek DC, Devries V, Nowak EC, Lu LF, Bennett KA, Scott ZA, Noelle RJ

**Year** 2008

Species Mouse

### Flow Cytometry (299)

Therapeutic advances in chronic disease

Prevention of EAE by tolerogenic vaccination with PEGylated antigenic peptides.

"Published figure using IFN gamma monoclonal antibody (Product # 12-7311-82) in Flow Cytometry"

Authors: Pfeil J,Simonetti M,Lauer U,von Thülen B,Durek P,Poulsen C,Pawlowska J,Kröger M,Krähmer R,Leenders F, Hoffmann U Hamann A

**Year** 2023

Cellular and molecular gastroenterology and hepatology

Klebsiella pneumoniae Induces Inflammatory Bowel Disease Through Caspase-11-Mediated IL18 in the Gut Epithelial Cells.

"Published figure using IFN gamma monoclonal antibody (Product # 12-7311-82) in Flow Cytometry" Authors: Zhang Q,Su X,Zhang C,Chen W,Wang Y,Yang X,Liu D,Zhang Y,Yang R

**Year** 2023

View more Flow references on thermofisher.com

#### More applications with references on thermofisher.com

ELISA (1) Neu (1) ISH (2)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization, Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT.

BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT.

BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS OURING THE WARRANTY PERIOD IS LIMITED. AND THE WARRANTY PERIOD IS