CD134 (OX40) Monoclonal Antibody (OX-86), Functional Grade, eBioscience™

Product	Dotails
FIOUUCI	Detalls

i ioddol Dolano	
Size	500 µg
Species Reactivity	Mouse
Published Species	Mouse
Host/Isotype	Rat / IgG1, kappa
Recommended Isotype Control	Rat IgG1 kappa Isotype Control (eBRG1), Functional Grade, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	OX-86
Conjugate	Functional Grade
Form	Liquid
Concentration	1 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	no preservative
Storage conditions	4° C
RRID	AB_469064

Applications	Tested Dilution	Publications
Immunohistochemistry (Frozen) (IHC (F))	-	2 Publications
Flow Cytometry (Flow)	0.5 µg/test	7 Publications
Neutralization (Neu)	-	1 Publication
Functional Assay (FN)	Assay-Dependent	2 Publications

Product Specific Information

Description: The OX-86 monoclonal antibody reacts with mouse CD134, also known as OX40. A member of the TNF receptor superfamily, CD134 is a 50 kDa type I membrane glycoprotein expressed by activated mouse T lymphocytes. Rat CD134 was initially identified as an activation marker only on activated CD4+ T cells. In contrast, mouse CD134 is expressed by both activated CD4+ and CD8+ T cells. The interaction of CD134 with CD252 (OX40 ligand) has been implicated in T cell-dependent humoral responses, regulation of primary T cell expansion, survival of T cells, size of the memory T cell pool, and regulation of tolerance in the CD4+ T cell compartment.

Applications Reported: The OX-86 antibody has been reported for use in flow cytometric analysis. It does not block binding of CD134 to its ligand, but has been reported for use in functional assays.

Applications Tested: The OX-86 antibody has been tested by flow cytometric analysis of activated mouse splenocytes. This can be used at less than or equal to 0.5 μ g per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Storage and handling: Use in a sterile environment.

1


Filtration: 0.2 µm post-manufacturing filtered.

Purity: Greater than 90%, as determined by SDS-PAGE.

Endotoxin Level: Less than 0.001 ng/ μ g antibody, as determined by LAL assay.

Aggregation: Less than 10%, as determined by HPLC.

Product Images For CD134 (OX40) Monoclonal Antibody (OX-86), Functional Grade, eBioscience™

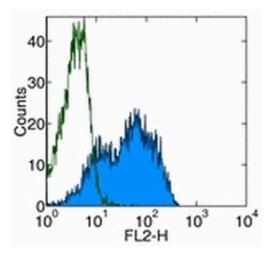

CD134 (OX40) Antibody (16-1341-85) in Flow

Staining of Con A stimulated mouse splenocytes with Anti-Mouse CD134 (OX40) PE.Appropriate isotype controls were used (open histogram).Total viable cells were used for analysis.

View more figures on thermofisher.com

□ 12 References

Immunohistochemistry (Frozen) (2)

Cancer research Adenovirus vector-mediated in vivo gene transfer of OX40 ligand to tumor cells enhances antitumor immunity of tumor-bearing hosts. Authors: Andarini S,Kikuchi T,Nukiwa M,Pradono P,Suzuki T,Ohkouchi S,Inoue A,Maemondo M,Ishii N,Saijo Y, Sugamura K,Nukiwa T	Year 2004
Journal of immunology (Baltimore, Md. : 1950) Regulation of T cell activation in vitro and in vivo by targeting the OX40- OX40 ligand interaction: amelioration of ongoing inflammatory bowel disease with an OX40-IgG fusion protein, but not with an OX40 ligand- IgG fusion protein. Authors: Higgins LM,McDonald SA,Whittle N,Crockett N,Shields JG,MacDonald TT	Year 1999
low Cytometry (7)	
Oncoimmunology	Year

Bendamustine with total body irradiation conditioning yields tolerant T- cells while preserving T-cell-dependent graft-versus-leukemia.	
"Published figure using CD134 (OX40) monoclonal antibody (Product # 16-1341-85) in Flow Cytometry" Authors: Stokes J,Hoffman EA,Molina MS,Kummet N,Simpson RJ,Zeng Y,Katsanis E	
	Year 2021

Notch4 signaling limits regulatory T-cell-mediated tissue repair and promotes severe lung inflammation in viral infections.

"Published figure using CD134 (OX40) monoclonal antibody (Product # 16-1341-85) in Flow Cytometry"

Authors: Harb H,Benamar M,Lai PS,Contini P,Griffith JW,Crestani E,Schmitz-Abe K,Chen Q,Fong J,Marri L,Filaci G,Del Zotto G,Pishesha N,Kolifrath S,Broggi A,Ghosh S,Gelmez MY,Oktelik FB,Cetin EA,Kiykim A,Kose M,Wang Z,Cui Y,Yu XG,Li JZ,Berra L,Stephen-Victor E,Charbonnier LM,Zanoni I,Ploegh H,Deniz G,De Palma R,Chatila TA

View more Flow references on thermofisher.com

More applications with references on thermofisher.com

Neu (1) FN (2)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by property trained individual. Unless otherwise stated in the Documentation, this warranty is provided herein is valid only when used by property trained individual. Unless otherwise stated in the Documentation, this warranty is provided herein is valid only when used by property trained individual. Unless otherwise stated in the Documentation, this warranty is provided herein is valid only when used by property trained individual. Unless otherwise stated in the Documentation, this warranty is provided herein is valid only when used by property trained individual. Unless otherwise stated in the Documentation, the warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample furnished to Buyer is provided therein is valid only MENDE TOR NON-CONFORMING PRODUCTS NOT MONE PROVING TER WARRANTY PERIOD IS NOT MARKEN PROVING TOR NON-CONFORMING PRODUCTS NOT NEE NON-CONFORMING PRODUCTS NOT AND CALUENCE, (III) MISUES, FAULT OR NEELIGUENCE OF OR BY BUYER, (III) DISO FTHE PRODUCTS IN ATMARKEN FOR WHICH THEY WERE NOT DE GONDALE NOT PROVIDENT IN A MANREN FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS IN ATMARKEN FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS IN ATMARKEN FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE FRODUCTS IN ATMARKEN FOR WHICH THEY WERE N