



CD8a Monoclonal Antibody (53-6.7), APC, eBioscience™

Product Details	
Size	100 μg
Species Reactivity	Mouse
Published Species	Mouse, Human
Host/Isotype	Rat / IgG2a, kappa
Recommended Isotype Control	Rat IgG2a kappa Isotype Control (eBR2a), APC, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	53-6.7
Conjugate	APC
Excitation/Emission Max	651/660 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_469335

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	1 Publication
Immunohistochemistry (Paraffin) (IHC (P))	-	1 Publication
Immunocytochemistry (ICC/IF)	-	1 Publication
Flow Cytometry (Flow)	0.125 µg/test	207 Publications
Functional Assay (FN)	-	3 Publications
Miscellaneous PubMed (Misc)	-	1 Publication

Product Specific Information

Description: The 53-6.7 monoclonal antibody reacts with the mouse CD8a molecule. CD8a is an approximately 32-34 kDa cell surface receptor expressed either as a heterodimer with the CD8 beta chain (CD8 alpha beta) or as a homodimer (CD8 alpha alpha). A majority of thymocytes and a subpopulation of mature alpha beta TCR T cells express CD8 alpha beta while gamma delta TCR T cells, a subpopulation of intestinal intraepithelial lymphocytes (IELs) and dendritic cells express CD8 alpha alpha. CD8 binds to MHC class I and through its association with protein tyrosine kinase p56lck plays a role in T cell development and activation of mature T cells.


Applications Reported: The 53-6.7 antibody has been reported for use in flow cytometric analysis.


Applications Tested: The 53-6.7 antibody has been tested by flow cytometric analysis of mouse thymocyte or splenocyte suspensions. This can be used at less than or equal to 0.125 μ g per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Excitation: 633-647 nm; Emission: 660 nm; Laser: Red Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For CD8a Monoclonal Antibody (53-6.7), APC, eBioscience™

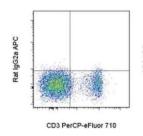

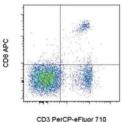

CD8a Antibody (17-0081-82) in Flow

Staining of mouse splenocytes with Anti-Mouse CD3 PerCP-eFluor® 710 (Product # 46-0032-82) and 0.06 μ g of Rat IgG2a K Isotype Control APC (Product # 17-4321-81) (left) or 0.06 μ g of Anti-Mouse CD8a APC (right). Cells in the lymphocyte gate were used for analysis.

□ 214 References

Immunohistochemistry (1)

European journal of immunology

Nitric oxide inhibits the accumulation of CD4+CD44hiTbet+CD69lo T cells in mycobacterial infection.

"17-0081 was used in Immunohistochemistry to support the hypothesis that effector T cells in mycobacterial granulomata are not a uniform effector population but exist in distinct subsets with differential susceptibility to the regulatory effects of nitric oxide."

Authors: Pearl JE, Torrado E, Tighe M, Fountain JJ, Solache A, Strutt T, Swain S, Appelberg R, Cooper AM

Year 2012

Species Mouse

Dilution 1:200

Immunohistochemistry (Paraffin) (1)

European journal of medical research

The inflammatory phase of fracture healing is influenced by oestrogen status in mice.

"17-0081 was used in Immunohistochemistry to indicate that oestrogen-deficiency influences the early inflammatory phase after fracture."

Authors: Haffner-Luntzer M, Fischer V, Prystaz K, Liedert A, Ignatius A

Year 2017

Species Mouse

Dilution 1:800

Immunocytochemistry (1)

Journal of immunology (Baltimore, Md.: 1950)

Programmed death-1 pathway in host tissues ameliorates Th17/Th1-mediated experimental chronic graft-versus-host disease.

"17-0081 was used in Flow cytometry/Cell sorting to clarify that programed death-ligand 1 deficiency in host tissues, not hematopoietic cells, is associated with exacerbated chronic graft-versus-host disease."

Authors: Fujiwara H,Maeda Y,Kobayashi K,Nishimori H,Matsuoka K,Fujii N,Kondo E,Tanaka T,Chen L,Azuma M,Yagita H,Tanimoto M

Year 2014

Species Mouse

Flow Cytometry (207)

Cancer immunology, immunotherapy: CII

Egr2 and 3 maintain anti-tumour responses of exhausted tumour infiltrating CD8+T cells.

"17-0081-82 was used in Flow cytometry/Cell sorting to conclude that Egr2 and 3 are important for maintaining anti-tumour responses of exhausted CD8 + TILs."

Authors: Symonds ALJ,Miao T,Busharat Z,Li S,Wang P

Year 2023

Species Mouse

View more Flow references on thermofisher.com

More applications with references on thermofisher.com

FN (3)

Misc (1)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production of source and a specifications and/or accompanying package inserts ("Documentation"), No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is immitted to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample invited to sample furnished to Buyer is merely illustrative of the general pand quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT.
BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS) AT SELLER'S SOLE OPTION. THERE IS NO DELIGATION TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS) AT SELLER'S SOLE OPTION. THERE IS NO DELIGATION TO REPAIR, REPLACE OR SOLE OPTION. THERE IS NO DELIGATION TO REPAIR, REPLACE OR SOLE OPTION. THERE IS NO DELIGATION TO REPAIR, REPLACE OR SOLE OPTION. THE PRODUCTS IN AMAINER FOR WHICH THEY WERE FOR TO TESICIARE FOR EACH SOLE OPTION. THERE IS NO DELIGATION TO REPAIR, REPLACE OR SOLE OPTION. THERE IS NO DELIGATION TO REPAIR, REPLACE OR SOLE OPTION. THERE IS NO DELIGATION TO REPAIR, REPLACE OR SOLE OPTION. THERE IS NO DELIGATION TO REPAIR, REPLACE OR SOLE OPTION. THERE IS NO DELIGATION TO REPAIR STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the doc