

# Phospho-Syk (Tyr348) Monoclonal Antibody (moch1ct), APC, eBioscience™

Product Details	
Size	100 Tests
Species Reactivity	Human, Mouse
Published Species	Human, Mouse
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), APC, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	moch1ct
Conjugate	APC
Excitation/Emission Max	651/660 nm
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.2% BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2573276

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 µL (0.06 µg)/test	4 Publications

## Product Specific Information

**Description:** This moch1ct monoclonal antibody recognizes human and mouse spleen tyrosine kinase (also known as SYK) when phosphorylated on tyrosine 348 (Y348). SYK is the founding member of the SYK family of kinases that also includes ZAP-70 (zeta-associated protein of 70 kD) and is expressed in hematopoietic cells, including B lymphocytes, immature (CD4, CD8 double-negative and double-positive) thymocytes, and myeloid cells, epithelial cell lines, and normal breast tissue. SYK is critical for B cell receptor (BCR) signaling and B cell development. Autophosphorylation of at Y348 is necessary for SYK to become fully catalytically active and creates a docking site for the SH2 domains of Vav, Grb2, p85 subunit of PI3 kinase, and PLC gamma.

Specificity of this moch1ct clone was confirmed by ELISA, flow cytometry, and western blotting.

**Applications Reported:** This moch1ct antibody has been reported for use in intracellular staining followed by flow cytometric analysis.

**Applications Tested:** This moch1ct antibody has been pre-titrated and tested by intracellular staining followed by flow cytometric analysis of stimulated mouse splenocytes. This can be used at 5 µL (0.06 µg) per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10<sup>5</sup> to 10<sup>8</sup> cells/test.


**Staining Protocol:** All protocols work well for this monoclonal antibody. Use of Protocol A: Two-step protocol: intracellular

(cytoplasmic) proteins allows for the greatest flexibility for detection of surface and intracellular (cytoplasmic) proteins. Use of Protocol B: One-step protocol: intracellular (nuclear) proteins is recommended for staining of transcription factors in conjunction with surface and phosphorylated intracellular (cytoplasmic) proteins. Protocol C: Two-step protocol: Fixation /Methanol allows for the greatest discrimination of phospho-specific signaling between unstimulated and stimulated samples, but with limitations on the ability to stain specific surface proteins (refer to "Clone Performance Following Fixation /Permeabilization" located in the Best Protocols Section under the Resources tab online). All Protocols can be found in the Flow Cytometry Protocols: "Staining Intracellular Antigens for Flow Cytometry Protocol" located in the Best Protocols Section under the Resources tab online.

Excitation: 633-647 nm; Emission: 660 nm; Laser: Red Laser.

Filtration: 0.2 µm post-manufacturing filtered.

**Product Images For Phospho-Syk (Tyr348) Monoclonal Antibody (moch1ct), APC, eBioscience™**

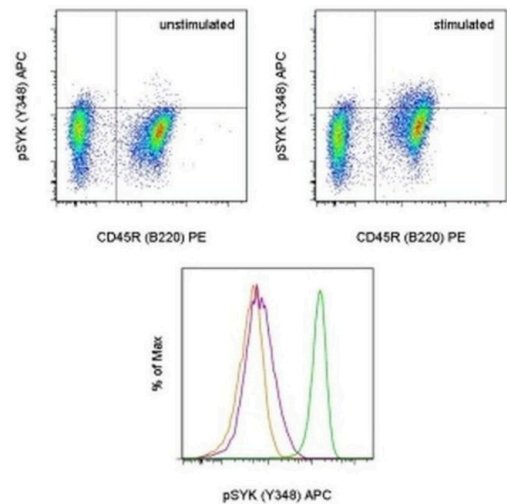

**Phospho-Syk (Tyr348) Antibody (17-9014-42) in Flow**  
TOP: Intracellular staining of unstimulated (left) or 2-minute F (ab')<sub>2</sub> Anti-Mouse IgM, u chain specific Functional Grade Purified-stimulated (Product # 16-5092-85) (right) mouse spleen cells with Anti-Mouse CD45R (B220) PE (Product # 12-0452-82) and Anti-Human/Mouse phospho-SYK (Y348) APC using the Intracellular Fixation and Permeabilization Buffer Set (Product # 88-8824-00) and protocol. BOTTOM: Intracellular staining of unstimulated (orange histogram), 2-minute F (ab')<sub>2</sub> Anti-Mouse IgM, u chain specific Functional Grade Purified-stimulated (Product # 16-5092-85) (purple histogram), or hydrogen peroxide-activated sodium pervanadate-treated mouse spleen cells (green histogram) with Anti-Human/Mouse phospho-SYK (Y348) APC. CD45R (B220)+ lymphocytes were used for analysis.

Flow Cytometry (4)

<p>Frontiers in immunology</p> <p><b>Complement Factor H Modulates Splenic B Cell Development and Limits Autoantibody Production.</b></p> <p>"17-9014-42 was used in Flow Cytometry to indicate that CFH, through its function as a complement repressor, acts as a negative regulator of BCR signaling and limits autoimmunity."</p> <p>Authors: Kiss MG,Ozsvár-Kozma M,Porsch F,Göderle L,Papac-Milievi N,Bartolini-Gritti B,Tsiantoulas D,Pickering MC,Binder CJ</p>	<p>Year 2020</p> <p>Species Mouse</p>
<p>Northern clinics of Istanbul</p> <p><b>IL-15 negatively regulates curdlan-induced IL-23 production by human monocyte-derived dendritic cells and subsequent Th17 response.</b></p> <p>"17-9014 was used in Flow cytometry/Cell sorting to assess the effects of long- and short-term IL-15 cytokine exposure of human monocyte-derived curdlan-matured dendritic cells (DCs) on the production of Th17 cell-polarizing cytokine IL-23 and subsequent Th17 cell activation."</p> <p>Authors: Eken A,Okus Z,Erdem S,Azizoglu ZB,Haliloglu Y,Bicer A,Gur TN,Yilmaz E,Karakukcu M,Altuntas HD,Canatan H</p>	<p>Year 2020</p> <p>Species Human</p>

[View more Flow references on thermofisher.com](#)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.