Performanc guarenteed

LAP (Latency Associated peptide) Monoclonal Antibody (FNLAP), APC, eBioscience[™]

Product Details	
Size	100 Tests
Species Reactivity	Human
Published Species	Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), APC, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	FNLAP
Conjugate	APC
Excitation/Emission Max	651/660 nm
Immunogen	E. coli derived hLAP-TGFb, mixed with VB3A9
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.2% BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2573316

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 μL (0.5 μg)/test	6 Publications

Product Specific Information

Description: The FNLAP monoclonal antibody reacts with human latency associated peptide (LAP, pro-TGF beta 1, LAP/TGF beta 1). Many different cells produce TGF beta and it mediates effects on the proliferation, differentiation and function of many cell types. TGF beta is synthesized as a precursor that contains LAP at the N-terminus and mature TGF beta at the C-terminus. Processing and cleavage of the precursor protein between amino acids 278 and 279 results in the formation of LAP dimers and TGF beta dimers that then non-covalently associate with each other to form the small latent TGF beta complex. LAP is secreted and can be found in the extracellular matrix. In addition, LAP can also be expressed on platelets and activated regulatory T cells. It is believed that this surface-expressed LAP is due to the binding of LAP to GARP (LRRC32), which is a transmembrane protein that is also found at high levels on platelets and activated regulatory T cells.

Applications Reported: This FNLAP antibody has been reported for use in flow cytometric analysis.


Applications Tested: This FNLAP antibody has been pre-titrated and tested by flow cytometric analysis of stimulated normal human peripheral blood cells using the Foxp3/Transcription Factor Staining Buffer Set (cat. 00-5523) and protocol. Please refer to Best Protocols: Protocol B: One step protocol for (nuclear) intracellular proteins located under the Resources Tab online. This can be used at 5 μ L (0.5 μ g) per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.

1

Excitation: 633-647 nm; Emission: 660 nm; Laser: Red Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For LAP (Latency Associated peptide) Monoclonal Antibody (FNLAP), APC, eBioscience™

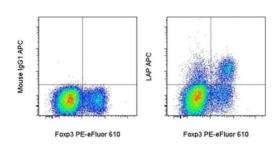

LAP (Latency Associated peptide) Antibody (17-9829-42) in Flow

Normal human peripheral blood cells were stimulated overnight with Human IL-2 Recombinant Protein (Product # 14-8029-81), Anti-Human CD3, and Anti-Human CD28 Functional Grade Purifieds (Product # 16-0037-81 and Product # 16-0289-81). These cells were then stained with Fixable Viability Dye eFluor® 450 (Product # 65-0863-14), followed by surface staining with Mouse IgG1 K Isotype Control APC (Product # 17-4714-81) (left) or Anti-Human LAP (Latency Associated Peptide) APC (right), then intracellular staining with Anti-Human Foxp3 PE-eFluor® 610 (Product # 61-4776-42) using the Foxp3/Transcription Factor Staining Buffer Set (Product # 00-5523-00) and protocol. Viable CD4+ lymphocytes were used for analysis.

View more figures on thermofisher.com

□ 6 References

Flow Cytometry (6)

eLife Dendritic cell Piezo1 directs the differentiation of T _H 1 and T _{reg} cells in	Year 2022
cancer.	
"Published figure using LAP (Latency Associated peptide) monoclonal antibody (Product # 17-9829-42) in Flow Cytometry"	
Authors: Wang Y,Yang H,Jia A,Wang Y,Yang Q,Dong Y,Hou Y,Cao Y,Dong L,Bi Y,Liu G	
JCI insight	Year 2019
Systematic testing and specificity mapping of alloantigen-specific chimeric antigen receptors in regulatory T cells.	
"Published figure using LAP (Latency Associated peptide) monoclonal antibody (Product # 17-9829-42) in Flow Cytometry"	
Authors: Dawson NA,Lamarche C,Hoeppli RE,Bergqvist P,Fung VC,McIver E,Huang Q,Gillies J,Speck M,Orban PC, Bush JW,Mojibian M,Levings MK	

View more Flow references on thermofisher.com

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is more of the product set of any product will conform to such model or sample function of asample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. No OTHER WARRANTES, EXPRESS OR INVELD, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES OF MERCHANTABILITY, FITNESS FOR ANY PARTICURCIS DURING THE WARRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES SOL ANY PARTICURCIS DURING THE WARRANTE REPRIND IS LIMITED TO REPARE, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS OURING THE WARRANT FES OL DEPTION DISAGENON TO REPERING PERIOD IS LIMITED TO REPARE, REPLACE OR OR BY BUYER, (III) USE OF THE PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE of THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDULGS IN THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE COMpanying the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commencial uses, in vitro diagnostic uses, arvivo or in vivo therapeutic uses, orany type of consumption to human or animals.

2