# TNF alpha Monoclonal Antibody (MAb11), PE-Cyanine7, eBioscience™

Prod	uct	Details
FIOU	uci	Delans

r roddor Details	
Size	100 µg
Species Reactivity	Human
Published Species	Non-human primate, Human, Rhesus monkey
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), PE-Cyanine7, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	MAb11
Conjugate	PE-Cyanine7
Excitation/Emission Max	569/780 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_469686

Applications	Tested Dilution	Publications
Immunocytochemistry (ICC/IF)	-	3 Publications
Flow Cytometry (Flow)	0.5 µg/test	41 Publications

# **Product Specific Information**

Description: The MAb11 antibody reacts with human tumor necrosis factor-alpha (TNF alpha), a 17 kDa cytokine produced by monocytes, macrophages, neutrophils, NK cells and CD4+ T cells. TN alpha has cytolytic activity against a range of tumor cells and is important in immune regulation. TNF alpha forms dimers and trimers and also exists as a 26 kDa membrane bound form. The MAb11 antibody is a neutralizing antibody.

Applications Reported: This MAb11 antibody has been reported for use in intracellular staining followed by flow cytometric analysis.

Applications Tested: This MAb11 antibody has been tested by intracellular staining and flow cytometric analysis of stimulated human peripheral blood cells using the Intracellular Fixation & Permeabilization Buffer Set (cat. 88-8824) and protocol. Please refer to Best Protocols: Protocol A: Two step protocol for (cytoplasmic) intracellular proteins located under the Resources Tab online.

This MAb11 antibody is offered in 2 formats:-  $\mu$ g size: This can be used at less than or equal to 0.5  $\mu$ g per test. A test is defined as the amount ( $\mu$ g) of antibody that will stain a cell sample in a final volume of 100  $\mu$ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest. - test size: has been pre-titrated and tested by intracellular staining and flow

1

cytometric analysis of stimulated human peripheral blood cells. This can be used at 5  $\mu$ L (0.25  $\mu$ g) per test. A test is defined as the amount ( $\mu$ g) of antibody that will stain a cell sample in a final volume of 100  $\mu$ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.


Light sensitivity: This tandem dye is sensitive photo-induced oxidation. Please protect this vial and stained samples from light.

Fixation: Samples can be stored in IC Fixation Buffer (cat. 00-8222) (100 µL cell sample + 100 µL IC Fixation Buffer) or 1-step Fix/Lyse Solution (cat. 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency /compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 488-561 nm; Emission: 775 nm; Laser: Blue Laser, Green Laser, Yellow-Green Laser.

Filtration: 0.2 µm post-manufacturing filtered.

# Product Images For TNF alpha Monoclonal Antibody (MAb11), PE-Cyanine7, eBioscience™

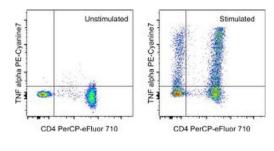

## TNF alpha Antibody (25-7349-82) in Flow

Normal human peripheral blood cells were unstimulated (left) or stimulated for 5 hours with Cell Stimulation Cocktail (plus protein transport inhibitors) (Product # 00-4975-03) (right). Cells were intracellularly stained with Anti-Human CD4 PerCP-eFluor® 710 (Product # 46-0047-42) and Anti-Human TNF alpha PE-Cyanine7 using the Intracellular Fixation & Permeabilization Buffer Set (Product # 88-8824-00) and protocol. Cells in the lymphocyte gate were used for analysis.

View more figures on thermofisher.com

# **44** References

# Immunocytochemistry (3)

#### Cell reports. Medicine

# Targeted T cell receptor gene editing provides predictable T cell product function for immunotherapy.

"25-7349-82 was used in Immunocytochemistry to conclude that increased T cell product homogeneity after targeted TCR gene editing correlates with predictable in vivo T cell responses, which represents a crucial aspect for clinical application in adoptive T cell immunotherapy."

Authors: Müller TR,Jarosch S,Hammel M,Leube J,Grassmann S,Bernard B,Effenberger M,Andrä I,Chaudhry MZ, Käuferle T,Malo A,Cicin-Sain L,Steinberger P,Feuchtinger T,Protzer U,Schumann K,Neuenhahn M,Schober K,Busch DH

#### Frontiers in immunology

# PD-1 Expression on *Mycobacterium tuberculosis*-Specific CD4 T Cells Is Associated With Bacterial Load in Human Tuberculosis.

"25-7349 was used in Immunocytochemistry to study the effect of PD-1 expression on bacterial load and functional capacity of Mycobacterial tuberculosis specific T cell responses."

Authors: Day CL, Abrahams DA, Bunjun R, Stone L, de Kock M, Walzl G, Wilkinson RJ, Burgers WA, Hanekom WA

### **Year** 2021

Species Human

Year 2019

Species Human

View more ICC/IF references on thermofisher.com

Flow	Cytometry	(41)
------	-----------	------

Nature communications PRDM1/BLIMP1 induces cancer immune evasion by modulating the USP22-SPI1-PD-L1 axis in hepatocellular carcinoma cells.	<b>Year</b> 2022
Published figure using TNF alpha monoclonal antibody (Product # 25-7349-82) in Flow Cytometry"	
Authors: Li Q,Zhang L,You W,Xu J,Dai J,Hua D,Zhang R,Yao F,Zhou S,Huang W,Dai Y,Zhang Y,Baheti T,Qian X,Pu L, Ku J,Xia Y,Zhang C,Tang J,Wang X	
	Year
Frontiers in immunology Safety and Immunogenicity of Malaria Vectored Vaccines Given with	Year 2022

"25-7349 was used in Flow cytometry/Cell sorting to demonstrate malaria vectored prime-boost vaccines coadministered with routine childhood immunizations were well tolerated. Potent humoral and cellular immunity induced by ChAd63 MVA ME-TRAP did not reduce the immunogenicity of co-administered EPI vaccines, supporting further evaluation of this regimen in infant populations."

Authors: Mensah VA,Roetynck S,Kanteh EK,Bowyer G,Ndaw A,Oko F,Bliss CM,Jagne YJ,Cortese R,Nicosia A,Roberts R,D'Alessio F,Leroy O,Faye B,Kampmann B,Cisse B,Bojang K,Gerry S,Viebig NK,Lawrie AM,Clarke E,Imoukhuede EB, Ewer KJ,Hill AVS,Afolabi MO

View more Flow references on thermofisher.com

Dilution

1:500

### More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in offect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications represent that any provided herein is valid only when used by property trained individuals. Unless otherwise stated in the Documentation, this warranty are rrow date of a simple furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample functions to general injustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample innovable WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFERNING PRODUCTS DURING THE WARRANT PERIOD IS LIMITED TO REPAR. REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS DURING THE WARRANT PERIOD IS LIMITED TO REPAR. REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OPTION. THERE IS NO BUILGATION TO REPAR. REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS DURING THE WARRANT PERIOD IS LIMITED TO REPAR. REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OPTION. THERE IS NO BUILGATION TO REPAR. REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OPTION. THERE IS NO BUILGATION TO REPAR. REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorizad commercial uses, in vitro diagnosciu. sex, or vivo or in vivo there