Performanc guarenteed

CD8a Monoclonal Antibody (RPA-T8), PE-Cyanine5.5, eBioscience™

Product Details	
Size	100 Tests
Species Reactivity	Human
Published Species	Rat, Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), PE-Cyanine5.5, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	RPA-T8
Conjugate	PE-Cyanine5.5
Excitation/Emission Max	564/692 nm
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.2% BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_11218701

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 μL (0.125 μg)/test	26 Publications

Product Specific Information

Description: The RPA-T8 monoclonal antibody reacts with the human CD8a molecule, an approximately 32-34 kDa cell surface receptor expressed either as a heterodimer with the CD8 beta chain (CD8 alpha/beta) or as a homodimer (CD8 alpha/alpha). A majority of thymocytes and a subpopulation of mature T cells and NK cells express CD8a. CD8 binds to MHC class I and through its association with protein tyrosine kinase p56lck plays a role in T-cell development and activation of mature T cells.

Applications Reported: This RPA-T8 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This RPA-T8 antibody has been pre-titrated and tested by flow cytometric analysis of normal human peripheral blood cells. This can be used at 5 μ L (0.125 μ g) per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.

Light sensitivity: This tandem dye is sensitive photo-induced oxidation. Please protect this vial and stained samples from light.


Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-822-49) (100 µL cell sample + 100 µL IC Fixation Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333-54) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone

1

specific performance should be determined empirically.

Excitation: 488-561 nm; Emission: 695 nm; Laser: Blue Laser, Green Laser, Yellow-Green Laser.

Product Images For CD8a Monoclonal Antibody (RPA-T8), PE-Cyanine5.5, eBioscience™

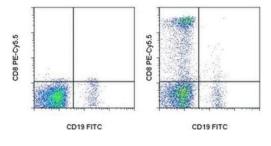

CD8a Antibody (35-0088-42) in Flow

Staining of normal human peripheral blood cells with Anti-Human CD19 FITC (Product # 11-0199-42) and Mouse IgG1 K Isotype Control PE-Cyanine5-5 (Product # 35-4714-82) (left) or Anti-Human CD8a (right) PE-Cyanine5-5. Cells in the lymphocyte gate were used for analysis.

View more figures on thermofisher.com

26 References

Flow Cytometry (26)

Frontiers in immunology Bispecific T-Cell Engaging Antibodies Against MUC16 Demonstrate	Year 2021
Efficacy Against Ovarian Cancer in Monotherapy and in Combination	
With PD-1 and VEGF Inhibition.	
"Published figure using CD8a monoclonal antibody (Product # 35-0088-42) in Flow Cytometry" Authors: Yeku OO,Rao TD,Laster I,Kononenko A,Purdon TJ,Wang P,Cui Z,Liu H,Brentjens RJ,Spriggs D	
Autrors. Teku OO,Rao TD,Laster I,Kononenko A,Purdon TJ,Wang P,Gui Z,Liu H,Bienijens KJ,Spinggs D	
Oncotarget	Year

Multipeptide stimulated PBMCs generate T_{FM}/T_{CM} for adoptive cell

therapy in multiple myeloma.

"Published figure using CD8a monoclonal antibody (Product # 35-0088-42) in Flow Cytometry" Authors: Vardam-Kaur T,Pathangey LB,McCormick DJ,Bergsagel PL,Cohen PA,Gendler SJ

View more Flow references on thermofisher.com

2021

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty pervised herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this varranty ger from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample furnished to Buyer is morely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample furnished. The WARRANTES, EXPRESS OR INNELD, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFERINGEMENT. BUYER'S EXCLUSIVE REMEMPT FOR NON-CONFERNING PRODUCTS DURING THE WARRANT PERIOD IS LIMITED TO REPAIR. REFLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS (S) AT SELLER'S OL DE OTION. THERE IS NO BULGATION TO REFRIR REFLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OPTION. THERE WARRANT PARTICULOR NEGLIGENCE OF OR BY BUYER'S SOLE OPTION. THERE WARRANT CONFORMATION TO REFLAR. REFLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER', (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (V) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product is intende