CD3 Monoclonal Antibody (SK7), Brilliant Ultra Violet[™] 496, eBioscience[™]

Product Details	
Size	100 Tests
Species Reactivity	Chimpanzee, Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), Brilliant Ultra Violet™ 496, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	SK7
Conjugate	Brilliant Ultra Violet™ 496
Excitation/Emission Max	348/493 nm
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2925315

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 μL (1.0 μg)/test	-

Product Specific Information

Description: The SK7 monoclonal antibody reacts with human and chimpanzee CD3e, a 20 kDa subunit of the TCR complex. Along with the other CD3 subunits gamma and delta, the epsilon chain is required for proper assembly, trafficking and surface expression of the TCR complex. CD3 is expressed by thymocytes in a developmentally regulated manner and by all mature T cells. The SK7 and UCHT1 monoclonal antibodies cross-block binding, suggesting recognition of overlapping epitope. In contrast clones OKT3 and SK7 see different epitopes.

The antibody SK7 recognizes chimpanzee CD3.

Applications Reported: This SK7 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This SK7 antibody has been pre-diluted and tested by flow cytometric analysis of normal human peripheral blood cells. This may be used at 5 μ L (1.0 μ g) per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells /test.

Brilliant Ultra Violet[™] 496 (BUV496) is a tandem dye that emits at 496 nm and is intended for use on cytometers equipped with an ultraviolet (355 nm) laser. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright, Brilliant Violet[™], Brilliant Ultra Violet[™], or other polymer dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401-42) or Brilliant Stain Buffer[™] (Product # 00-4409-75) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super

Bright Staining Buffer or Brilliant Stain Buffer for more information.


Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.

Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222-49) (100 µL of cell sample + 100 µL of IC Fixation Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333-54) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone-specific performance should be determined empirically.

Excitation: 350 nm; Emission: 496 nm; Laser: Ultraviolet Laser.

BRILLIANT ULTRA VIOLET[™] is a trademark or registered trademark of Becton, Dickinson and Company or its affiliates, and is used under license. Powered by Sirigen.[™]

Product Images For CD3 Monoclonal Antibody (SK7), Brilliant Ultra Violet[™] 496, eBioscience[™]

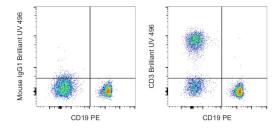

CD3 Antibody (364-0036-42) in Flow

Normal human peripheral blood cells were stained with CD19 Monoclonal Antibody, PE (Product # 12-0199-42) and Mouse IgG1 kappa Isotype Control, Brilliant Ultra Violet 496 (Product # 364-4714-81) (left) or CD3 Monoclonal Antibody, Brilliant Ultra Violet 496 (right). Viable cells in the lymphocyte gate were used for analysis, as determined by 7-AAD (Product # 00-6993-50).

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and general to perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, this warranty is provided herein is valid only when used by properly trained individuals. Duless otherwise stated in the Documentation, this warranty a long year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTES, EXPRESS OR INNELD, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES OF MERCHANTABULTY, FITNESS FOR ANY PARTICURS DURING THE WARRANT PERIOD IS LIMITED TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OFTION. THERE IS NO OBLICATION TO REPARIN, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OFTION. THERE IS NO BULGATION TO FERRAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLICENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product to human or animals.