



CD170 (Siglec F) Monoclonal Antibody (1RNM44N), Brilliant Ultra Violet™ 563, eBioscience™

Product Details	
Size	25 μg
Species Reactivity	Mouse
Host/Isotype	Rat / IgG2a, kappa
Recommended Isotype Control	Rat IgG2a kappa Isotype Control (eBR2a), Brilliant Ultra Violet™ 563, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	1RNM44N
Conjugate	Brilliant Ultra Violet™ 563
Excitation/Emission Max	349/563 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2925388

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	0.5 µg/test	-

Product Specific Information

Description: This 1RNM44N monoclonal antibody recognizes mouse CD170, also known as Siglec F. Siglec F is a cell surface lectin belonging to the Ig superfamily that binds glycoconjugates containing sialic acids that are commonly found on various cell types. The cytoplasmic domain of Siglec F contains an immunoreceptor tyrosine-based inhibitory motif (ITIM) that initiates intracellular signaling upon ligand binding. Siglec F is expressed mostly on eosinophils and alveolar macrophages and lower levels of this receptor have also been reported on immature myelomonocytic cells. mouse Siglec F is a functional ortholog of human Siglec 8, however, unlike human Siglec 8, mouse Siglec F is not expressed on mast cells.

This 1RNM44N antibody will recognize a formaldehyde-fixed epitope.

Applications Reported: This 1RNM44N antibody has been reported for use in flow cytometric analysis.

Applications Tested: This 1RNM44N antibody has been tested by flow cytometric analysis of mouse thioglycolate-elicited peritoneal exudate cells. This may be used at less than or equal to 0.5 μ g per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Brilliant Ultra Violet[™] 563 (BUV563) is a tandem dye that emits at 564 nm and is intended for use on cytometers equipped with an ultraviolet (355 nm) laser. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright, Brilliant Violet™, Brilliant Ultra Violet™, or other polymer dye-conjugated antibodies in a

staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401-42) or Brilliant Stain Buffer™ (Product # 00-4409-75) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer or Brilliant Stain Buffer for more information.


Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.

Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222-49) ($100~\mu$ L of cell sample + $100~\mu$ L of IC Fixation Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333-54) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone-specific performance should be determined empirically.

Excitation: 350 nm; Emission: 564 nm; Laser: Ultraviolet Laser.

BRILLIANT ULTRA VIOLET™ is a trademark or registered trademark of Becton, Dickinson and Company or its affiliates, and is used under license. Powered by Sirigen™.

Product Images For CD170 (Siglec F) Monoclonal Antibody (1RNM44N), Brilliant Ultra Violet™ 563, eBioscience™

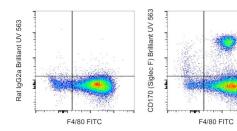

CD170 (Siglec F) Antibody (365-1702-80) in Flow

BALB/c mouse thioglycolate-elicited peritoneal exudate cells were stained with F4 /80 Monoclonal Antibody, FITC (Product # 11-4801-82) and 0.25 µg of Rat IgG2a kappa Isotype Control, Brilliant Ultra Violet 563 (Product # 365-4321-81) (left) or 0.25 µg of CD170 (Siglec F) Monoclonal Antibody, Brilliant Ultra Violet 563 (right). Total viable cells were used for analysis, as determined by 7-AAD (Product # 00-6993-50).