CD279 (PD-1) Monoclonal Antibody (J43), Brilliant Ultra Violet[™] 805, eBioscience[™]

Product Details	
Size	25 µg
Species Reactivity	Mouse
Host/Isotype	Armenian hamster / IgG
Recommended Isotype Control	Armenian Hamster IgG Isotype Control (eBio299Arm), Brilliant Ultra Violet™ 805, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	J43
Conjugate	Brilliant Ultra Violet™ 805
Excitation/Emission Max	349/804 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2896178

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	0.5 µg/test	-

Product Specific Information

Description: The J43 monoclonal antibody reacts with mouse PD-1 (programmed death-1), a 55 kDa member of the Ig superfamily. PD-1 contains the immunoreceptor tyrosine-based inhibitory motif (ITIM) and plays a key role in peripheral tolerance and autoimmune disease in mice. PD-1 is expressed mainly on activated T and B lymphocytes. Two novel B7 Family members have been identified as PD-1 ligands, PD-L1 (B7-H1) and PD-L2 (B7-DC). Evidence reported to date suggests overlapping functions for these ligands and their constitutive expression on some normal tissues and upregulation on activated antigen-presenting cells. It is reported that J43 inhibits the binding of mouse PD-L1-Ig and mouse PD-L2-Ig to PD-1/BHK transfected cells. When administrated in vivo, both intact and Fab of J43 are reported to enhance contact hypersensitivity and exacerbate acute GVHD similar to transfer of PD-1-deficient cells. Injection of J43 also exacerbates EAE and NOD diabetes as do specific antibodies to mouse PD-L1 and PD-L2.

Applications Reported: This J43 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This J43 antibody has been tested by flow cytometric analysis of stimulated mouse splenocytes. This may be used at less than or equal to 0.5 μ g per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Brilliant Ultra Violet[™] 805 (BUV805) is a tandem dye that emits at 797 nm and is intended for use on cytometers equipped with an ultraviolet (355 nm) laser. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright, Brilliant Violet™, Brilliant Ultra Violet™, or other polymer dye-conjugated antibodies in a

staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) or Brilliant Stain Buffer (Product # 00-4409-75) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer or Brilliant Stain Buffer for more information.

Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.


Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222) (100 μ L of cell sample + 100 μ L of IC Fixation Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Our internal testing suggests that Brilliant Ultra Violet[™] 805 (BUV805) is compatible with short-term methanol-based fixation, but should not be stored in buffers containing methanol for longer than one hour.

Excitation: 355 nm; Emission: 797 nm; Laser: Ultraviolet Laser.

BRILLIANT ULTRA VIOLET[™] is a trademark or registered trademark of Becton, Dickinson and Company or its affiliates, and is used under license. Powered by Sirigen[™]

Product Images For CD279 (PD-1) Monoclonal Antibody (J43), Brilliant Ultra Violet™ 805, eBioscience™

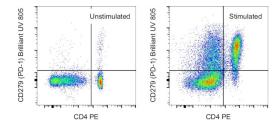

CD279 (PD-1) Antibody (368-9985-80) in Flow

BALB/c mouse splenocytes were unstimulated (left) or stimulated for 72 hours with CD3e and CD28 Monoclonal Antibodies, Functional Grade (Product # 16-0031-85) and (Product # 16-0281-85) (right). Cells were then stained with CD4 Monoclonal Antibody, PE (Product # 12-0042-82) and 0.25 µg of CD279 (PD-1) Monoclonal Antibody, Brilliant Ultra Violet 805 (BUV805). Viable cells in the lymphocyte gate were used for analysis, as determined by 7-AAD (Product # 00-6993-50).

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein its valid only when used by property trained individuals. Unless otherwise stated in the Documentation, this warranty are ror date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not textend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample furnices Does Developed textend to anyone other than the Buyer. Any model textend to anyone REFUND FOR THE WARRANTIES, EXPRESS OR INNERNICE DIRCLOING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULTS DURING THE WARRANTIES IN CENTRE VARCHANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES SOF MERCHANTABILITY, FITNESS FOR ANY PARTICUCTS DURING THE WARRANT FERION IS LIMITED TO REFLAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS DURING THE WARRANT FERION IS LIMITED TO REFLAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOL CENTION. THERE IS NO BUILSTATION TO FORCE MAJEURE, (III) USE OF ANY PARTICULTY, SOL CENTRO TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (III) USES, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product is intended for research only and is not to be used for