



CD4 Monoclonal Antibody (SK3 (SK-3)), Brilliant Violet™ 421, eBioscience™

Product Details	
Size	25 Tests
Species Reactivity	Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), Brilliant Violet™ 421, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	SK3 (SK-3)
Conjugate	Brilliant Violet™ 421
Excitation/Emission Max	406/423 nm
Form	Liquid
Concentration	5 μL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2925486

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 μL (0.125 μg)/test	-

Product Specific Information

Description: The SK3 monoclonal antibody reacts with human CD4, a 59-kDa cell surface receptor expressed by a majority of thymocytes, a subpopulation of mature T helper cells, and at low levels on monocytes. CD4 is a receptor for the human immunodeficiency virus (HIV). SK3 blocks HIV binding and mixed lymphocyte reaction. The SK3 and RPA-T4 monoclonal antibodies do not cross-block binding, suggesting recognition of distinct epitopes.

Applications Reported: This SK3 (SK-3) antibody has been reported for use in flow cytometric analysis.

Applications Tested: This SK3 (SK-3) antibody has been pre-diluted and tested by flow cytometric analysis of normal human peripheral blood cells. This may be used at 5 µL (0.125 µg) per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.


Brilliant Violet™ 421 (BV421) is a dye that emits at 423 nm and is intended for use on cytometers equipped with a violet (405 nm) laser. Please make sure that your instrument is capable of detecting this fluorochrome.


When using two or more Super Bright, Brilliant Violet™, Brilliant Ultra Violet™, or other polymer dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401-42) or Brilliant Stain Buffer™ (Product # 00-4409-75) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer or Brilliant Stain Buffer for more information.

Excitation: 407 nm; Emission: 423 nm; Laser: Violet Laser.

BRILLIANT VIOLETTM is a trademark or registered trademark of Becton, Dickinson and Company or its affiliates, and is used under license. Powered by SirigenTM.

Product Images For CD4 Monoclonal Antibody (SK3 (SK-3)), Brilliant Violet™ 421, eBioscience™

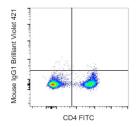

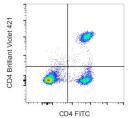

CD4 Antibody (404-0047-41) in Flow

Normal human peripheral blood cells were stained with CD3 Monoclonal Antibody, FITC (Product # 11-0037-42) and Mouse IgG1 kappa Isotype Control, Brilliant Violet 421 (Product # 404-4714-81) (left) or CD4 Monoclonal Antibody, Brilliant Violet 421 (right). Viable cells in the lymphocyte gate were used for analysis, as determined by 7-AAD (Product # 00-6993-50).

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation specifications and/use accompanying package inserts ("Documentation,"), to claim of substituted to proper the product is usually accompanying package inserts ("Documentation,"), to claim of substituted to proper the product is usually accompanying package inserts. In the product is usually accompanying package inserts ("Documentation,"), to claim of substituted to the product is usually accompanying package inserts. In the product is usually accompanying package inserts in the product is usually accompanying package inserts. In the product is usually accompanying package in the product is usually accompanying package in the product is used to accompanying the product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT.

BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERFIGIO IS LIMITED, ACCEDING TO BE PERIAR REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS, AT SELECTION, THERE IS NO GOLIGATION TO REPLAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS, AT SELECTION, THE RESON OF THE PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, or vivo or in vivo therapeutic uses,