



CD14 Monoclonal Antibody (61D3), Brilliant Violet™ 421, eBioscience™

Product Details		
Size	25 Tests	
Species Reactivity	Human	
Host/Isotype	Mouse / IgG1, kappa	
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), Brilliant Violet™ 421, eBioscience™	
Class	Monoclonal	
Туре	Antibody	
Clone	61D3	
Conjugate	Brilliant Violet™ 421	
Excitation/Emission Max	406/423 nm	
Form	Liquid	
Concentration	5 μL/Test	
Purification	Affinity chromatography	
Storage buffer	PBS, pH 7.2, with BSA	
Contains	0.09% sodium azide	
Storage conditions	4° C, store in dark, DO NOT FREEZE!	
RRID	AB_2921018	

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 μL (0.5 μg)/test	-

Product Specific Information

Description: The 61D3 monoclonal antibody reacts with human CD14, a 53-55 kDa GPI-linked glycoprotein. CD14 is expressed on monocytes, interfollicular macrophages and some dendritic cells. Complexes of LPS and LBP (LPS-Binding Protein) bind with high affinity to monocytes through the surface CD14.

Applications Reported: This 61D3 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This 61D3 antibody has been pre-diluted and tested by flow cytometric analysis of normal human peripheral blood cells. This may be used at 5 µL (0.5 µg) per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10^5 to 10^8 cells /test.


Brilliant Violet™ 421 (BV421) is a dye that emits at 423 nm and is intended for use on cytometers equipped with a violet (405 nm) laser. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright, Brilliant Violet™, Brilliant Ultra Violet™, or other polymer dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401-42) or Brilliant Stain Buffer™ (Product # 00-4409-75) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer or Brilliant Stain Buffer for more information.

Excitation: 407 nm; Emission: 423 nm; Laser: Violet Laser.

BRILLIANT VIOLET™ is a trademark or registered trademark of Becton, Dickinson and Company or its affiliates, and is used under license. Powered by Sirigen.™

Product Images For CD14 Monoclonal Antibody (61D3), Brilliant Violet™ 421, eBioscience™

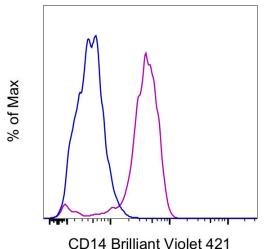

CD14 Antibody (404-0149-41) in Flow

Normal human peripheral blood cells were stained with Mouse IgG1 kappa Isotype Control, Brilliant Violet 421 (Product # 404-4714-81) (blue histogram) or CD14 Monoclonal Antibody, Brilliant Violet 421 (purple histogram). Cells in the monocyte gate were used for analysis.

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production of source and a specifications and/or accompanying package inserts ("Documentation"), No claim of suitability for use in applications regulated by EDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, in the applications regulated by EDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, and is a supplementation of the provided is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample is merely illustrative of the general pand quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICUS OR, NON INFRINGEMENT.

BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING HE WARRANTY PERFIOI DIS LIMITED. OR PEAPIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS, AT SELLER'S SOLE OPTION, THERE IS NO BUIGATION TO REPLIA, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS, AT SELLER'S SOLE OPTION, THERE IS NO. IN PROPERLY SOLE OPTION. THERE IS NO. IN PROPERLY SOLE OPTION. THE PRODUCTS AS THE RESULT OF (i) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (ii) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (iii) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (iv) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, includi