



CD127 Monoclonal Antibody (A7R34), PerCP-Cyanine5.5, eBioscience™

Product Details	
Size	25 μg
Species Reactivity	Mouse
Published Species	Mouse
Host/Isotype	Rat / IgG2a, kappa
Recommended Isotype Control	Rat IgG2a kappa Isotype Control (eBR2a), PerCP-Cyanine5.5, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	A7R34
Conjugate	PerCP-Cyanine5.5
Excitation/Emission Max	489/679 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_906212

Applications	Tested Dilution	Publications
Immunohistochemistry (Frozen) (IHC (F))	-	1 Publication
Immunocytochemistry (ICC/IF)	-	1 Publication
Flow Cytometry (Flow)	0.25 μg/test	33 Publications

Product Specific Information

Description: The A7R34 monoclonal antibody reacts with mouse CD127, the high affinity alpha subunit of the mouse IL-7 receptor. IL-7 receptor alpha chain is expressed by immature B cells in the bone marrow, double-negative (CD4-CD8-), single-positive (CD4+ and CD8+), but not double-positive (CD4+CD8+) thymocytes. In the periphery, mature T cells express CD127 at low level. A7R34 inhibits binding of IL-7 to its receptor and has been used in in vivo and in vitro studies to elucidate the role of IL-7 in T and B cell development and activation. Binding of A7R34 blocks the binding of SB/199, another antibody which recognizes mouse CD127.


Applications Reported: This A7R34 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This A7R34 antibody has been tested by flow cytometric analysis of mouse splenocytes. This can be used at less than or equal to $0.25~\mu g$ per test. A test is defined as the amount (μg) of antibody that will stain a cell sample in a final volume of $100~\mu L$. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Excitation: 488 nm; Emission: 695 nm; Laser: Blue Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For CD127 Monoclonal Antibody (A7R34), PerCP-Cyanine5.5, eBioscience™

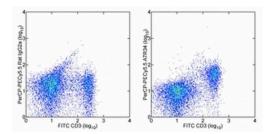

CD127 Antibody (45-1271-80) in Flow

Staining of C57BL/6 splenocytes with Anti-Mouse CD3e FITC (Product # 11-0031-82) and 0.5 µg of Rat IgG2a K Isotype Control PerCP-Cyanine5-5 (Product # 45-4321-80) (left) or 0.125 µg of Anti-Mouse CD127 PerCP-Cyanine5-5 (right). Cells in the lymphocyte gate were used for analysis.

View more figures on thermofisher.com

□ 35 References

Immunohistochemistry (Frozen) (1)

The American journal of pathology

Foxp3-expressing CD103+ regulatory T cells accumulate in dendritic cell aggregates of the colonic mucosa in murine transfer colitis.

Authors: Leithäuser F, Meinhardt-Krajina T, Fink K, Wotschke B, Möller P, Reimann J

Year 2006

Immunocytochemistry (1)

Journal of immunology (Baltimore, Md.: 1950)

LTbetaR signaling induces cytokine expression and up-regulates lymphangiogenic factors in lymph node anlagen.

Authors: Vondenhoff MF, Greuter M, Goverse G, Elewaut D, Dewint P, Ware CF, Hoorweg K, Kraal G, Mebius RE

Year 2009

Flow Cytometry (33)

ANJ, Libby P, Shi GP, Guo J

Advanced science (Weinheim, Baden-Wurttemberg, Germany)

Group 2 Innate Lymphoid Cells Protect Mice from Abdominal Aortic Aneurysm Formation via IL5 and Eosinophils.

"Published figure using CD127 monoclonal antibody (Product # 45-1271-82) in Flow Cytometry"

Authors: Zhang Y,Liu T,Deng Z,Fang W,Zhang X,Zhang S,Wang M,Luo S,Meng Z,Liu J,Sukhova GK,Li D,McKenzie

Year 2023

Cell reports

PI3K coordinates transcriptional, chromatin, and metabolic changes to promote effector CD8⁺ T cells at the expense of central memory.

"Published figure using CD127 monoclonal antibody (Product # 45-1271-82) in Flow Cytometry"

Authors: Cannons JL, Villarino AV, Kapnick SM, Preite S, Shih HY, Gomez-Rodriguez J, Kaul Z, Shibata H, Reilley JM, Huang B, Handon R, McBain IT, Gossa S, Wu T, Su HC, McGavern DB, O'Shea JJ, McGuire PJ, Uzel G, Schwartzberg PL

Year 2021

View more Flow references on thermofisher.com

More applications with references on thermofisher.com

For Research Use Only, Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production of specifications and/or accompanying package inserts ("Documentation"). No claim of studiability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is important to the product is used to a sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample in. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT.

BUYER'S EXCLUSIVE REMEDY FOR NON-CORNORNING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED. A PERIOD IS LIMITED. A PERION ES LIMITED. A PERIOD FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, or why or in vito therapeutic uses, or any type of consumption to human or animals.