CD29 (Integrin beta 1) Monoclonal Antibody (eBioHMb1-1 (HMb1-1)), PerCP-eFluor™ 710, eBioscience™

Size	100 µg
Species Reactivity	Mouse, Rat
Published Species	Mouse, Human
Host/Isotype	Armenian hamster / IgG
Recommended Isotype Control	Armenian Hamster IgG Isotype Control (eBio299Arm), PerCP-eFluor™ 710, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	eBioHMb1-1 (HMb1-1)
Conjugate	PerCP-eFluor™ 710
Excitation/Emission Max	482/708 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_10670099

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	0.5 μg/test	55 Publications

Product Specific Information

Description: The eBioHMb1-1 monoclonal antibody reacts with mouse and rat CD29 (integrin beta 1), a 110-120 kDa member of the beta integrin family expressed by leukocytes, endothelial, smooth muscle and epithelial cells. CD29 binds non-covalently with the alpha integrins CD49a-f to form the VLA-1 through VLA-6 complexes, as well as with CD51. These alpha-beta integrin heterodimers are capable of mediating a variety of cellular responses including adhesion, trafficking, proliferaton and differentiation. All integrins which include CD29 bind to extracellular matrix proteins including collagen, laminin, fibronectin and vitronectin, whereas some CD29-containing integrins can also interact with cellular receptors such as VCAM-1 and MadCAM-1.

Applications Reported: This eBioHMb1-1 (HMb1-1) antibody has been reported for use in flow cytometric analysis.

Applications Tested: This eBioHMb1-1 (HMb1-1) antibody has been tested by flow cytometric analysis of mouse bone marrow cells. This can be used at less than or equal to 0.5 μ g per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells /test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

PerCP-eFluor® 710 emits at 710 nm and is excited with the blue laser (488 nm); it can be used in place of PerCP-Cyanine5.5. We recommend using a 710/50 bandpass filter, however, the 695/40 bandpass filter is an acceptable alternative. Please make sure that your instrument is capable of detecting this fluorochrome.


Fixation: Samples can be stored in IC Fixation Buffer (cat. 00-8222) (100 µL cell sample + 100 µL IC Fixation Buffer) or 1-step

Fix/Lyse Solution (cat. 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency /compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 488 nm; Emission: 710 nm; Laser: Blue Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For CD29 (Integrin beta 1) Monoclonal Antibody (eBioHMb1-1 (HMb1-1)), PerCP-eFluor™ 710, eBioscience™

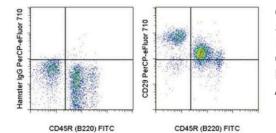

CD29 (Integrin beta 1) Antibody (46-0291-82) in Flow

Staining of C57BI/6 bone marrow cells with Anti-Human/Mouse CD45R (B220) FITC (Product # 11-0452-82) and 0.25 µg of Armenian Hamster IgG Isotype Control PerCP-eFluor® 710 (Product # 46-4888-82) (left) or 0.25 µg of Anti-Mouse/Rat CD29 (Integrin beta 1) PerCP-eFluor® 710 (right). Cells in the small /lymphocyte scatter population were used for analysis.

View more figures on thermofisher.com

□ 55 References

Flow Cytometry (55)

FEBS open bio A comparative study of mouse bone marrow mesenchymal stem cells isolated using three easy-to-perform approaches. "Published figure using CD29 (Integrin beta 1) monoclonal antibody (Product # 46-0291-82) in Flow Cytometry" Authors: Lu Y,Han Y,Zhou L,Shi G,Bai L,Wang K,Qin C	Year 2022
Heliyon Single-cell transcriptomics reveals variable trajectories of CSPCs in the progression of osteoarthritis.	Year 2022
"Published figure using CD29 (Integrin beta 1) monoclonal antibody (Product # 46-0291-82) in Flow Cytometry" Authors: Qi L,Wang J,Chen X,Ding Y,Ling B,Wang W,Xu J,Xue Z	

View more Flow references on thermofisher.com

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in offect at the time of saids as set forth in the Production documentation, specifications and/or sample functional companying published by trained individuals. Unless otherwise stated in the Documentation, this warrant is initiated to experif trained individuals. Unless otherwise stated in the Documentation, this warrant is used by roported or sample functional or sample functional provided herein is valid only when used by provide to the individuals. Unless otherwise stated in the Documentation, this warrant is used to example functional or sample functional provide to the period or sample functional or sample functional provide the entry is used or sample functional provide to the general type and quality of goods and does not represent that any Product will conform to such model or sample functional. PRESS OR INNERNES, EXPRESS OR INNERNES, EXPRESS OR INNERNE, EXPRESS OR

2