

CD73 Monoclonal Antibody (eBioTY/11.8 (TY/11.8)), PerCP-eFluor™ 710, eBioscience™

Product Details	
Size	100 µg
Species Reactivity	Mouse
Published Species	Mouse
Host/Isotype	Rat / IgG1
Recommended Isotype Control	Rat IgG1 kappa Isotype Control (eBRG1), PerCP-eFluor™ 710, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	eBioTY/11.8 (TY/11.8)
Conjugate	PerCP-eFluor™ 710
Excitation/Emission Max	482/708 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_10853356

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	0.125 µg/test	14 Publications

Product Specific Information

Description: eBioTY/11.8 recognizes the with CD73 a 69-kDa GPI-anchored cell-surface protein with ecto-5'-nucleotidase activity. Expression on myeloid cells (CD11b) is restricted to the bone marrow. In human CD73 can be induced to secrete a soluble form with IL-2 suggesting a role in mediating activation signals. Differences between human and mouse CD73 have been reported. BALB/c mice have more CD4+CD73+ than CD8+CD73+ while the reciprocal is documented in humans.

CD73 is expressed on a subset of lymphocytes and increases during lymphocyte maturation. Recently it has been found that memory CD4 T cells express and are similar to the uncommitted primed precursor helper cells (Thpp) that can differentiate into Th1 or Th2 cells. Furthermore CD73 has been found on regulatory T cells.

Applications Reported: This eBioTY/11.8 (TY/11.8) antibody has been reported for use in flow cytometric analysis.

Applications Tested: This eBioTY/11.8 (TY/11.8) antibody has been tested by flow cytometric analysis of mouse splenocytes. This can be used at less than or equal to 0.125 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.


PerCP-eFluor® 710 emits at 710 nm and is excited with the blue laser (488 nm); it can be used in place of PerCP-Cyanine5.5. We recommend using a 710/50 bandpass filter, however, the 695/40 bandpass filter is an acceptable alternative. Please make sure that your instrument is capable of detecting this fluorochrome.

Fixation: Samples can be stored in IC Fixation Buffer (cat. 00-8222) (100 µL cell sample + 100 µL IC Fixation Buffer) or 1-step Fix/Lyse Solution (cat. 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency /compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 488 nm; Emission: 710 nm; Laser: Blue Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For CD73 Monoclonal Antibody (eBioTY/11.8 (TY/11.8)), PerCP-eFluor™ 710, eBioscience™

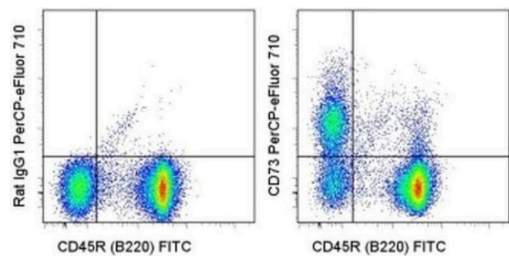

CD73 Antibody (46-0731-82) in Flow
Staining of C57Bl/6 splenocytes with Anti-Human/Mouse CD45R (B220) FITC (Product # 11-0452-82) and 0.06 µg of Rat IgG1 K Isotype Control PerCP-eFluor® 710 (Product # 46-4301-80) (left) or 0.06 µg of Anti-Mouse CD73 PerCP-eFluor® 710 (right). Total viable cells were used for analysis.

View more figures on thermofisher.com

14 References

Flow Cytometry (14)

Frontiers in endocrinology	Year 2022
Preventing type 1 diabetes in late-stage pre-diabetic NOD mice with insulin: A central role for alum as adjuvant.	Species Mouse
"46-0731-82 was used in Flow cytometry/Cell sorting to support the use of alum as adjuvant to optimise the efficacy of antigen-specific immunotherapy in future trials."	
Authors: Martens PJ,Ellis D,Buggeman Y,Viaene M,Laureys J,Teyton L,Mathieu C,Gysemans C	
BMC musculoskeletal disorders	Year 2021
Age-related trabecular bone loss is associated with a decline in serum Galectin-1 level.	
"Published figure using CD73 monoclonal antibody (Product # 46-0731-82) in Flow Cytometry"	
Authors: Xu W,Ni C,Wang Y,Zheng G,Zhang J,Xu Y	

View more Flow references on thermofisher.com

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.