

Phospho-STAT5 (Tyr694) Monoclonal Antibody (SRBCZX), PerCP-eFluor 710, eBioscience™

Product Details	
Size	100 Tests
Species Reactivity	Human, Mouse
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), PerCP-eFluor 710, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	SRBCZX
Conjugate	PerCP-eFluor™ 710
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.1% gelatin, 0.2% BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2573860

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 µL (0.06 µg)/test	5 Publications

Product Specific Information

Description: This SRBCZX monoclonal antibody recognizes signal transducer and activator of transcription 5 (STAT5) when phosphorylated on tyrosine 694. STAT proteins are activated by ligand binding to receptors, such as cytokine receptors, that associate with Janus kinase (JAK) family members. Following their phosphorylation by JAKs, STAT proteins translocate to the nucleus where they bind to DNA and regulate transcription of specific genes in a cell type- and cytokine-specific manner. In response to cytokines that signal through the common gamma chain such as IL-2, IL-7, and IL-15, STAT5 is phosphorylated on tyrosine 694 by JAK1 and JAK3. Cytokines such as IL-3, IL-5, and GM-CSF that signal via the common beta chain induce STAT5 phosphorylation on tyrosine 694 by JAK 2. Phosphorylation of STAT5 on tyrosine 694 is essential for STAT5 dimer formation, nuclear translocation, and DNA binding activity.

Specificity of this SRBCZX clone was determined by ELISA and flow cytometry.

Applications Reported: This SRBCZX antibody has been reported for use in intracellular staining followed by flow cytometric analysis.

Applications Tested: This SRBCZX antibody has been pre-titrated and tested by intracellular staining followed by flow cytometric analysis of stimulated normal human peripheral blood cells. This can be used at 5 µL (0.06 µg) per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test.

Staining Protocol: We recommend using Protocol C: Two-step protocol: Fixation/Methanol. Protocol A: Two-step protocol: intracellular (cytoplasmic) proteins and Protocol B: One-step protocol: intracellular (nuclear) proteins cannot be used. All Protocols can be found in the Flow Cytometry Protocols: "Staining Intracellular Antigens for Flow Cytometry Protocol" located in the Best Protocols Section under the Resources tab online.

PerCP-eFluor® 710 emits at 710 nm and is excited with the blue laser (488 nm); it can be used in place of PerCP-Cyanine5.5. We recommend using a 710/50 bandpass filter, however, the 695/40 bandpass filter is an acceptable alternative. Please make sure that your instrument is capable of detecting this fluorochrome.


Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.

Fixation: Samples can be stored in IC Fixation Buffer (cat. 00-8222) (100 µL of cell sample + 100 µL of IC Fixation Buffer) or 1-step Fix/Lyse Solution (cat. 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency /compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 488 nm; Emission: 710 nm; Laser: Blue Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For Phospho-STAT5 (Tyr694) Monoclonal Antibody (SRBCZX), PerCP-eFluor 710, eBioscience™

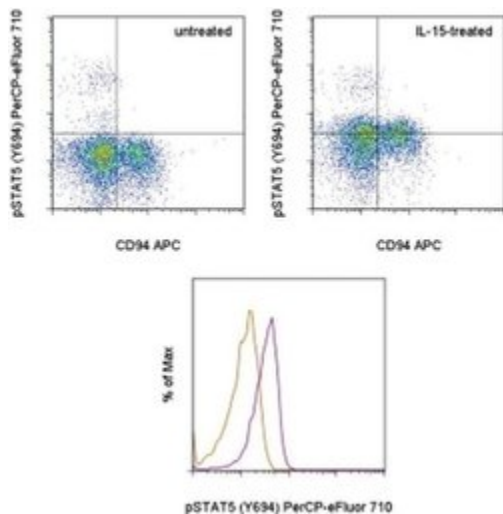

Phospho-STAT5 (Tyr694) Antibody (46-9010-42) in Flow

TOP: Intracellular staining of untreated (left) or 15-minute IL-15-treated (right) normal human peripheral blood cells with Anti-Canine/Human CD94 APC (Product # 17-5094-42) and Anti-Human/Mouse phospho-STAT5 (Y694) PerCP-eFluor® 710. Cells in the lymphocyte gate were used for analysis. BOTTOM: Intracellular staining of untreated (orange histogram) or 15-minute IL-15-treated (purple histogram) normal human peripheral blood cells with Anti-Human/Mouse phospho-STAT5 (Y694) PerCP-eFluor® 710. CD94+ cells in the lymphocyte gate were used for analysis. In both panels, cells were stained using the Fixation/Methanol protocol.

[View more figures on thermofisher.com](https://www.thermofisher.com)

5 References

Flow Cytometry (5)

eLife

Receptor-mediated dimerization of JAK2 FERM domains is required for JAK2 activation.

"Published figure using Phospho-STAT5 (Tyr694) monoclonal antibody (Product # 46-9010-42) in Flow Cytometry"

Authors: Ferrao RD,Wallweber HJ,Lupardus PJ

Species
Not Applicable

Dilution
Not Cited

Year
2018

Scientific reports

Erythropoietin enhances Kupffer cell number and activity in the challenged liver.

"Published figure using Phospho-STAT5 (Tyr694) monoclonal antibody (Product # 46-9010-42) in Flow Cytometry"

Authors: Gilboa D,Haim-Ohana Y,Deshet-Unger N,Ben-Califa N,Hiram-Bab S,Reuveni D,Zigmond E,Gassmann M, Gabet Y,Varol C,Neumann D

Species
Not Applicable

Dilution
Not Cited

Year
2017

[View more Flow references on thermofisher.com](#)

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.