



c-MAF Monoclonal Antibody (sym0F1), PerCP-eFluor™ 710, eBioscience™

Product Details		
Size	100 Tests	
Species Reactivity	Human, Mouse	
Published Species	Mouse	
Host/Isotype	Mouse / IgG2b, kappa	
Recommended Isotype Control	Mouse IgG2b kappa Isotype Control (eBMG2b), PerCP-eFluor™ 710, eBioscience™	
Class	Monoclonal	
Туре	Antibody	
Clone	sym0F1	
Conjugate	PerCP-eFluor™ 710	
Excitation/Emission Max	482/708 nm	
Immunogen	Full-length E. coli-produced c-Maf	
Form	Liquid	
Concentration	5 μL/Test	
Purification	Affinity chromatography	
Storage buffer	PBS, pH 7.2, with 0.2% BSA	
Contains	0.09% sodium azide	
Storage conditions	4° C, store in dark, DO NOT FREEZE!	
RRID	AB_2573908	

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 μL (0.0075 μg)/test	4 Publications

Product Specific Information

Description: The sym0F1 monoclonal antibody reacts with human and mouse c-Maf, a 42 kDa basic leucine zipper transcription factor shown to be involved in the neural, ocular and hematopoietic systems. In hematopoietic cells, it was first shown to be crucial for IL-4 expression in Th2 and was the first transcription factor believed to be Th subset-specific. More recent evidence shows that specific phospho-tyrosine residues lead to upregulation of IL-4. c-Maf has also been shown to be important to differentiation and function in both Th17 and Tfh cells. It drives expression of IL-21 in both cell types, while promoting expression of IL-23R in Th17 and CXCR5 in Tfh as well.

Applications Reported: This sym0F1 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This sym0F1 antibody has been pre-titrated and tested by intracellular staining followed by flow cytometric analysis of Th17-polarized mouse splenocytes using the Foxp3/Transcription Factor Buffer Set and protocol. Please see Best Protocols Section (Staining intracellular Antigens for Flow Cytometry) for staining protocol (refer to Protocol B: One-step protocol for intracellular (nuclear) proteins. This can be used at 5 μ L (0.0075 μ g) per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.

PerCP-eFluor® 710 emits at 710 nm and is excited with the blue laser (488 nm); it can be used in place of PerCP-Cyanine5.5.


We recommend using a 710/50 bandpass filter, however, the 695/40 bandpass filter is an acceptable alternative. Please make sure that your instrument is capable of detecting this fluorochrome.


Fixation: Samples can be stored in IC Fixation Buffer (cat. 00-8222) (100 μ L cell sample + 100 μ L IC Fixation Buffer) or 1-step Fix/Lyse Solution (cat. 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency /compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 488 nm; Emission: 710 nm; Laser: Blue Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For c-MAF Monoclonal Antibody (sym0F1), PerCP-eFluor™ 710, eBioscience™

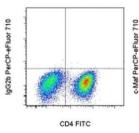

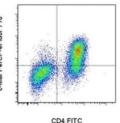

c-MAF Antibody (46-9855-42) in Flow

Staining of 3-day Th17-polarized mouse splenocytes with Anti-Mouse CD4 FITC (Product # 11-0042-82) and Mouse IgG2b K Isotype Control PerCP-eFluor® 710 (Product # 46-4732-82) (left) or Anti-Human/Mouse c-Maf PerCP-eFluor® 710 (right) using the Foxp3/Transcription Factor Buffer Set and protocol. Cells in the lymphocyte gate were used for analysis.

View more figures on thermofisher.com

□ 4 References

Flow Cytometry (4)

Journal of cell science

Year 2022

GM-CSF-activated human dendritic cells promote type 1 T follicular helper cell polarization in a CD40-dependent manner.

"Published figure using c-MAF monoclonal antibody (Product # 46-9855-42) in Flow Cytometry"

Authors: Korniotis S,Saichi M,Trichot C,Hoffmann C,Amblard E,Viguier A,Grondin S,Noel F,Mattoo H,Soumelis V

JCI insight

Year 2018

C5aR1 regulates T follicular helper differentiation and chronic graftversus-host disease bronchiolitis obliterans.

"Published figure using c-MAF monoclonal antibody (Product # 46-9855-42) in Flow Cytometry"

 $\label{thm:prop:matching} \mbox{Authors: Verghese DA,Chun N,Paz K,Fribourg M,Woodruff TM,Flynn R,Hu Y,Xiong H,Zhang W,Yi Z,Du J,Blazar BR, Heeger PS$

View more Flow references on thermofisher.com

For Research Use Only, Not for use in diagnostic procedures. Not for resale without express authorization, Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production aspectations and care accompanying package inserts (Documentation). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is important to the product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT.

BUYER'S EXCLUSIVE REMEDY FOR NON-CORNORMOR PRODUCTS DURING THE WARRANTY PERFIGIO IS LIMITED. ACCEDENT OF OR REFUND FOR THE NON-CONFORMING PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, or vivy or in vivo therapeutic uses, or any type of consumption to hourand or animal to home an animal to the buse of the support of the product of the major of consumption to hourand or animal to how therapeutic uses, or vivo or in vivo therapeutic uses, or vivo or in vivo therapeutic uses, or any type of consump