



# CD3 Monoclonal Antibody (UCHT1), APC-eFluor™ 780, eBioscience™

<b>Product Details</b>	
Size	100 Tests
Species Reactivity	Human
Published Species	Mouse, Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), APC-eFluor™ 780, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	UCHT1
Conjugate	APC-eFluor™ 780
Excitation/Emission Max	756/785 nm
Form	Liquid
Concentration	5 μL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.2% BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_1272042

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 μL (0.5 μg)/test	31 Publications

#### **Product Specific Information**

Description: The UCHT1 monoclonal antibody reacts with human CD3e, a 20 kDa subunit of the TCR complex. Along with the other CD3 subunits gamma and delta, the epsilon chain is required for proper assembly, trafficking and surface expression of the TCR complex. CD3 is expressed by thymocytes in a developmentally regulated manner and by all mature T cells. Crosslinking of TCR via immobilized UCHT1 initiates an intracellular biochemical pathway resulting in cellular activation and proliferation.

Applications Reported: This UCHT1 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This UCHT1 antibody has been pre-titrated and tested by flow cytometric analysis of normal human peripheral blood cells. This can be used at 5  $\mu$ L (0.5  $\mu$ g) per test. A test is defined as the amount ( $\mu$ g) of antibody that will stain a cell sample in a final volume of 100  $\mu$ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells /test.

APC-eFluor 780 emits at 780 nm and is excited with the Red laser (633 nm). Please make sure that your instrument is capable of detecting this fluorochome.

Light sensitivity: This tandem is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.


Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222) (100 µL cell sample + 100 µL IC Fixation Buffer) or 1-


step Fix/Lyse Solution (Product # 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 633-647 nm; Emission: 780 nm; Laser: Red Laser.

Filtration: 0.2 µm post-manufacturing filtered.


# Product Images For CD3 Monoclonal Antibody (UCHT1), APC-eFluor™ 780, eBioscience™

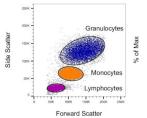


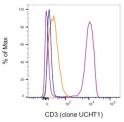

# CD3 Antibody (47-0038-42)

Staining of human peripheral blood cells. As expected based on known relative expression patterns, CD3 clone UCHT1 stains a subset of lymphocytes (T cells) and does not stain monocytes and granulocytes. Details: Normal human whole blood was surface stained with CD3 (clone UCHT1). After staining, red blood cells were lysed using 1-step Fix/Lyse Buffer. Cells in the lymphocyte (purple histogram), monocyte (orange histogram), or granulocyte (blue histogram) gates were used for analysis. {RE}

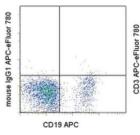

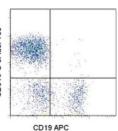

### CD3 Antibody (47-0038-42) in Flow

Staining of normal human peripheral blood cells with Anti-Human CD19 APC (Product # 17-0199-42) and Mouse IgG1 K Isotype Control APC-eFluor® 780 (Product # 47-4714-82) (left) or Anti-Human CD3 APC-eFluor® 780 (right). Cells in the lymphocyte gate were used for analysis.

#### **□ 31 References**

# Flow Cytometry (31)

#### Cell reports

# Engineered human cytokine/antibody fusion proteins expand regulatory T cells and confer autoimmune disease protection.

"47-0038-42 was used in Flow cytometry/Cell sorting to provide a roadmap for IC design and establish a Treg-biased immunotherapy that could be clinically translated for autoimmune disease treatment."

Authors: VanDyke D,Iglesias M,Tomala J,Young A,Smith J,Perry JA,Gebara E,Cross AR,Cheung LS,Dykema AG, Orcutt-Jahns BT,Henclová T,Golias J,Balolong J,Tomasovic LM,Funda D,Meyer AS,Pardoll DM,Hester J,Issa F,Hunter CA,Anderson MS,Bluestone JA,Raimondi G,Spangler JB

**Year** 2022

Species Mouse

Dilution 1:50

### **NPJ** vaccines

# Vaccine genetics of IGHV1-2 VRC01-class broadly neutralizing antibody precursor naïve human B cells.

"47-0038-42 was used in Flow cytometry/Cell sorting to find that IGHV1-2 alleles traditionally thought to be incompatible with VRC01-class responses are relatively common in various human populations and that germline variation within IGHV1-2 associates with gene usage frequencies in the naive BCR repertoire."

 $\label{prop:control} \mbox{Authors: Lee JH,Toy L,Kos JT,Safonova Y,Schief WR,Havenar-Daughton C,Watson CT,Crotty S} \\$ 

**Year** 2021

Species Human

Dilution 1:100

View more Flow references on thermofisher.com

#### More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization, Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications required by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT.

BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED. A REFUND FOR THE NON-CONFORMING PRODUCTS, SI TESLELER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS, SI TESLELER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS, SI TESLELER'S SOLE OPTION. THERE IS NO GREAT REPORT OF THE PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro di