

# TER-119 Monoclonal Antibody (TER-119), eFluor™ 450, eBioscience™

Product Details	
Size	25 µg
Species Reactivity	Mouse
Published Species	Fish, Mouse, Human
Host/Isotype	Rat / IgG2b, kappa
Recommended Isotype Control	Rat IgG2b kappa Isotype Control (eB149/10H5), eFluor™ 450, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	TER-119
Conjugate	eFluor™ 450
Excitation/Emission Max	405/445 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_1518809

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	6 Publications
Immunocytochemistry (ICC/IF)	-	1 Publication
Flow Cytometry (Flow)	0.5 µg/test	77 Publications

## Product Specific Information

**Description:** The TER-119 monoclonal antibody reacts with mouse erythroid cells from early proerythroblast to mature erythrocyte stages. The TER-119 antigen is present in yolk sac, fetal and newborn liver, but is not expressed by cells carrying BFU-E and CFU-E activities. Several erythroleukemia cell lines tested so far are negative for expression of TER-119 antigen even after dimethylsulfoxide stimulation. Biochemical and molecular analysis of the TER-119 antigen indicate that this molecule is associated with the surface glycoprotein A, but is not a typical glycoprotein.

**Applications Reported:** This TER-119 antibody has been reported for use in flow cytometric analysis.


**Applications Tested:** This TER-119 antibody has been tested by flow cytometric analysis of mouse bone marrow cells. This can be used at less than or equal to 0.5 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10<sup>5</sup> to 10<sup>8</sup> cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

eFluor® 450 is an alternative to Pacific Blue®. eFluor® 450 emits at 445 nm and is excited with the Violet laser (405 nm). Please make sure that your instrument is capable of detecting this fluorochrome.

Excitation: 405 nm; Emission: 445 nm; Laser: Violet Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For TER-119 Monoclonal Antibody (TER-119), eFluor™ 450, eBioscience™

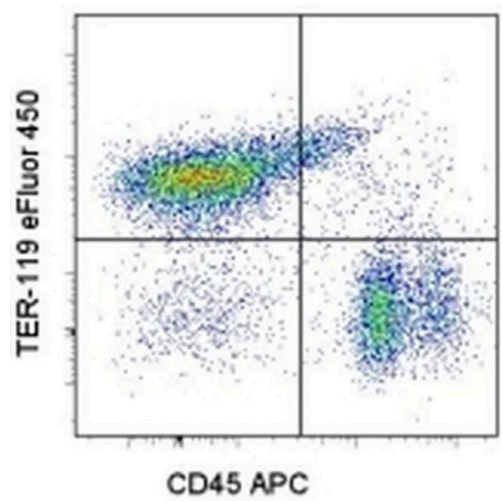

**TER-119 Antibody (48-5921-80) in Flow**  
Staining of C57BL/6 bone marrow cells with Anti-Mouse CD45 APC (Product # 17-0451-82) and 0.25 µg of Anti-Mouse TER-119 eFluor® 450. Total viable cells were used for analysis.

View more figures on [thermofisher.com](https://thermofisher.com)

## Immunohistochemistry (6)

<p><b>The Journal of clinical investigation</b></p> <p><b>Ribonuclease inhibitor 1 regulates erythropoiesis by controlling GATA1 translation.</b></p> <p>"Published figure using TER-119 monoclonal antibody (Product # 48-5921-82) in Immunocytochemistry"</p> <p>Authors: Chennupati V,Veiga DF,Maslowski KM,Andina N,Tardivel A,Yu EC,Stilinovic M,Simillion C,Duchosal MA,Quadroni M,Roberts I,Sankaran VG,MacDonald HR,Fasel N,Angelillo-Scherrer A,Schneider P,Hoang T,Allam R</p>	<p><b>Year</b> 2018</p>
<p><b>Development (Cambridge, England)</b></p> <p><b>Transient loss of venous integrity during developmental vascular remodeling leads to red blood cell extravasation and clearance by lymphatic vessels.</b></p> <p>"Published figure using TER-119 monoclonal antibody (Product # 48-5921-82) in Immunofluorescence"</p> <p>Authors: Zhang Y,Daubel N,Stritt S,Mäkinen T</p>	<p><b>Year</b> 2018</p>

[View more IHC references on thermofisher.com](#)

## Immunocytochemistry (1)

<p><b>The Journal of clinical investigation</b></p> <p><b>Ribonuclease inhibitor 1 regulates erythropoiesis by controlling GATA1 translation.</b></p> <p>"Published figure using TER-119 monoclonal antibody (Product # 48-5921-82) in Immunocytochemistry"</p> <p>Authors: Chennupati V,Veiga DF,Maslowski KM,Andina N,Tardivel A,Yu EC,Stilinovic M,Simillion C,Duchosal MA,Quadroni M,Roberts I,Sankaran VG,MacDonald HR,Fasel N,Angelillo-Scherrer A,Schneider P,Hoang T,Allam R</p>	<p><b>Year</b> 2018</p>
--	-----------------------------

## Flow Cytometry (77)

<p><b>PloS one</b></p> <p><b>Engineered red blood cells carrying PCSK9 inhibitors persistently lower LDL and prevent obesity.</b></p> <p>"Published figure using TER-119 monoclonal antibody (Product # 48-5921-82) in Flow Cytometry"</p> <p>Authors: Deshycka R,Sudaryo V,Huang NJ,Xie Y,Smeding LY,Choi MK,Ploegh HL,Lodish HF,Pishesha N</p>	<p><b>Year</b> 2021</p>
--	-----------------------------

[View more Flow references on thermofisher.com](#)

## More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.