IFN gamma Monoclonal Antibody (4S.B3), eFluor™ 450, eBioscience™

Size	100 Tests
Species Reactivity	Human
Published Species	Human, Mouse, Rhesus monkey
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), eFluor™ 450, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	4S.B3
Conjugate	eFluor™ 450
Excitation/Emission Max	405/445 nm
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.2% BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2043866

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	1 Publication
Flow Cytometry (Flow)	5 μL (0.25 μg)/test	66 Publications
ELISA (ELISA)	-	1 Publication

Product Specific Information

Description: The 4S.B3 monoclonal antibody reacts with interferon-gamma (IFN gamma). Human IFN gamma is a 17 kDa factor produced by activated T and NK cells and is an anti-viral and anti-parasitic cytokine. IFN gamma in synergy with other cytokines, such as TNF alpha, inhibits proliferation of normal and transformed cells. Immunomodulatory effects of IFN gamma are exerted on a wide range of cell types expressing the high affinity receptors for IFN gamma. Glycosylation of IFN gamma does not affect its biological activity.

Applications Reported: This 4S.B3 antibody has been reported for use in intracellular staining followed by flow cytometric analysis.


Applications Tested: This 4S.B3 antibody has been pre-titrated and tested by intracellular staining and flwo cytometric analysis of stimulated normal human peripheral blood cells. This can be used at 5 μ L (0.25 μ g) per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.

eFluor® 450 is an alternative to Pacific Blue®. eFluor® 450 emits at 445 nm and is excited with the Violet laser (405 nm). Please make sure that your instrument is capable of detecting this fluorochome.

Excitation: 405 nm; Emission: 445 nm; Laser: Violet Laser.


Filtration: 0.2 µm post-manufacturing filtered.

Product Images For IFN gamma Monoclonal Antibody (4S.B3), eFluor™ 450, eBioscience™

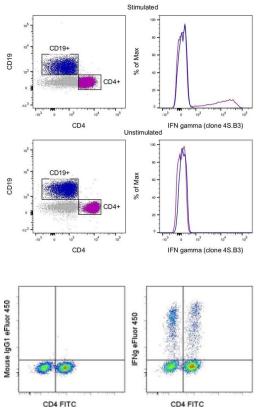

IFN gamma Antibody (48-7319-42)

Intracellular staining of stimulated human peripheral blood cells. As expected based on known expression patterns, IFN gamma clone 4S.B3 stains a subset of CD4+ T cells only after stimulation and does not stain CD19+ B cells regardless of stimulation. Details: Normal human peripheral blood cells were cultured in the presence of Protein Transport Inhibitors (500X) (Unstimulated, bottom row) or Cell Stimulation Cocktail (plus protein transport inhibitors, 500X) for 5 hours (Stimulated, top row). Cells were fixed and permeabilized with the IC Fixation & Permeabilization Buffer Set and protocol followed by intracellular staining with CD19 (clone SJ25C1), CD4 (clone RPA-T4), and IFN gamma (clone 4S.B3). Cells in the CD19+ (blue histogram) or CD4+ (purple histogram) gates were used for analysis. {TM}

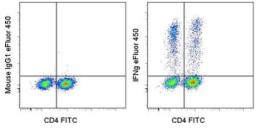

IFN gamma Antibody (48-7319-42) in Flow

Intracellular staining of 5-hour PMA/ionomycin stimulated normal human peripheral blood cells (in the presence of Brefeldin A) with Anti-Human CD4 FITC (Product # 11-0049-42) and Mouse IgG1 K Isotype Control eFluor® 450 (Product # 48-4714-82) (left) or Anti-Human IFN gamma eFluor® 450 (right). Cells in the lymphocyte gate were used for analysis.

View more figures on thermofisher.com

G 68 References

Immunohistochemistry (1)

Journal of translational medicine

Niraparib exhibits a synergistic anti-tumor effect with PD-L1 blockade by inducing an immune response in ovarian cancer.

"Published figure using IFN gamma monoclonal antibody (Product # 48-7319-42) in Flow Cytometry" Authors: Meng J,Peng J,Feng J,Maurer J,Li X,Li Y,Yao S,Chu R,Pan X,Li J,Zhang T,Liu L,Zhang Q,Yuan Z,Bu H,Song K,Kong B

Flow Cytometry (66)

Cell & bioscience GPR116 receptor regulates the antitumor function of NK cells via Gq /HIF1/NF-B signaling pathway as a potential immune checkpoint.	Year 2023
"Published figure using IFN gamma monoclonal antibody (Product # 48-7319-42) in Flow Cytometry"	
Authors: Guo D,Jin C,Gao Y,Lin H,Zhang L,Zhou Y,Yao J,Duan Y,Ren Y,Hui X,Ge Y,Yang R,Jiang W	
Cancer immunology research	Year
CD8+ T cell-Dependent Remodeling of the Tumor Microenvironment	2023

Overcomes Chemoresistance. "Published figure using IFN gamma monoclonal antibody (Product # 48-7319-42) in Flow Cytometry" Authors: Lao L,Zeng W,Huang P,Chen H,Jia Z,Wang P,Huang D,Chen J,Nie Y,Yang L,Wu W,Liu J

View more Flow references on thermofisher.com

Year 2021

Year 2021

ELISA (1)

Frontiers in immunology A Method of Assessment of Human Natural Killer Cell Phenotype and Function in Whole Blood.

"Published figure using IFN gamma monoclonal antibody (Product # 48-7319-42) in ELISA" Authors: Market M,Tennakoon G,Ng J,Scaffidi M,de Souza CT,Kennedy MA,Auer RC

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation in Suitability for use in applications regulated by FDA is made. The warranty provided here in is vaile only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warrant is imited to one year from date of shipment when the Product is subjected to normal, proper and individuals. Unless otherwise stated and use in applications are extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTES, EXPRESS OR INNERLO, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINCEMENT. BUYER'S EXCLUSIVE REMEUP FOR NON-CONFORMING PRODUCTS DURING THE WARRANT PERIOD IS LIMITED TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OPTION. THERE IS NO BUILGATION TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OPTION. THERE IS NO BUILGATION TO REPAR, REPLACE OR REFUND FOR THE RODUCTS AS THE RESULT OF (I) ACCLIEDENT, DISASTER OR EVENT OF FORCE MAJEURE, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. UNder softenwise extra or in the documentation accompanying the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vin vinterapurpose, including without limita