

CD20 Monoclonal Antibody (L26), eFluor 660, eBioscience

Product Details	
Size	100 µg
Species Reactivity	Cynomolgus Monkey, Human, Rhesus monkey
Published Species	Human
Host/Isotype	Mouse / IgG2a, kappa
Recommended Isotype Control	Mouse IgG2a kappa Isotype Control (eBM2a), eFluor 660, eBioscience
Class	Monoclonal
Type	Antibody
Clone	L26
Conjugate	eFluor® 660
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.1% gelatin
Contains	0.09% sodium azide
Storage Conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_11150959

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	3 Publications
Immunohistochemistry (Paraffin) (IHC (P))	20 µg/mL	-
Miscellaneous PubMed (Misc)	-	1 Publication
Immunocytochemistry (ICC/IF)	20 µg/mL	1 Publication

Product Specific Information

Description: The monoclonal antibody L26 recognizes human CD20, a two transmembrane spanning protein found on precursor and mature B lymphocytes. Expression is lost upon differentiation to plasma cells. The levels of CD20 are critical in evaluating malignant B lymphocytes including Non-Hodgkin's Lymphoma, and Acute B cell, Chronic Lymphocytic and Hairy-Cell Leukemias. CD20 plays a role in B cell differentiation and activation. The structure of the protein shows a small extracellular loop with the carboxy and amino termini in the cytoplasmic compartment. The L26 antibody recognizes an epitope in the cytoplasmic domain and therefore requires permeabilization of the tissue for staining applications. Unlike many anti-CD20 antibodies that recognize the extracellular domain, the epitope recognized by L26 will not be blocked by therapeutic antibodies to CD20, like Rituximab.

Applications Reported: This L26 antibody has been reported for use in immunohistochemical staining of formalin-fixed paraffin embedded tissue sections.


Applications Tested: This L26 antibody has been tested by immunohistochemistry on formalin-fixed paraffin embedded human tonsils using IHC Antigen Retrieval Solution - Low pH (cat. 00-4955).

eFluor® 660 is a replacement for Alexa Fluor® 647. eFluor® 660 emits at 659 nm and is excited with the red laser (633 nm). Please make sure that your instrument is capable of detecting this fluorochrome. This can be used at less than or equal to 20 µg /mL. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Excitation: 633-647 nm; Emission: 668 nm; Laser: Red Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For CD20 Monoclonal Antibody (L26), eFluor 660, eBioscience

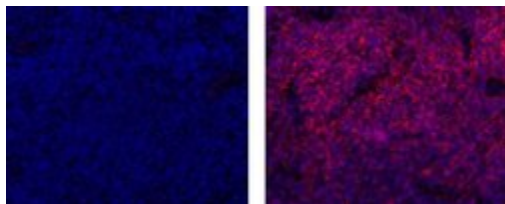

CD20 Antibody (50-0202-82) in IHC (P)

Immunohistochemistry on formalin-fixed paraffin embedded human tonsil, using 20 µg/mL of Mouse IgG2a K Isotype Control eFluor® 660 (left) or 20 µg/mL Anti-Human CD20 eFluor® 660 (right). Nuclei are stained with DAPI.

[View more figures on thermofisher.com](https://www.thermofisher.com)

5 References

Immunohistochemistry (3)

JCI insight

Heterogeneous fibroblasts underlie age-dependent tertiary lymphoid tissues in the kidney.

"Published figure using CD20 monoclonal antibody (Product # 50-0202-82) in Immunofluorescence"

Authors: Sato Y, Mii A, Hamazaki Y, Fujita H, Nakata H, Masuda K, Nishiyama S, Shibuya S, Haga H, Ogawa O, Shimizu A, Narumiya S, Kaisho T, Arita M, Yanagisawa M, Miyasaka M, Sharma K, Minato N, Kawamoto H, Yanagita M

Species
Not Applicable

Dilution
Not Cited

Year
2016

The Journal of clinical investigation

Accumulation of follicular CD8+ T cells in pathogenic SIV infection.

"Published figure using CD20 monoclonal antibody (Product # 50-0202-82) in Immunohistochemistry"

Authors: Ferrando-Martinez S, Moysi E, Pegu A, Andrews S, Nganou Makamdop K, Ambrozak D, McDermott AB, Palesch D, Paiardini M, Pavlakis GN, Brenchley JM, Douek D, Mascola JR, Petrovas C, Koup RA

Species
Not Applicable

Dilution
Not Cited

Year
2018

[View more IHC references on thermofisher.com](#)

Miscellaneous PubMed (1)

The Journal of clinical investigation

Accumulation of follicular CD8+ T cells in pathogenic SIV infection.

"Published figure using CD20 monoclonal antibody (Product # 50-0202-82) in Immunohistochemistry"

Authors: Ferrando-Martinez S, Moysi E, Pegu A, Andrews S, Nganou Makamdop K, Ambrozak D, McDermott AB, Palesch D, Paiardini M, Pavlakis GN, Brenchley JM, Douek D, Mascola JR, Petrovas C, Koup RA

Species
Not Applicable

Dilution
Not Cited

Year
2018

Immunocytochemistry (1)

The Journal of clinical investigation

Accumulation of follicular CD8+ T cells in pathogenic SIV infection.

"Published figure using CD20 monoclonal antibody (Product # 50-0202-82) in Immunohistochemistry"

Authors: Ferrando-Martinez S, Moysi E, Pegu A, Andrews S, Nganou Makamdop K, Ambrozak D, McDermott AB, Palesch D, Paiardini M, Pavlakis GN, Brenchley JM, Douek D, Mascola JR, Petrovas C, Koup RA

Species
Not Applicable

Dilution
Not Cited

Year
2018

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION FOR REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.