Musashi-2 Monoclonal Antibody (C1), eFluor™ 660, eBioscience™

Product Details

Size	100 µg
Species Reactivity	Human, Mouse
Host/Isotype	Mouse / IgG3, kappa
Class	Monoclonal
Туре	Antibody
Clone	C1
Conjugate	eFluor™ 660
Excitation/Emission Max	651/668 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2574345

Applications	Tested Dilution	Publications
Immunohistochemistry (Paraffin) (IHC (P))	5 μg/mL	-
Immunocytochemistry (ICC/IF)	Assay-Dependent	-

Product Specific Information

Description: The monoclonal antibody C1 recognizes human and mouse Musashi-2 (Msi2), a highly conserved RNA-binding protein. Musashi-2 contains two RNA-recognition domains and is expressed as one of two isoforms of 36 and 37kDa. Musashi 2 with its family member Musashi-1, are expressed at high levels in stem and tumor cells. During neural development, Musashi-2 is expressed in progenitor cells in the ventricular and subventricular zones, while postnatally, Musashi-2 expression is limited to multiple subsets of neurons and astrocytes in addition to adult stem cell populations. Functionally, Musashi-2 associates with Sox2 and is required for self-renewal of embryonic stem cells. Much like its role in maintaining stemness in neural progenitors, Musashi-2 is expressed in hematopoetic stem cells and plays a similar role in maintaining self-renewal. Musashi-2 has been found to be overexpressed in human myeloid leukemia and the levels of expression correlate with decreased survival. Musashi-2 may serve as a prognostic marker for acute myeloid leukemia (AML).

Applications Reported: This C1 antibody has been reported for use in immunohistochemical staining of formalin-fixed paraffin embedded tissue sections, microscopy, and immunocytochemistry.

Applications Tested: This C1 antibody has been tested by immunohistochemistry of formalin-fixed paraffin embedded tissue using high pH antigen retrieval and can be used at less than or equal to 5 µg/mL. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.


eFluor® 660 is a replacement for Alexa Fluor® 647. eFluor® 660 emits at 659 nm and is excited with the red laser (633 nm). Please make sure that your instrument is capable of detecting this fluorochome.

1

Excitation: 633-647 nm; Emission: 668 nm; Laser: Red Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For Musashi-2 Monoclonal Antibody (C1), eFluor™ 660, eBioscience™

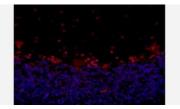

Musashi-2 Antibody (50-9677-82) in IHC (P)

Immunohistochemistry of formalin-fixed paraffin embedded human cerebellar tissue stained with 5 μ g/mL of Anti-Human/Mouse Musashi-2 eFluor® 660. Nuclei are stained with DAPI.

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty pervided herein is valid only when used by property trained individuals. Unless otherwise stated in the Documentation, this warranty is provided herein is valid only when used by property trained individuals. Unless otherwise stated in the Documentation, this warranty des not be trained to use other than the buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample furnished to Buyer is morely illustrative of the general type and quality of goods and does not superise that any Product will conform to such model or sample. NO OTHER WARRANTES, EXPRESS OR INNERDING WITHOUT LIMITATION, IMPLIED WARRANTES OF MERCHANTABILITY, FITHESS FOR ANY PARTICULA PURPOSE, OR NON INFRINCEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTED TO EXPRESS OR LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCTS) AS THE RESULT OF () ACCIDENT TO FORCE MAILEVILLE (I) MISUSE, FAULT OR NEGLEGATE OR REFUND FOR THE NON-CONFORMING PRODUCTS) AS THE RESULT OF () ACCIDENT TO FORCE MAILEVILLE (I) MISUSE, FAULT OR NEGLEGATE OR BY DUTER, (III) USE OF THE PRODUCTS) IN AMARKER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, are vinv or in