



CD86 (B7-2) Monoclonal Antibody (IT2.2), Alexa Fluor™ 488, eBioscience™

Product Details	
Size	100 Tests
Species Reactivity	Human
Published Species	Cynomolgus monkey, Human
Host/Isotype	Mouse / IgG2b, kappa
Recommended Isotype Control	Mouse IgG2b kappa Isotype Control (eBMG2b), Alexa Fluor™ 488, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	IT2.2
Conjugate	Alexa Fluor™ 488
Excitation/Emission Max	499/520 nm
Form	Liquid
Concentration	5 μL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.2% BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_10717088

Applications	Tested Dilution	Publications
Immunocytochemistry (ICC/IF)	-	1 Publication
Flow Cytometry (Flow)	5 μL (0.25 μg)/test	14 Publications
Neutralization (Neu)	-	1 Publication

Product Specific Information


Description: The IT2.2 monoclonal antibody reacts with human CD86, an ~80 kDa surface receptor also known as B7-2. CD86 and CD80 are members of the B7 family of costimulatory molecules. CD86 is expressed at low levels on B cells, macrophages, and dendritic cells and is upregulated on B cells through a variety of surface stimuli including the BCR complex, CD40 and some cytokine receptors. In addition to CD80 (B7-1), CD86 is a counter-receptor for the T cell surface molecules CD28 and CD152 (CTLA-4). The interaction of CD86 with its ligands plays a critical role in T-B crosstalk, T cell costimulation, autoantibody production and Th2-mediated Ig production. The kinetics of upregulation of CD86 upon stimulation supports its major contribution during the primary phase of an immune response.

Applications Reported: This IT2.2 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This IT2.2 antibody has been pre-titrated and tested by flow cytometric analysi sof normal human perioheral blood cells. This can be used at 5 μ L (0.25 μ g) per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test

Excitation: 488 nm; Emission: 519 nm; Laser: Blue Laser.

Product Images For CD86 (B7-2) Monoclonal Antibody (IT2.2), Alexa Fluor™ 488, eBioscience™

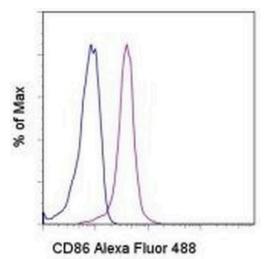

CD86 (B7-2) Antibody (53-0869-42) in Flow Staining of normal human peripheral blood cells with Mouse IgG2b K Isotype Control Alexa Fluor® 488 (Product # 53-4732-80) (blue histogram) or Anti-

Control Alexa Fluor® 488 (Product # 53-4732-80) (blue histogram) or Anti-Human CD86 (B7-2) Alexa Fluor® 488 (purple histogram). Cells in the monocyte gate were used for analysis.

View more figures on thermofisher.com

□ 16 References

Immunocytochemistry (1)

Biochemistry and biophysics reports

Pentraxin-3 regulates the inflammatory activity of macrophages.

"Published figure using CD86 (B7-2) monoclonal antibody (Product # 53-0869-42) in Immunocytochemistry" Authors: Shiraki A,Kotooka N,Komoda H,Hirase T,Oyama JI,Node K

Year 2016

Species Human

Flow Cytometry (14)

PeerJ

SARS-CoV-2 Delta (B.1.617.2) variant replicates and induces syncytia formation in human induced pluripotent stem cell-derived macrophages.

"Published figure using CD86 (B7-2) monoclonal antibody (Product # 53-0869-42) in Flow Cytometry"

 $\label{thm:model} Authors: Thaweer attanasinp T, Wanitchang A, Saenboonrueng J, Srisutthisamphan K, Wanasen N, Sungsuwan S, Jongkaewwattana A, Chailangkarn T$

Year 2023

Cells

TcpC Inhibits M1 but Promotes M2 Macrophage Polarization via Regulation of the MAPK/NF-B and Akt/STAT6 Pathways in Urinary Tract Infection.

"Published figure using CD86 (B7-2) monoclonal antibody (Product # 53-0869-42) in Flow Cytometry" Authors: Fang J,Ou Q,Wu B,Li S,Wu M,Qiu J,Cen N,Hu K,Che Y,Ma Y,Pan J

Year 2022

View more Flow references on thermofisher.com

Neutralization (1)

PLoS pathogens

Pomalidomide restores immune recognition of primary effusion lymphoma through upregulation of ICAM-1 and B7-2.

"Published figure using CD86 (B7-2) monoclonal antibody (Product # 53-0869-42) in Neutralization" Authors: Shrestha P,Davis DA,Jaeger HK,Stream A,Aisabor AI,Yarchoan R

Year 2021

More applications with references on thermofisher.com

For Research Use Only, Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications and/or accompanying package inserts ("Documentation"), No claim of suitability for use in application, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRRINGEMENT.

BUYER'S EXCLUSIVE EXEMBED FOR NON-CORNORNOR PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRRINGEMENT.

BUYER'S EXCLUSIVE EXEMBED FOR NON-CORNORN PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED, AND ADDITIONAL PURPOSE, OR NON INFRRINGEMENT.

BUYER'S EXCLUSIVE EXEMBED FOR NON-CORNORN PRODUCTS SOLE OF THE NON-CONFORMING PRODUCTS, AT SELECTION THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS, AT SELECTION THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS, AT SELECTION THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS AT THE NON-CONFORMING PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IM