Performand

CD20 Monoclonal Antibody (2H7), Alexa Fluor™ 700, eBioscience™

r roddor Botano	
Size	100 Tests
Species Reactivity	Human
Published Species	Human, Rhesus monkey
Host/Isotype	Mouse / IgG2b, kappa
Recommended Isotype Control	Mouse IgG2b kappa Isotype Control (eBMG2b), Alexa Fluor™ 700, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	2H7
Conjugate	Alexa Fluor™ 700
Excitation/Emission Max	696/719 nm
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.2% BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_1963574

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	4 Publications
Immunocytochemistry (ICC/IF)	-	1 Publication
Flow Cytometry (Flow)	5 µL (0.06 µg)/test	11 Publications

Product Specific Information

Description: The 2H7 monoclonal antibody reacts with human CD20, a 33-36 kDa transmembrane protein. CD20 is expressed by developing B cells as well as mature B cells but not plasma cells. CD20 has been detected at low levels on a small subset of mature T cells. It is suggested that CD20 plays a role in B-cell activation.

Applications Reported: This 2H7 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This 2H7 antibody has been pre-titrated and tested by flow cytometric analysis of normal human peripheral blood cells. This can be used at 5 μ L (0.06 μ g) per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.


Alexa Fluor® 700 emits at 723 nm and can be excited with the red laser (633 nm). Most instruments will require a 685 LP mirror and 710/20 filter. Please make sure that your instrument is capable of detecting this fluorochrome.

1

Excitation: 633-647 nm; Emission: 723 nm; Laser: Red Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For CD20 Monoclonal Antibody (2H7), Alexa Fluor™ 700, eBioscience™

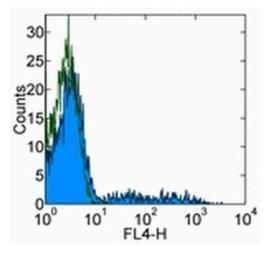

CD20 Antibody (56-0209-42) in Flow

Staining of normal human peripheral blood cells with staining buffer (autofluorescence) (open histogram) or Anti-Human CD20 Alexa Fluor® 700 (filled histogram). Cells in the lymphocyte gate were used for analysis.

View more figures on thermofisher.com

□16 References

Immunohistochemistry (4)

Frontiers in immunology	Year
L-33 Alarmin and Its Active Proinflammatory Fragments Are Released in Small Intestine in Celiac Disease.	2021
Published figure using CD20 monoclonal antibody (Product # 56-0209-42) in Immunohistochemistry"	
Authors: Perez F,Ruera CN,Miculan E,Carasi P,Dubois-Camacho K,Garbi L,Guzman L,Hermoso MA,Chirdo FG	
Stem cell research & therapy	Year
Relative safety of various spermatogenic stem cell purification methods	2019
for application in spermatogenic stem cell transplantation.	
Published figure using CD20 monoclonal antibody (Product # 56-0209-42) in Immunocytochemistry"	
Authors: Tian J,Ma K,Pei CB,Zhang SH,Li X,Zhou Y,Yan B,Wang HY,Ma LH	
View more	e IHC references on thermofisher.
nmunocytochemistry (1)	
Stem cell research & therapy	Year
Relative safety of various spermatogenic stem cell purification methods	2019
for application in spermatogenic stem cell transplantation.	
Published figure using CD20 monoclonal antibody (Product # 56-0209-42) in Immunocytochemistry"	
Authors: Tian J,Ma K,Pei CB,Zhang SH,Li X,Zhou Y,Yan B,Wang HY,Ma LH	

Flow Cytometry (11)

Aging and disease	Year 2022
A Soluble NK-CAR Mediates the Specific Cytotoxicity of NK Cells toward	
the Target CD20 ⁺ Lymphoma Cells.	
"Published figure using CD20 monoclonal antibody (Product # 56-0209-42) in Flow Cytometry"	
Authors: Liu R,Luo Q,Luo W,Wan L,Zhu Q,Yin X,Lu X,Song Z,Wei L,Xiang Z,Zou Y	

View more Flow references on thermofisher.com

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product till conform to such model or sample. No OTHER MARANTIES, EXPRESS OR IMPLED, ARRANTIES OF MERCHANTABILITY, THTRESS FOR ANY PARTICLAR PURPOSE, OR NON INFRINCEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS UNITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCTS) AS THE RESULT OF (i) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (ii) MISUSF, FAULT OR NEGLIGENCE OF OR BY BUYER, (iii) USE OF THE RODUCTS) A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (iv) IMPROPER STORAGE AND HANDLING OF THE RODUCTS. UNlish other warranty accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to numar or animals.