Performance guarenteed

Ki-67 Monoclonal Antibody (SolA15), Alexa Fluor™ 700, eBioscience™

Product Details	
Size	100 µg
Species Reactivity	Dog, Cynomolgus monkey, Human, Mouse, Non-human primate, Rat
Published Species	Rat, Mouse, Human
Host/Isotype	Rat / IgG2a, kappa
Recommended Isotype Control	Rat IgG2a kappa Isotype Control (eBR2a), Alexa Fluor™ 700, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	SolA15
Conjugate	Alexa Fluor™ 700
Excitation/Emission Max	696/719 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2637480

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	70 Publications
Immunohistochemistry (Paraffin) (IHC (P))	-	3 Publications
Immunohistochemistry (PFA fixed) (IHC (PFA))	-	1 Publication
Immunohistochemistry (Frozen) (IHC (F))	-	1 Publication
Immunocytochemistry (ICC/IF)	-	46 Publications
Flow Cytometry (Flow)	0.25 µg/test	59 Publications
Functional Assay (FN)	-	1 Publication

Product Specific Information

Description: The monoclonal antibody SolA15 recognizes mouse and rat Ki-67, a 300 kDa nuclear protein. Ki-67 is present during all active phases of the cell cycle (G1, S, G2, and mitosis), but is absent from resting cells (G0). Ki-67 is detected within the nucleus during interphase but redistributes to the chromosomes during mitosis. Ki-67 is used as a marker for determining the growth fraction of a given population of cells. In studies of tumor cells, the "Ki-67 labeling index" refers to the number of Ki-67 positive cells within the population and this is used to predict outcome of particular cancer types. Ki-67 has been shown to interact with the DNA-bound protein chromobox protein homolog 3 (CBX3) (heterochromatin).

The SoIA15 antibody also recognizes human, non-human primate and canine Ki-67.

Applications Reported: This SolA15 antibody has been reported for use in intracellular staining followed by flow cytometric

1

analysis.


Applications Tested: This SolA15 antibody has been tested by intracellular staining and flow cytometric analysis of stimulated mouse splenocytes using the Foxp3/Transcription Factor Staining Buffer Set (cat. 00-5523) and protocol. Please refer to Best Protocols: Protocol B: One step protocol for (nuclear) intracellular proteins located under the Resources Tab online. This antibody can be used at less than or equal to 0.25 μ g per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells /test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Alexa Fluor® 700 emits at 723 nm and can be excited with the red laser (633 nm). Most instruments will require a 685 LP mirror and 710/20 filter. Please make sure that your instrument is capable of detecting this fluorochrome.

Excitation: 633-647 nm; Emission: 723 nm; Laser: Red Laser.

Filtration: 0.2 µm post-manufacturing filtered.

Product Images For Ki-67 Monoclonal Antibody (SoIA15), Alexa Fluor™ 700, eBioscience™

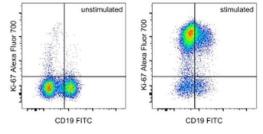

Ki-67 Antibody (56-5698-82) in Flow

Mouse splenocytes were unstimulated (left) or stimulated for 2 days with Anti-Mouse CD3e Functional Grade Purified (Product # 16-0031-82) (right). Cells were surface stained with Anti-Mouse CD19 FITC (Product # 11-0193-82) then fixed and permeabilized with the Foxp3/Trasncription Factor Staining Buffer Set (Product # 00-5523-00) and intracellularly stained with 0.125 µg of Anti-Mouse /Rat Ki-67 Alexa Fluor® 700. Total viable cells, as determined by Fixable Viability Dye eFluor® 450 (Product # 65-0863-14), were used for analysis.

View more figures on thermofisher.com

□ 181 References

Immunohistochemistry (70)

Scientific reports Application of an instructive hydrogel accelerates re-epithelialization of Renografted human skin wounds.	Year 2022
Published figure using Ki-67 monoclonal antibody (Product # 56-5698-82) in Immunohistochemistry" Authors: Sparks HD,Mandla S,Vizely K,Rosin N,Radisic M,Biernaskie J	
ournal of orthopaedic surgery and research	Year
GTSE1 is possibly involved in the DNA damage repair and cisplatin esistance in osteosarcoma.	2021
Published figure using Ki-67 monoclonal antibody (Product # 56-5698-82) in Immunohistochemistry" .uthors: Xie C,Xiang W,Shen H,Shen J	

Immunohistochemistry (Paraffin) (3)


Frontiers in oncology	Year 2020
Radiation Induced Metabolic Alterations Associate With Tumor	
Aggressiveness and Poor Outcome in Glioblastoma.	
Published figure using Ki-67 monoclonal antibody (Product # 56-5698-82) in Immunocytochemistry"	
Authors: Gupta K,Vuckovic I,Zhang S,Xiong Y,Carlson BL,Jacobs J,Olson I,Petterson XM,Macura SI,Sarkaria J,Burns IC	
Science	Year
Misershieta, Course of Duringer, Course and Mound Hashing, and Museum	2020

Microbiota-Sourced Purines Support Wound Healing and Mucous Barrier Function.

"Published figure using Ki-67 monoclonal antibody (Product # 56-5698-82) in Immunohistochemistry (Paraffin)" Authors: Lee JS,Wang RX,Goldberg MS,Clifford GP,Kao DJ,Colgan SP

View more IHC (P) references on thermofisher.com

More applications with references on thermofisher.com

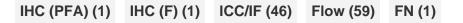

For Research Use Only, Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production occumentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications represent material to apprete or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, this warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation or sample control or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. No OTHER WARRANTES, EXPRESS OR INVELD, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES OF MERCHANTABILITY, FITNESS FOR ANY PARTICLAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEMEY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTES NO E DEFINIC, PERIOD IS LIMITED TO REPAIR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLED OFTICUTS ONE TO FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLICENCE OF OR BY BUYER'S SOLE OFTICUES AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLICENCE OF OR BY BUYER', (III) USE OF THE PRODUCTS. UNING STATE WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. UNless otherwise expressly stated on the Product on the documentation accompanying the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorizad commercial uses, in vitro diagnostic uses, or vivo frow o