

CD19 Monoclonal Antibody (SJ25C1), Alexa Fluor™ 532, eBioscience™

Product Details	
Size	25 Tests
Species Reactivity	Human
Published Species	Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), Alexa Fluor™ 532, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	SJ25C1
Conjugate	Alexa Fluor™ 532
Excitation/Emission Max	534/553 nm
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.2% BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_11219076

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	1 Publication
Flow Cytometry (Flow)	5 µL (0.125 µg)/test	11 Publications

Product Specific Information

Description: The SJ25C1 monoclonal antibody reacts with human CD19, a 95 kDa transmembrane glycoprotein. CD19 is expressed by B cells during all stages of development excluding the terminally differentiated plasma cells. Follicular dendritic cells also express this molecule. CD19, along with CD21, CD81, Leu13, and MHC class II, form a multimolecular complex that associates with the BCR. Signaling through CD19 induces tyrosine phosphorylation, calcium flux, and proliferation of B cells. The SJ25C1 antibody and the HIB19 monoclonal antibody recognize overlapping epitopes.

Applications Reported: This SJ25C1 antibody has been reported for use in flow cytometric analysis.


Applications Tested: This SJ25C1 antibody has been pre-titrated and tested by flow cytometric analysis of normal human peripheral blood cells. This can be used at 5 µL (0.125 µg) per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test.

Alexa Fluor® 532 is excited with the Green laser (532 nm) and emits at 561 nm. This cannot be used with the Yellow-Green laser (561 nm). We recommend using a 560/14 band pass filter. Please make sure that your instrument is capable of detecting

this fluorochrome.

Excitation: 532 nm; Emission: 561 nm; Laser: Green Laser

Product Images For CD19 Monoclonal Antibody (SJ25C1), Alexa Fluor™ 532, eBioscience™

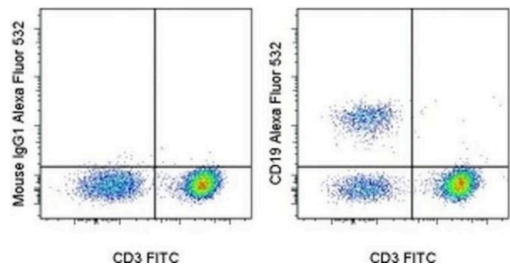

CD19 Antibody (58-0198-41) in Flow
Staining of normal human peripheral blood cells with Anti-Human CD3 FITC (Product # 11-0038-42) and Mouse IgG1 K Isotype Control Alexa Fluor® 532 (Product # 58-4714-80) (left) or Anti-Human CD19 Alexa Fluor® 532 (right). Cells in the lymphocyte gate were used for analysis.

[View more figures on thermofisher.com](#)

12 References

Immunohistochemistry (1)

<p>Nature communications</p> <p>Innate-like self-reactive B cells infiltrate human renal allografts during transplant rejection.</p> <p>"Published figure using CD19 monoclonal antibody (Product # 58-0198-41) in Flow Cytometry"</p> <p>Authors: Asano Y,Daccache J,Jain D,Ko K,Kinloch A,Veselits M,Wolfgeher D,Chang A,Josephson M,Cunningham P, Tambur A,Khan AA,Pillai S,Chong AS,Clark MR</p>	<p>Year</p> <p>2021</p>
---	-------------------------

Flow Cytometry (11)

<p>Biomolecules</p> <p>Tyrosine Kinase Inhibitors Target B Lymphocytes.</p> <p>"Published figure using CD19 monoclonal antibody (Product # 58-0198-41) in Flow Cytometry"</p> <p>Authors: Upfold NLE,Petakh P,Kamyshnyi A,Oksenyich V</p>	<p>Year</p> <p>2023</p>
--	-------------------------

<p>Indian journal of pharmacology</p> <p>Evaluation of peripheral lymphocyte subsets in acne vulgaris patients before and after systemic isotretinoin treatment.</p> <p>"Published figure using CD19 monoclonal antibody (Product # 58-0198-41) in Flow Cytometry"</p> <p>Authors: Gunaydin SD,Tezcan I</p>	<p>Year</p> <p>2022</p>
--	-------------------------

[View more Flow references on thermofisher.com](#)

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON-INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.