



## CD161 Monoclonal Antibody (HP-3G10), PE-eFluor™ 610, eBioscience™

<b>Product Details</b>	
Size	100 Tests
Species Reactivity	Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), PE-eFluor™ 610, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	HP-3G10
Conjugate	PE-eFluor™ 610
Excitation/Emission Max	565/606 nm
Form	Liquid
Concentration	5 μL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with 0.2% BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2574585

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 μL (0.5 μg)/test	1 Publication

#### **Product Specific Information**

Description: This HP-3G10 monoclonal antibody reacts with CD161 (also known as NKR-P1A), a member of the C-type lectin superfamily. The human homologue of NK1.1 in C57BL/6 mice, CD161 is expressed on natural killer cells and is upregulated in response to IL-12. CD161 can also be detected on various T cell subsets, including memory/effector CD4+ T cells, CD8+ T cells, gamma delta TCR T cells, and a subset of CD3+ thymocytes. Finally, CD161 expression has been demonstrated on human Th17 CD4+ T cells. The function of this receptor is unclear although studies suggest a possible stimulatory role. Nevertheless, Lectin-like transcript-1 (LLT1), which is also known as osteoclast inhibitory lectin, has been identified as the ligand for CD161.

Applications Reported: This HP-3G10 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This HP-3G10 antibody has been pre-titrated and tested by flow cytometric analysis of normal human peripheral blood cells. This can be used at 5 µL (0.5 µg) per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10^5 to 10^8 cells /test.

PE-eFluor® 610 can be excited with laser lines from 488-561 nm and emits at 607 nm. We recommend using a 610/20 band pass filter (equivalent to PE-Texas Red®). Please make sure that your instrument is capable of detecting this fluorochome.


Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.


Fixation: Samples can be stored in IC Fixation Buffer (cat. 00-8222) (100  $\mu$ L of cell sample + 100  $\mu$ L of IC Fixation Buffer) or 1-step Fix/Lyse Solution (cat. 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency /compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 488-561 nm; Emission: 607 nm; Laser: Blue Laser, Green Laser, Yellow-Green Laser.

Filtration: 0.2 µm post-manufacturing filtered.

#### Product Images For CD161 Monoclonal Antibody (HP-3G10), PE-eFluor™ 610, eBioscience™

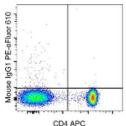

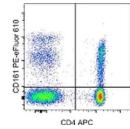

#### CD161 Antibody (61-1619-42) in Flow

Staining of normal human peripheral blood cells with Anti-Human CD4 APC (Product # 17-0049) and Mouse IgG1 K Isotype Control PE-eFluor® 610 (Product # 61-4714) (left) or Anti-Human CD161 PE-eFluor® 610 (right). Cells in the lymphocyte gate were used for analysis.

#### **□ 1 Reference**

### Flow Cytometry (1)

Frontiers in immunology

**Year** 2021

# Biological Therapy in Inflammatory Bowel Disease Patients Partly Restores Intestinal Innate Lymphoid Cell Subtype Equilibrium.

"Published figure using CD161 monoclonal antibody (Product # 61-1619-42) in Flow Cytometry"

Authors: Creyns B,Jacobs I,Verstockt B,Cremer J,Ballet V,Vandecasteele R,Vanuytsel T,Ferrante M,Vermeire S,Van Assche G,Ceuppens JL,Breynaert C

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production obcumentation specifications and/or accompanying package inserts ("Documentation"), No claim of suitability for use in applications regulated by EPA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, in the product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT.

BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCTS) AT SELLER'S SOLD OR SULTEN, DISASTER OR EVENT OF FORCE MAJEURE, (ii) INSUESE, FAULT OR NEGLIGENCE OF OR BY BUYER, (iii) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (iv) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, or vivy for er) of consumption to hourand or animal to hour and an animal research only and is not to be used for any other purpose, including without limitation to human or animal diagnostic uses, or vivy for er) of consumption to hourand or animal research.