. . . . .

## CD5 Monoclonal Antibody (L17F12), Super Bright<sup>™</sup> 436, eBioscience<sup>™</sup>

Product Details	
Size	25 Tests
Species Reactivity	Human
Host/Isotype	Mouse / IgG2a, kappa
Recommended Isotype Control	Mouse IgG2a kappa Isotype Control (eBM2a), Super Bright™ 436, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	L17F12
Conjugate	Super Bright™ 436
Excitation/Emission Max	413/431 nm
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2762493

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 μl (0.25 μg)/test	-

## **Product Specific Information**

Description: The L17F12 monoclonal antibody reacts with human CD5, a 67 kDa protein expressed by a majority of thymocytes and mature T cells and a subset of B cells. Signaling through the CD5 molecule activates T cells and binding of CD5 to its ligand on B cells, CD72, and plays an important role in T-B interaction and proliferation.

Applications Reported: This L17F12 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This L17F12 antibody has been pre-diluted and tested by flow cytometric analysis of normal human peripheral blood cells. This may be used at 5  $\mu$ L (0.25  $\mu$ g) per test. A test is defined as the amount ( $\mu$ g) of antibody that will stain a cell sample in a final volume of 100  $\mu$ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.


Super Bright 436 can be excited with the violet laser line (405 nm) and emits at 436 nm. We recommend using a 450/50 bandpass filter, or equivalent. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer for more information.

Excitation: 405 nm; Emission: 436 nm; Laser: Violet Laser

Super Bright Polymer Dyes are sold under license from Becton, Dickinson and Company.

## Product Images For CD5 Monoclonal Antibody (L17F12), Super Bright<sup>™</sup> 436, eBioscience<sup>™</sup>

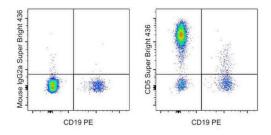

## CD5 Antibody (62-0058-41) in Flow

Normal human peripheral blood cells were stained with CD19 Monoclonal Antibody, PE (Product # 12-0199-42) and Mouse IgG2a kappa Isotype Control, Super Bright 436 (Product # 62-4724-82) (left) or CD5 Monoclonal Antibody, Super Bright 436 (right). Cells in the lymphocyte gate were used for analysis.

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regress nuthorization. Product sare warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, this warranty go no extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. No OTHER WARRANTES, EXPRESS OR INNELD, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES OF MERCHANTABILITY, FITNESS FOR ANY PARTICURS DE NON INFRINCEMENT. BUYER'S EXCLUSIVE EMEMORY FOR NON-CONFORMING PRODUCTS JUNISH THE WARRANT PERIOD IS LIMITED TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS OURING THE WARRANT FESOLE OFTION. THERE IS NO BULGATION TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OFTION. THERE IS NO BULGATION TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OFTION. THERE IS NO BULGATION TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OFTION. THERE IS NO BULGATION TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OFTION. THERE IS NO BULGATION TO REPAR, REPLACE OR REFUND FOR THE NON-CONFORMING PRODUCTS IN SALES ON ON UNPROPER STORAGE AND HANDILING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, arvivo or in vivo therapseutic uses, orany ty