CD27 Monoclonal Antibody (O323), Super Bright[™] 436, eBioscience[™]

Product Details	
Size	100 Tests
Species Reactivity	Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), Super Bright [™] 436, eBioscience [™]
Class	Monoclonal
Туре	Antibody
Clone	O323
Conjugate	Super Bright [™] 436
Excitation/Emission Max	413/431 nm
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2637433

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 μL (0.25 μg)/test	11 Publications

Product Specific Information

Description: The O323 monoclonal antibody reacts with human CD27, a lymphocyte-specific member of the TNFR superfamily. CD27 is expressed by a subset of thymocytes and virtually all mature T cells and is upregulated upon T-cell stimulation. CD27 binds to CD70, and through this interaction, plays an important role in T cell-B cell interaction.

Applications Reported: This O323 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This O323 antibody has been pre-titrated and tested by flow cytometric analysis of normal human peripheral blood cells. This can be used at 5 μ L (0.25 μ g) per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.


Super Bright 436 can be excited with the violet laser line (405 nm) and emits at 436 nm. We recommend using a 450/50 bandpass filter, or equivalent. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer for more information.

Excitation: 405 nm; Emission: 436 nm; Laser: Violet Laser

Super Bright Polymer Dyes are sold under license from Becton, Dickinson and Company.

Product Images For CD27 Monoclonal Antibody (O323), Super Bright[™] 436, eBioscience[™]

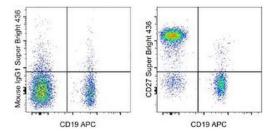

CD27 Antibody (62-0279-42) in Flow

Staining of normal human peripheral blood cells with Anti-Human CD19 APC (Product # 17-0199-42) and Mouse IgG1 K Isotype Control Super Bright 436 (Product # 62-4714-82) (left) or Anti-Human CD27 Super Bright 436 (right). Cells in the lymphocyte gate were used for analysis.

View more figures on thermofisher.com

□ 11 References

Flow Cytometry (11)

Cell discovery More than one antibody of individual B cells revealed by single-cell immune profiling.	Year 2023
"Published figure using CD27 monoclonal antibody (Product # 62-0279-42) in Flow Cytometry"	
Authors: Shi Z,Zhang Q,Yan H,Yang Y,Wang P,Zhang Y,Deng Z,Yu M,Zhou W,Wang Q,Yang X,Mo X,Zhang C,Huang J,Dai H,Sun B,Zhao Y,Zhang L,Yang YG,Qiu X	
Frontiers in immunology	Year
Oxygen level is a critical regulator of human B cell differentiation and	2023
IgG class switch recombination.	
"Published figure using CD27 monoclonal antibody (Product # 62-0279-42) in Flow Cytometry"	

View more Flow references on thermofisher.com

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuation. Unless otherwise stated in the Documentation, this warranty gere from tase of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample furnished to Buyer is provided therein is valid only when used by PDA is made. NO THER WARRANTES, EXPRESS OR INNELED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEMPT FOR NON-CONFORMING PRODUCTS DURING THE WARRANTES IN CEPTIAL REPLACE OR REFLIND. FOR THE NON-CONFORMING PRODUCTS SOLEN OF THE WARRANTE IN CEPTIAL REPLACE OR REFLIND. FOR THE NON-CONFORMING PRODUCTS SOLEN OF THE ANY PARTICULAR ULRAPURPER, CEPTIAC OR REFLIND. FOR THE NON-CONFORMING PRODUCTS SOLEN OF THE ANY OF IN CONFORMING PRODUCTS DURING TO TO REPRESS OF ANY EXCLUSIVE REMEMPT. BUYER'S SOLE OFTION. THERE IS NO BILLICATION TO REPRESS OF ANY EXCLUSIVE REMEMPT. THE NON-CONFORMING PRODUCTS SOLEN OFTION. THERE IS NO BILLICATION TO REPRESS OF ANY EXCLUSIVE REMEMPT. THE NON-CONFORMING PRODUCTS AS THE RESULT OF (I) ACCLIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER', (III) USE OF THE PRODUCTS. UNless otherwise expressly stated on the Product or in the documentation accompanying the Product, the