

CD69 Monoclonal Antibody (H1.2F3), Super Bright™ 436, eBioscience™

Product Details	
Size	100 µg
Species Reactivity	Mouse
Host/Isotype	Armenian hamster / IgG
Recommended Isotype Control	Armenian Hamster IgG Isotype Control (eBio299Arm), Super Bright™ 436, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	H1.2F3
Conjugate	Super Bright™ 436
Excitation/Emission Max	413/431 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2688109

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	0.5 µg/test	27 Publications

Product Specific Information

Description: The H1.2F3 monoclonal antibody reacts with mouse CD69, also known as very early activation antigen (VEA). CD69 is approximately 35 kDa and is expressed on the surface as a disulfide-linked dimer. While a small subset of lymphocytes in the thymus, spleen and lymph nodes express this antigen, activation of both T and B cells rapidly upregulates the surface expression of CD69, suggesting a role for CD69 in lymphocyte development and activation.

Applications Reported: This H1.2F3 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This H1.2F3 antibody has been tested by flow cytometric analysis of stimulated mouse splenocytes. This may be used at less than or equal to 0.5 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.


Super Bright 436 can be excited with the violet laser line (405 nm) and emits at 436 nm. We recommend using a 450/50 bandpass filter, or equivalent. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer for more information.

Excitation: 405 nm; Emission: 436 nm; Laser: Violet Laser

Super Bright Polymer Dyes are sold under license from Becton, Dickinson and Company.

Product Images For CD69 Monoclonal Antibody (H1.2F3), Super Bright™ 436, eBioscience™

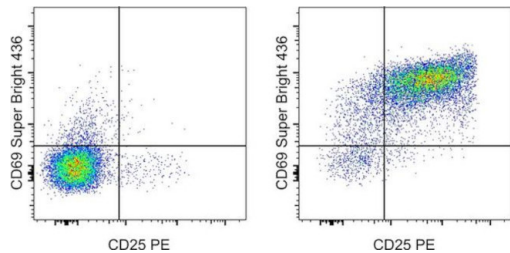

CD69 Antibody (62-0691-82) in Flow
BALB/c mouse splenocytes were unstimulated (left) or stimulated overnight with Con A (Product # 00-4978-03) (right). Cells were then stained with CD25 Monoclonal Antibody, PE (Product # 12-0251-82) and 0.25 µg of CD69 Monoclonal Antibody, Super Bright 436. Cells in the lymphocyte gate were used for analysis.

View more figures on thermofisher.com

27 References

Flow Cytometry (27)

Theranostics	Year 2023
HSF1 promotes CD69⁺ Treg differentiation to inhibit colitis progression.	
"Published figure using CD69 monoclonal antibody (Product # 62-0691-82) in Flow Cytometry"	
Authors: Yu L,Zhou B,Zhu Y,Li L,Zhong Y,Zhu L,Wang H,Chen H,Xu J,Guo T,Feng L,Wang X,Cai Z,Wang J,Jin H	
PloS one	Year 2023
Novel Bruton's tyrosine kinase inhibitor TAS5315 suppresses the progression of inflammation and joint destruction in rodent collagen-induced arthritis.	
"Published figure using CD69 monoclonal antibody (Product # 62-0691-82) in Flow Cytometry"	
Authors: Akasaka D,Iguchi S,Kaneko R,Yoshiga Y,Kajiwara D,Nakachi Y,Noma N,Tanaka K,Shimizu A,Hosoi F	

View more Flow references on thermofisher.com

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.