



CD127 Monoclonal Antibody (A7R34), Super Bright™ 436, eBioscience™

Product Details	
Size	100 μg
Species Reactivity	Mouse
Host/Isotype	Rat / IgG2a, kappa
Recommended Isotype Control	Rat IgG2a kappa Isotype Control (eBR2a), Super Bright™ 436, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	A7R34
Conjugate	Super Bright [™] 436
Excitation/Emission Max	413/431 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2744786

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	1.0 µg/test	7 Publications

Product Specific Information

Description: The A7R34 monoclonal antibody reacts with mouse CD127, the high affinity alpha subunit of the mouse IL-7 receptor. IL-7 receptor alpha chain is expressed by immature B cells in the bone marrow, double-negative (CD4-CD8-), singlepositive (CD4+ and CD8+), but not double-positive (CD4+CD8+) thymocytes. In the periphery, mature T cells express CD127 at low level. A7R34 inhibits binding of IL-7 to its receptor and has been used in in vivo and in vitro studies to elucidate the role of IL-7 in T and B cell development and activation. Binding of A7R34 blocks the binding of SB/199, another antibody which recognizes mouse CD127.

Applications Reported: This A7R34 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This A7R34 antibody has been tested by flow cytometric analysis of mouse splenocytes. This can be used at less than or equal to 1 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 μL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.


Super Bright 436 can be excited with the violet laser line (405 nm) and emits at 436 nm. We recommend using a 450/50 bandpass filter, or equivalent. Please make sure that your instrument is capable of detecting this fluorochrome.


When using two or more Super Bright dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer for more information.

Excitation: 405 nm; Emission: 436 nm; Laser: Violet Laser

Super Bright Polymer Dyes are sold under license from Becton, Dickinson and Company.

Product Images For CD127 Monoclonal Antibody (A7R34), Super Bright™ 436, eBioscience™

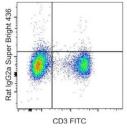

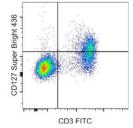

CD127 Antibody (62-1271-82) in Flow

Swiss Webster mouse splenocytes were stained with CD3e Monoclonal Antibody, FITC (Product # 11-0031-82) and 1.0 µg of Rat IgG2a kappa Isotype Control, Super Bright 436 (Product # 62-4321-82) (left) or 1.0 µg of CD127 Monoclonal Antibody, Super Bright 436 (right). Cells in the lymphocyte gate were used for analysis.

View more figures on thermofisher.com

□ 7 References

Flow Cytometry (7)

Advanced science (Weinheim, Baden-Wurttemberg, Germany)

Year 2023

Group 2 Innate Lymphoid Cells Protect Mice from Abdominal Aortic Aneurysm Formation via IL5 and Eosinophils.

"Published figure using CD127 monoclonal antibody (Product # 62-1271-82) in Flow Cytometry"

Authors: Zhang Y,Liu T,Deng Z,Fang W,Zhang X,Zhang S,Wang M,Luo S,Meng Z,Liu J,Sukhova GK,Li D,McKenzie ANJ,Libby P,Shi GP,Guo J

Cell reports

Year 2021

PI3K coordinates transcriptional, chromatin, and metabolic changes to promote effector CD8⁺ T cells at the expense of central memory.

"Published figure using CD127 monoclonal antibody (Product # 62-1271-82) in Flow Cytometry"

Authors: Cannons JL,Villarino AV,Kapnick SM,Preite S,Shih HY,Gomez-Rodriguez J,Kaul Z,Shibata H,Reilley JM, Huang B,Handon R,McBain IT,Gossa S,Wu T,Su HC,McGavern DB,O'Shea JJ,McGuire PJ,Uzel G,Schwartzberg PL

View more Flow references on thermofisher.com

More applications with references on thermofisher.com

For Research Use Only, Not for use in diagnostic procedures. Not for resale without express authorization, Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"), No claim of sultability for use in applications, and/or accompanying package inserts ("Documentation"), No claim of sultability for use in application, this warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated to in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTIBAILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT.

BUYER'S EXCLUSIVE EXEMBLY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITACE OR REFUND FOR THE NON-CONFORMING PRODUCTS SOLE OF THE PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, or vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.