

# HLA-DR Monoclonal Antibody (LN3), Super Bright™ 436, eBioscience™

Product Details	
Size	100 Tests
Species Reactivity	Human
Host/Isotype	Mouse / IgG2b, kappa
Recommended Isotype Control	Mouse IgG2a kappa Isotype Control (eBM2a), Super Bright™ 436, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	LN3
Conjugate	Super Bright™ 436
Excitation/Emission Max	413/431 nm
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4°C, store in dark, DO NOT FREEZE!
RRID	AB_2744822

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	3 Publications
Immunohistochemistry (Frozen) (IHC (F))	-	1 Publication
Immunocytochemistry (ICC/IF)	-	1 Publication
Flow Cytometry (Flow)	5 µL (0.125 µg)/test	10 Publications

## Product Specific Information

**Description:** The LN3 mAb reacts with the human major histocompatibility complex (MHC) class II, HLA-DR. HLA-DR is expressed on the surface of human antigen presenting cells (APC) including B cells, monocytes, macrophages, DCs, and activated T cells. HLA-DR is a heterodimeric transmembrane protein composed of alpha and beta subunits and plays an important role in the presentation of peptides to CD4+ T lymphocytes.

**Applications Reported:** This LN3 antibody has been reported for use in flow cytometric analysis.

**Applications Tested:** This LN3 antibody has been pre-titrated and tested by flow cytometric analysis of normal human peripheral blood cells. This can be used at 5 µL (0.125 µg) per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10<sup>5</sup> to 10<sup>8</sup> cells/test.

Super Bright 436 can be excited with the violet laser line (405 nm) and emits at 436 nm. We recommend using a 450/50 bandpass filter, or equivalent. Please make sure that your instrument is capable of detecting this fluorochrome.


When using two or more Super Bright dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright

Complete Staining Buffer (Product # SB-4401) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer for more information.

Excitation: 405 nm; Emission: 436 nm; Laser: Violet Laser

Super Bright Polymer Dyes are sold under license from Becton, Dickinson and Company.

## Product Images For HLA-DR Monoclonal Antibody (LN3), Super Bright™ 436, eBioscience™

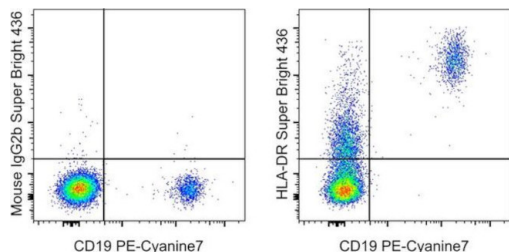

### HLA-DR Antibody (62-9956-42) in Flow

Normal human peripheral blood cells were stained with CD19 Monoclonal Antibody, PE-Cyanine7 (Product # 25-0199-82) and Mouse IgG2b kappa Isotype Control, Super Bright 436 (Product # 62-4732-82) (left) or HLA-DR Monoclonal Antibody, Super Bright 436 (right). Cells in the lymphocyte gate were used for analysis.

[View more figures on thermofisher.com](https://www.thermofisher.com)

## 15 References

### Immunohistochemistry (3)

Antioxidants (Basel, Switzerland)

#### Dimethyl Fumarate, an Approved Multiple Sclerosis Treatment, Reduces Brain Oxidative Stress in SIV-Infected Rhesus Macaques: Potential Therapeutic Repurposing for HIV Neuroprotection.

"Published figure using HLA-DR monoclonal antibody (Product # 62-9956-42) in Immunohistochemistry"

Authors: Garcia-Mesa Y, Xu HN, Vance P, Gruenewald AL, Garza R, Midkiff C, Alvarez-Hernandez X, Irwin DJ, Gill AJ, Kolson DL

Year  
2021

Nature communications

#### Gene expression and functional deficits underlie TREM2-knockout microglia responses in human models of Alzheimer's disease.

"Published figure using HLA-DR monoclonal antibody (Product # 62-9956-42) in Immunohistochemistry"

Authors: McQuade A, Kang YJ, Hasselmann J, Jairaman A, Sotelo A, Coburn M, Shabestari SK, Chadarevian JP, Fote G, Tu CH, Danhash E, Silva J, Martinez E, Cotman C, Prieto GA, Thompson LM, Steffan JS, Smith I, Davtyan H, Cahalan M, Cho H, Blurton-Jones M

Year  
2020

[View more IHC references on thermofisher.com](#)

### Immunohistochemistry (Frozen) (1)

Acta neuropathologica communications

#### Brain macrophages acquire distinct transcriptomes in multiple sclerosis lesions and normal appearing white matter.

"Published figure using HLA-DR monoclonal antibody (Product # 62-9956-42) in Immunohistochemistry (Frozen)"

Authors: Miedema A, Gerrits E, Brouwer N, Jiang Q, Kracht L, Meijer M, Nutma E, Peferoen-Baert R, Pijnacker ATE, Wesseling EM, Wijering MHC, Gabius HJ, Amor S, Eggen BJL, Kooistra SM

Year  
2022

### Immunocytochemistry (1)

Oncoimmunology

#### The role of irreversible electroporation in promoting M1 macrophage polarization via regulating the HMGB1-RAGE-MAPK axis in pancreatic cancer.

"Published figure using HLA-DR monoclonal antibody (Product # 62-9956-42) in Immunocytochemistry"

Authors: He C, Sun S, Zhang Y, Xie F, Li S

Year  
2021

### More applications with references on thermofisher.com

## Flow (10)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON-INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.