CD27 Monoclonal Antibody (O323), Super Bright[™] 600, eBioscience[™]

Product Details	
Size	100 Tests
Species Reactivity	Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), Super Bright [™] 600, eBioscience [™]
Class	Monoclonal
Туре	Antibody
Clone	O323
Conjugate	Super Bright™ 600
Excitation/Emission Max	414/601 nm
Form	Liquid
Concentration	5 µL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2662723

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 μL (0.25 μg)/test	4 Publications

Product Specific Information

Description: The O323 monoclonal antibody reacts with human CD27, a lymphocyte-specific member of the TNFR superfamily. CD27 is expressed by a subset of thymocytes and virtually all mature T cells and is upregulated upon T-cell stimulation. CD27 binds to CD70, and through this interaction, plays an important role in T cell-B cell interaction.

Applications Reported: This O323 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This O323 antibody has been pre-titrated and tested by flow cytometric analysis of normal human peripheral blood cells. This can be used at 5 μ L (0.25 μ g) per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.

Super Bright 600 is a tandem dye that can be excited with the violet laser line (405 nm) and emits at 600 nm. We recommend using a 610/20 bandpass filter. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer for more information.

Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Protect this vial and stained samples from light.

Fixation: Samples can be stored in IC Fixation Buffer (cat. 00-8222) (100 µL of cell sample + 100 µL of IC Fixation Buffer) or 1-


1

step Fix/Lyse Solution (cat. 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency /compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 405 nm; Emission: 600 nm; Laser: Violet Laser

Super Bright Polymer Dyes are sold under license from Becton, Dickinson and Company.

Product Images For CD27 Monoclonal Antibody (O323), Super Bright[™] 600, eBioscience[™]

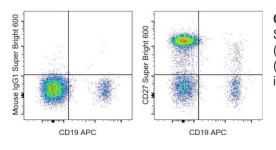

CD27 Antibody (63-0279-42) in Flow

Staining of normal human peripheral blood cells with Anti-Human CD19 APC (Product # 17-0199-42) and Mouse IgG1 K Isotype Control Super Bright 600 (Product # 63-4714-82) (left) or Anti-Human CD27 Super Bright 600 (right). Cells in the lymphocyte gate were used for analysis.

View more figures on thermofisher.com

Cell discovery	Year
More than one antibody of individual B cells revealed by single-cell	2023
mmune profiling.	
Published figure using CD27 monoclonal antibody (Product # 63-0279-42) in Flow Cytometry"	
Authors: Shi Z,Zhang Q,Yan H,Yang Y,Wang P,Zhang Y,Deng Z,Yu M,Zhou W,Wang Q,Yang X,Mo X,Zhang C,Huang I,Dai H,Sun B,Zhao Y,Zhang L,Yang YG,Qiu X	
Frontiers in immunology	Year
Oxygen level is a critical regulator of human B cell differentiation and	2023
gG class switch recombination.	
Published figure using CD27 monoclonal antibody (Product # 63-0279-42) in Flow Cytometry"	

Authors: Koers J,Marsman C,Steuten J,Tol S,Derksen NIL,Ten Brinke A,van Ham SM,Rispens T

View more Flow references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production or supplications and/or accompanying participations in effect at the time of sale, as set forth in the Decumentation, this warrant is specifications and/or accompanying part from date of shipment when the Product subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is bujected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample furnished. Buyer is Buyer AP model Studies and Buyer AP model Studies and Buyer AP model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not sepresent that any Product will conform to such model or sample furnished. Product subjected to normal, proper and intended usage. This warranty does not textend to anyone REFUND FOR THE WARRANTIES FOR ANY PARTICULTS DURING THE WARRANTES, EXPRESS OR INHELD, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTES SOF MERCHANTABILITY, FITNESS FOR ANY PARTICUTS DURING THE WARRANT PERFIDIO IS LIMITED TO REPAR, REPLACE CREATERY FOR THE NON-CONFORMING PRODUCTS SOLE OPTION. THERE IS NO BUIGATION TO REPAR, REPLACE OR REFUND FOR THE RONDLY of () ACCULTENT, DISASTER OR EVENT OF FORCE MAJEURE, (III) USE OF THE PRODUCTS NA MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (V) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. UN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (V) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. UN a mannee, set vico or in vico therapeutic uses, or any type of consumption by or application to

2