

CD45R (B220) Monoclonal Antibody (RA3-6B2), Super Bright™ 600, eBioscience™

Product Details	
Size	100 µg
Species Reactivity	Human, Mouse
Published Species	Fish, Human, Mouse
Host/Isotype	Rat / IgG2a, kappa
Recommended Isotype Control	Rat IgG2a kappa Isotype Control (eBR2a), Super Bright™ 600, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	RA3-6B2
Conjugate	Super Bright™ 600
Excitation/Emission Max	414/601 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Protein G
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2637457

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	30 Publications
Immunohistochemistry (Paraffin) (IHC (P))	-	2 Publications
Immunohistochemistry (Frozen) (IHC (F))	-	1 Publication
Immunocytochemistry (ICC/IF)	-	5 Publications
Flow Cytometry (Flow)	0.5 µg/test	62 Publications
Miscellaneous PubMed (Misc)	-	1 Publication

Product Specific Information

Description: The RA3-6B2 monoclonal antibody reacts with exon A-restricted isoform of mouse CD45, a 220 kDa surface molecule. CD45R/B220 epitope is mainly expressed by the B cell lineage from early Pro-B to mature B cells. However, some activated T cells, lymphokine activated killer cells (LAK), NK cell progenitors in the bone marrow, and T cells of the lpr/lpr mutant mouse also express this antigen.

Applications Reported: This RA3-6B2 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This RA3-6B2 antibody has been tested by flow cytometric analysis of mouse splenocytes. This can be used at less than or equal to 0.5 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Super Bright 600 is a tandem dye that can be excited with the violet laser line (405 nm) and emits at 600 nm. We recommend using a 610/20 bandpass filter. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer for more information.

Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Protect this vial and stained samples from light.

Fixation: Samples can be stored in IC Fixation Buffer (cat. 00-8222) (100 µL of cell sample + 100 µL of IC Fixation Buffer) or 1-step Fix/Lyse Solution (cat. 00-5333) for up to 3 days in the dark at 2-8°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorochrome performance after fixation can be made, but clone specific performance should be determined empirically.


Excitation: 405 nm; Emission: 600 nm; Laser: Violet Laser

Super Bright Polymer Dyes are sold under license from Becton, Dickinson and Company.

Product Images For CD45R (B220) Monoclonal Antibody (RA3-6B2), Super Bright™ 600, eBioscience™


CD45R (B220) Antibody (63-0452-82)

Staining of mouse splenocytes and bone marrow cells. As expected based on known relative expression patterns, CD45R (B220) clone RA3-6B2 stains subsets in both the splenocytes gate and the bone marrow lymphoid gate but not the bone marrow myeloid gate. Details: Balb/c bone marrow cells (left) and splenocytes (middle) were surface stained with CD45R (B220, clone RA3-6B2) followed by staining with 7-AAD. Viable bone marrow cells in the lymphoid (blue histogram) and myeloid (purple histogram) gates and viable splenocytes (orange histogram) were used for analysis. {RE}

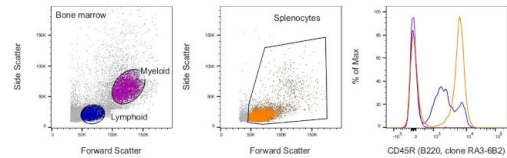

CD45R (B220) Antibody (63-0452-82) in Flow

Staining of C57Bl/6 splenocytes with Anti-Mouse CD3 APC (Product # 17-0032-82) and 0.25 µg of Rat IgG2a K Isotype Control Super Bright 600 (Product # 63-4321-82) (left) or 0.25 µg of Anti-Human/Mouse CD45R (B220) Super Bright 600 (right). Total viable cells were used for analysis.

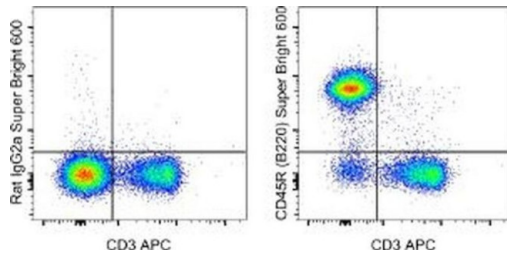

View more figures on thermofisher.com

Immunohistochemistry (30)

<p>Frontiers in immunology</p> <p>Role of CXCL5 in Regulating Chemotaxis of Innate and Adaptive Leukocytes in Infected Lungs Upon Pulmonary Influenza Infection.</p> <p>"Published figure using CD45R (B220) monoclonal antibody (Product # 63-0452-82) in Immunohistochemistry"</p> <p>Authors: Guo L,Li N,Yang Z,Li H,Zheng H,Yang J,Chen Y,Zhao X,Mei J,Shi H,Worthen GS,Liu L</p>	<p>Year</p> <p>2022</p>
<p>Cancers</p> <p>Targeting Ovarian Carcinoma with TSP-1:CD47 Antagonist TAX2 Activates Anti-Tumor Immunity.</p> <p>"Published figure using CD45R (B220) monoclonal antibody (Product # 63-0452-82) in Immunohistochemistry"</p> <p>Authors: Jeanne A,Sarazin T,Charlé M,Moali C,Fichel C,Boulagnon-Rombi C,Callewaert M,Andry MC,Diesis E,Delolme F,Rioul D,Dedieu S</p>	<p>Year</p> <p>2021</p>

[View more IHC references on thermofisher.com](#)

Immunohistochemistry (Paraffin) (2)

<p>Pathogens (Basel, Switzerland)</p> <p>Pituitary Adenylate Cyclase-Activating Polypeptide Alleviates Intestinal, Extra-Intestinal and Systemic Inflammatory Responses during Acute <i>Campylobacter jejuni</i>-induced Enterocolitis in Mice.</p> <p>"Published figure using CD45R (B220) monoclonal antibody (Product # 63-0452-82) in Immunohistochemistry (Paraffin)"</p> <p>Authors: Heimesaat MM,Mousavi S,Kløve S,Genger C,Weschka D,Tamas A,Reglodi D,Bereswill S</p>	<p>Year</p> <p>2020</p>
<p>European journal of microbiology & immunology</p> <p>Anti-inflammatory Effects of the Octapeptide NAP in Human Microbiota-Associated Mice Suffering from Subacute Ileitis.</p> <p>"Published figure using CD45R (B220) monoclonal antibody (Product # 63-0452-82) in Immunohistochemistry (Paraffin)"</p> <p>Authors: Escher U,Giladi E,Dunay IR,Bereswill S,Gozes I,Heimesaat MM</p>	<p>Year</p> <p>2018</p>

More applications with references on thermofisher.com

- IHC (F) (1)
- ICC/IF (5)
- Flow (62)
- Misc (1)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON-INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (i) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (ii) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (iii) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (iv) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.