

CD45.2 Monoclonal Antibody (104), Super Bright™ 600, eBioscience™

| Product Details | |
|-----------------------------|--|
| Size | 100 µg |
| Species Reactivity | Mouse |
| Host/Isotype | Mouse / IgG2a, kappa |
| Recommended Isotype Control | Mouse IgG2a kappa Isotype Control (eBM2a), Super Bright™ 600, eBioscience™ |
| Class | Monoclonal |
| Type | Antibody |
| Clone | 104 |
| Conjugate | Super Bright™ 600 |
| Excitation/Emission Max | 414/601 nm |
| Form | Liquid |
| Concentration | 0.2 mg/mL |
| Purification | Affinity chromatography |
| Storage buffer | PBS, pH 7.2, with BSA |
| Contains | 0.09% sodium azide |
| Storage conditions | 4°C, store in dark, DO NOT FREEZE! |
| RRID | AB_2662428 |

| Applications | Tested Dilution | Publications |
|----------------------------|-----------------|-----------------|
| Immunohistochemistry (IHC) | - | 1 Publication |
| Flow Cytometry (Flow) | 1 µg/test | 14 Publications |

Product Specific Information

Description: The 104 monoclonal antibody reacts with the mouse CD45 molecule, the leukocyte common antigen (LCA) in CD45.2-expressing mouse strains. The strains that express CD45.2 include the most commonly used mouse strains C57BL/6, BALB/c, C58, DBA/1, DBA/2, C3H/He, CBA, 129, A and AKR. CD45.2 is expressed by all leukocytes in these strains.

Applications Reported: This 104 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This 104 antibody has been tested by flow cytometric analysis of mouse splenocytes. This can be used at less than or equal to 1 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Super Bright 600 is a tandem dye that can be excited with the violet laser line (405 nm) and emits at 600 nm. We recommend using a 610/20 bandpass filter. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer for more information.


Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Protect this vial and stained samples from light.

Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222) (100 μ L of cell sample + 100 μ L of IC Fixation Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333) for up to 3 days in the dark at 2-8°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorochrome performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 405 nm; Emission: 600 nm; Laser: Violet Laser

Super Bright Polymer Dyes are sold under license from Becton, Dickinson and Company.

Product Images For CD45.2 Monoclonal Antibody (104), Super Bright™ 600, eBioscience™

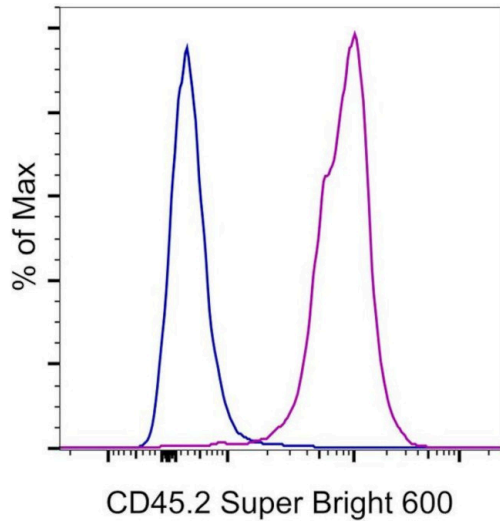

CD45.2 Antibody (63-0454-82) in Flow

Staining of SJL (blue histogram) and BALB/c (purple histogram) splenocytes with 0.5 μ g of Anti-Mouse CD45.2 Super Bright 600. Cells in the lymphocyte gate were used for analysis.

[View more figures on thermofisher.com](https://www.thermofisher.com)

Immunohistochemistry (1)

Molecular medicine reports

Aggravation of hepatic ischemiareperfusion injury with increased inflammatory cell infiltration is associated with the TGF/Smad3 signaling pathway.

"Published figure using CD45.2 monoclonal antibody (Product # 63-0454-82) in Immunohistochemistry"

Authors: Li H, Shen X, Tong Y, Ji T, Feng Y, Tang Y, Mai R, Ye J, Que T, Luo X

Year
2021

Flow Cytometry (14)

Frontiers in immunology

Ddb1 Is Essential for the Expansion of CD4⁺ Helper T Cells by Regulating Cell Cycle Progression and Cell Death.

"Published figure using CD45.2 monoclonal antibody (Product # 63-0454-82) in Flow Cytometry"

Authors: Yang L, Chen W, Li L, Xiao Y, Fan S, Zhang Q, Xia T, Li M, Hong Y, Zhao T, Li Q, Liu WH, Xiao N

Year
2021

Cell reports

Differentiation of fetal hematopoietic stem cells requires ARID4B to restrict autocrine KITLG/KIT-Src signaling.

"Published figure using CD45.2 monoclonal antibody (Product # 63-0454-82) in Flow Cytometry"

Authors: Young IC, Wu B, Andricovich J, Chuang ST, Li R, Tzatsos A, Wu RC, Wu MY

Year
2021

[View more Flow references on thermofisher.com](#)

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.