



CD223 (LAG-3) Monoclonal Antibody (3DS223H), Super Bright[™] 600, eBioscience[™]

Product Details	
Size	100 Tests
Species Reactivity	Human
Host/Isotype	Mouse / IgG1, kappa
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), Super Bright™ 600, eBioscience™
Class	Monoclonal
Туре	Antibody
Clone	3DS223H
Conjugate	Super Bright [™] 600
Excitation/Emission Max	414/601 nm
Immunogen	hCD223-FC PFA 4%
Form	Liquid
Concentration	5 μL/Test
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2688188

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	5 μL (0.06 μg)/test	5 Publications

Product Specific Information

Description: This 3DS223H monoclonal antibody recognizes human CD223 also known as Lymphocyte Activation Gene 3 (LAG-3). LAG-3 is a 70-kDa surface glycoprotein belonging to the Ig superfamily with homology to CD4. LAG-3 binds to MHC class II with higher affinity than CD4 and is thought to be involved in the negative regulation of T cell activation and homeostatic proliferation. Surface expression of LAG-3 has been reported on activated T cells (including regulatory T cells) and NK cells. CD8+ T cells usually express LAG-3 at significantly higher levels than CD4+ T cells. Coexpression of LAG-3 and CD49b has been proposed to identify human and mouse Type 1 regulatory T cells (Tr1 cells).

This 3DS223H antibody will recognize a formaldehyde-fixed epitope.

Applications Reported: This 3DS223H antibody has been reported for use in flow cytometric analysis.

Applications Tested: This 3DS223H antibody has been pre-diluted and tested by flow cytometric analysis of normal human peripheral blood cells. This may be used at 5 μ L (0.06 μ g) per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.

Super Bright 600 is a tandem dye that can be excited with the violet laser line (405 nm) and emits at 600 nm. We recommend using a 610/20 bandpass filter. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer for more information.


Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.


Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222) (100 µL of cell sample + 100 µL of IC Fixation Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 405 nm; Emission: 600 nm; Laser: Violet Laser

Super Bright Polymer Dyes are sold under license from Becton, Dickinson and Company.

Product Images For CD223 (LAG-3) Monoclonal Antibody (3DS223H), Super Bright™ 600, eBioscience™

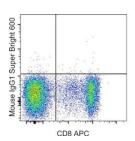

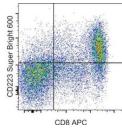

CD223 (LAG-3) Antibody (63-2239-42) in Flow

Normal human peripheral blood cells were stimulated for 3 days with Human IL-2 Recombinant Protein (Product # 14-8029-81), Anti-Human CD3 (Product # 16-0037-81), and Anti-Human CD28 Functional Grade Purifieds (Product # 16-0289-81). These cells were then surface stained with Anti-Human CD8a APC (Product # 17-0088-42) and Mouse IgG1 K Isotype Control Super Bright 600 (Product # 63-4714-82) (left) or Anti-Human CD223 (LAG-3) Super Bright 600 (right). Cells in the lymphocyte gate were used for analysis.

View more figures on thermofisher.com

□ 5 References

Flow Cytometry (5)

Frontiers in oncology

Low-density lipoprotein balances T cell metabolism and enhances response to anti-PD-1 blockade in a HCT116 spheroid model.

"Published figure using CD223 (LAG-3) monoclonal antibody (Product # 63-2239-42) in Flow Cytometry"

Authors: Babl N,Hofbauer J,Matos C,Voll F,Menevse AN,Rechenmacher M,Mair R,Beckhove P,Herr W,Siska PJ, Renner K,Kreutz M,Schnell A

Frontiers in immunology

Bispecific T-Cell Engaging Antibodies Against MUC16 Demonstrate Efficacy Against Ovarian Cancer in Monotherapy and in Combination With PD-1 and VEGF Inhibition.

"Published figure using CD223 (LAG-3) monoclonal antibody (Product # 63-2239-42) in Flow Cytometry" Authors: Yeku OO,Rao TD,Laster I,Kononenko A,Purdon TJ,Wang P,Cui Z,Liu H,Brentjens RJ,Spriggs D

Year 2023

Year 2021

View more Flow references on thermofisher.com

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization, Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package intentition?). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid not herein in valid not herein in valid not herein is valid not herein in valid not herein in valid not herein is valid not herein in valid