

TER-119 Monoclonal Antibody (TER-119), Super Bright™ 600, eBioscience™

Product Details	
Size	100 µg
Species Reactivity	Mouse
Host/Isotype	Rat / IgG2b, kappa
Recommended Isotype Control	Rat IgG2b kappa Isotype Control (eB149/10H5), Super Bright™ 600, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	TER-119
Conjugate	Super Bright™ 600
Excitation/Emission Max	414/601 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2688117

Applications	Tested Dilution	Publications
Western Blot (WB)	-	1 Publication
Flow Cytometry (Flow)	0.5 µg/test	15 Publications

Product Specific Information

Description: The TER-119 monoclonal antibody reacts with mouse erythroid cells from early proerythroblast to mature erythrocyte stages. The TER-119 antigen is present in yolk sac, fetal and newborn liver, but is not expressed by cells carrying BFU-E and CFU-E activities. Several erythroleukemia cell lines tested so far are negative for expression of TER-119 antigen even after dimethylsulfoxide stimulation. Biochemical and molecular analysis of the TER-119 antigen indicate that this molecule is associated with the surface glycoprotein A, but is not a typical glycoprotein.^{^M}

^{^M}

Applications Reported: This TER-119 antibody has been reported for use in flow cytometric analysis.^{^M}

^{^M}

Applications Tested: This TER-119 antibody has been tested by flow cytometric analysis of mouse bone marrow cells. This may be used at less than or equal to 0.25 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.^{^M}

^{^M}

Super Bright 600 is a tandem dye that can be excited with the violet laser line (405 nm) and emits at 600 nm. We recommend using a 610/20 bandpass filter. Please make sure that your instrument is capable of detecting this fluorochrome.^{^M}

^{^M}

When using two or more Super Bright dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer for more information.^{^M}

^M

Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.^M

^M

Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222) (100 μ L of cell sample + 100 μ L of IC Fixation Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.^M


^M

Excitation: 405 nm; Emission: 600 nm; Laser: Violet Laser^M

^M

Super Bright Polymer Dyes are sold under license from Becton, Dickinson and Company.

Product Images For TER-119 Monoclonal Antibody (TER-119), Super Bright™ 600, eBioscience™

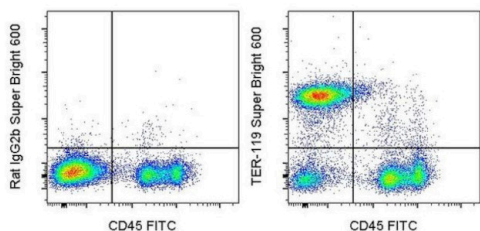

TER-119 Antibody (63-5921-82) in Flow

C57BL/6 mouse bone marrow cells were stained with CD45 Monoclonal Antibody, FITC (Product # 11-0451-82) and 0.25 μ g of Rat IgG2b kappa Isotype Control, Super Bright 600 (Product # 63-4031-82) (left) or 0.25 μ g of TER-119 Monoclonal Antibody, Super Bright 600 (right). Total viable cells were used for analysis.

[View more figures on thermofisher.com](https://www.thermofisher.com)

16 References

Western Blot (1)

Life science alliance

SnRNA sequencing defines signaling by RBC-derived extracellular vesicles in the murine heart.

"Published figure using TER-119 monoclonal antibody (Product # 63-5921-82) in Western Blot"

Authors: Valkov N, Das A, Tucker NR, Li G, Salvador AM, Chaffin MD, Pereira De Oliveira Junior G, Kur I, Gokulnath P, Ziegler O, Yeri A, Lu S, Khamesra A, Xiao C, Rodosthenous R, Srinivasan S, Toxavidis V, Tigges J, Laurent LC, Momma S, Kitchen R, Ellinor P, Ghiran I, Das S

Year
2021

Flow Cytometry (15)

PloS one

Engineered red blood cells carrying PCSK9 inhibitors persistently lower LDL and prevent obesity.

"Published figure using TER-119 monoclonal antibody (Product # 63-5921-82) in Flow Cytometry"

Authors: Deshycka R, Sudaryo V, Huang NJ, Xie Y, Smeding LY, Choi MK, Ploegh HL, Lodish HF, Pishesha N

Year
2021

Frontiers in immunology

Sex Matters: Physiological Abundance of Immuno-Regulatory CD71+ Erythroid Cells Impair Immunity in Females.

"Published figure using TER-119 monoclonal antibody (Product # 63-5921-82) in Flow Cytometry"

Authors: Mashhour S, Koleva P, Huynh M, Okoye I, Shahbaz S, Elahi S

Year
2021

[View more Flow references on thermofisher.com](#)

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.