CD326 (EpCAM) Monoclonal Antibody (1B7), Super Bright[™] 600, eBioscience[™]

Product Details

Size	25 Tests	
Species Reactivity	Human	
Host/Isotype	Mouse / IgG1, kappa	
Recommended Isotype Control	Mouse IgG1 kappa Isotype Control (P3.6.2.8.1), Super Bright [™] 600, eBioscience [™]	
Class	Monoclonal	
Туре	Antibody	
Clone	1B7	
Conjugate	Super Bright™ 600	
Form	Liquid	
Concentration	5 µL/Test	
Storage conditions	4° C, store in dark, DO NOT FREEZE!	
RRID	AB_2802444	

Applications	Tested Dilution	Publications
Immunocytochemistry (ICC/IF)	-	0 Publication
Flow Cytometry (Flow)	5 μL (0.5 μg)/Test	1 Publication

Product Specific Information

Description: EpCAM (Epithelial cell adhesion molecule, CD326, KSA, TROP1) is a 40 kD cell-surface adhesion molecule participating in homophilic, calcium-independent cell-cell interactions. EpCAM is a type-I transmembrane protein, and is expressed primarily on the basolateral surface of most epithelia. Although normal epithelia express low levels of EpCAM, increased expression has been correlated with increased proliferation and progression to a mesenchymal phenotype. EpCAM has also been used as a diagnostic marker on circulating metastatic carcinoma cells, while cancer cells of non-epithelial origin do not express EpCAM.

Applications Reported: This 1B7 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This 1B7 antibody has been pre-diluted and tested by flow cytometric analysis of A549 cells. This may be used at 5 μ L (0.5 μ g) per test. A test is defined as the amount (μ g) of antibody that will stain a cell sample in a final volume of 100 μ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test.

Super Bright 600 is a tandem dye that can be excited with the violet laser line (405 nm) and emits at 600 nm. We recommend using a 610/20 bandpass filter. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer for more information.

Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.


Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222-49) (100 µL of cell sample + 100 µL of IC Fixation Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333-57) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be

made, but clone specific performance should be determined empirically.

Excitation: 405 nm; Emission: 600 nm; Laser: Violet Laser

Super Bright Polymer Dyes are sold under license from Becton, Dickinson and Company.

Product Images For CD326 (EpCAM) Monoclonal Antibody (1B7), Super Bright[™] 600, eBioscience[™]

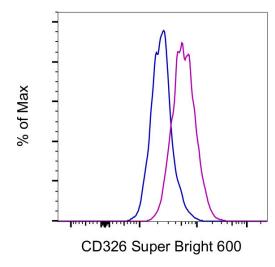

CD326 (EpCAM) Antibody (63-9326-41) in Flow

A549 cells were stained with Mouse IgG1 kappa Isotype Control, Super Bright 600 (Product # 63-4714-82) (blue histogram) or CD326 Monoclonal Antibody, Super Bright 600 (purple histogram). Total viable cells were used for analysis, as determined by 7-AAD (Product # 00-6993-50).

□ 1 Reference

Flow Cytometry (1)

Cell death & disease FZD5 contributes to TNBC proliferation, DNA damage repair and stemness. **Year** 2020

"Published figure using CD326 (EpCAM) monoclonal antibody (Product # 63-9326-42) in Flow Cytometry" Authors: Sun Y,Wang Z,Na L,Dong D,Wang W,Zhao C

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. No OTHER WARRANTIES, EXPRESS OR IMPLED, ARRANTES OF MERCHANTABILITY. FITNESS FOR ANY PARTICILAR PURPAGE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS QUINNG THE WARRANTY PERIOD IS UNITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCTS) ANY PARTICILAR PURPAGE. (III) MISUEF, APLICE OR REFUND FOR THE NON-CONFORMING PRODUCTS) AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (III) MISUEF, APLICE OR OR BY BUYER, (III) USE OF THE NON-CONFORMING PRODUCTS). IN MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HADDLING OF THE PRODUCTS SA THE RESULT OF (II) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (III) MISUEF, APLICE OR OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FWHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HADDLING OF THE PRODUCTS. IN A MANNER FWHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HADDLING OF THE PRODUCTS. IN MAINTES APPLICE OR THE HOULD'S IN A MANNER FWHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HADDLING OF THE PRODUCTS. IN A MANNER FWHICH THEY WER