

CD4 Monoclonal Antibody (GK1.5), Super Bright™ 702, eBioscience™

Product Details	
Size	25 µg
Species Reactivity	Mouse
Host/Isotype	Rat / IgG2b, kappa
Recommended Isotype Control	Rat IgG2b kappa Isotype Control (eB149/10H5), Super Bright™ 702, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	GK1.5
Conjugate	Super Bright™ 702
Excitation/Emission Max	413/702 nm
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2722962

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	5 Publications
Immunohistochemistry (Frozen) (IHC (F))	-	1 Publication
Immunocytochemistry (ICC/IF)	-	2 Publications
Flow Cytometry (Flow)	0.25 µg/test	78 Publications

Product Specific Information

Description: The GK1.5 monoclonal antibody reacts with the mouse CD4 molecule, a 55 kDa cell surface receptor expressed by a majority of thymocytes, subpopulation of mature T cells and dendritic cells. CD4 binds to MHC class II on the surface of antigen presenting cells and plays an important role both in T cell development and in optimal functioning of mature T cells. In T cells, CD4 associates with protein tyrosine kinase p56lck through its cytoplasmic tail. Binding of GK1.5 is blocked by RM4-5.

Applications Reported: This GK1.5 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This GK1.5 antibody has been tested by flow cytometric analysis of mouse splenocytes. This can be used at less than or equal to 0.25 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Super Bright 702 is a tandem dye that can be excited with the violet laser line (405 nm) and emits at 702 nm. We recommend using a 710/50 bandpass filter. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer for more information.™

™

Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.™

™

Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222) (100 µL of cell sample + 100 µL of IC Fixation Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.™


™

Excitation: 405 nm; Emission: 702 nm; Laser: Violet Laser™

™


Super Bright Polymer Dyes are sold under license from Becton, Dickinson and Company.

Product Images For CD4 Monoclonal Antibody (GK1.5), Super Bright™ 702, eBioscience™

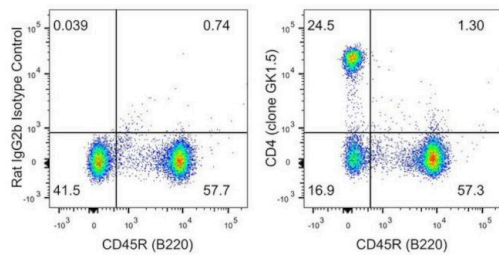

CD4 Antibody (67-0041-80)

Staining of mouse splenocytes. As expected based on known relative expression patterns, CD4 clone GK1.5 stains about 25% of total lymphocytes and does not stain CD45R (B220)+ B cells. Details: Balb/c splenocytes were surface stained with CD45R (B220, clone RA3-6B2) and Rat IgG2b Isotype Control or CD4 (clone GK1.5). Cells in the lymphocyte gate were used for analysis. {RE}

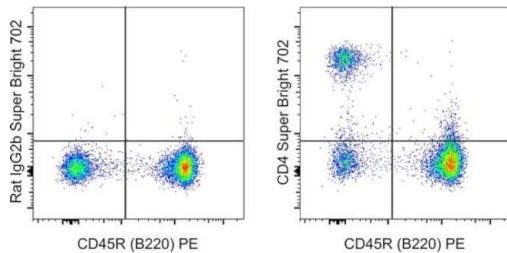

CD4 Antibody (67-0041-80) in Flow

C57BL/6 mouse splenocytes were stained with CD45R (B220) Monoclonal Antibody, PE (Product # 12-0452-82) and 0.125 µg of Rat IgG2b kappa Isotype Control, Super Bright 702 (Product # 67-4031-82) (left) or 0.125 µg of CD4 Monoclonal Antibody, Super Bright 702 (right). Total viable cells were used for analysis.

[View more figures on thermofisher.com](https://www.thermofisher.com)

86 References

Immunohistochemistry (5)

Cell reports

Btla signaling in conventional and regulatory lymphocytes coordinately tempers humoral immunity in the intestinal mucosa.

"Published figure using CD4 monoclonal antibody (Product # 67-0041-82) in Immunohistochemistry"

Authors: Stienne C, Virgen-Slane R, Elmén L, Vény M, Huang S, Nguyen J, Chappell E, Balmert MO, Shui JW, Hurchla MA, Kronenberg M, Peterson SN, Murphy KM, Ware CF, Šedý JR

Year
2022

Cancers

Targeting Ovarian Carcinoma with TSP-1:CD47 Antagonist TAX2 Activates Anti-Tumor Immunity.

"Published figure using CD4 monoclonal antibody (Product # 67-0041-82) in Flow Cytometry"

Authors: Jeanne A, Sarazin T, Charlé M, Moali C, Fichel C, Boulagnon-Rombi C, Callewaert M, Andry MC, Diesis E, Delolme F, Rioult D, Dedieu S

Year
2021

[View more IHC references on thermofisher.com](#)

Immunohistochemistry (Frozen) (1)

International journal of molecular sciences

Positive Effects of Oral Antibiotic Administration in Murine Chronic Graft-Versus-Host Disease.

"Published figure using CD4 monoclonal antibody (Product # 67-0041-82) in Immunocytochemistry"

Authors: Sato S, Shimizu E, He J, Ogawa M, Asai K, Yazu H, Rusch R, Yamane M, Yang F, Fukuda S, Kawakami Y, Tsubota K, Ogawa Y

Year
2021

Immunocytochemistry (2)

Nature nanotechnology

Cell-mimicking nanodecoys neutralize SARS-CoV-2 and mitigate lung injury in a non-human primate model of COVID-19.

"Published figure using CD4 monoclonal antibody (Product # 67-0041-82) in Immunocytochemistry"

Authors: Li Z, Wang Z, Dinh PC, Zhu D, Popowski KD, Lutz H, Hu S, Lewis MG, Cook A, Andersen H, Greenhouse J, Pessaint L, Lobo LJ, Cheng K

Year
2021

More applications with references on thermofisher.com

Flow (78)

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.