



# CD19 Monoclonal Antibody (eBio1D3 (1D3)), Super Bright™ 702, eBioscience™

Size25 μgSpecies ReactivityMouseHost/IsotypeRat / IgG2a, kappaRecommended Isotype ControlRat IgG2a kappa Isotype ControClassMonoclonalTypeAntibodyCloneeBio1D3 (1D3)ConjugateSuper Bright™ 702Excitation/Emission Max413/702 nmFormLiquid	ol (eBR2a), Super Bright™ 702, eBioscience™
Host/Isotype Rat / IgG2a, kappa  Recommended Isotype Control  Class Monoclonal  Type Antibody  Clone eBio1D3 (1D3)  Conjugate Super Bright™ 702  Excitation/Emission Max	ol (eBR2a), Super Bright™ 702, eBioscience™
Recommended Isotype Control  Class  Monoclonal  Type  Antibody  Clone  eBio1D3 (1D3)  Conjugate  Excitation/Emission Max  Rat IgG2a kappa Isotype Contro  Monoclonal  Antibody  EBio1D3 (1D3)  Super Bright™ 702  413/702 nm	ol (eBR2a), Super Bright™ 702, eBioscience™
Control  Class Monoclonal  Type Antibody  Clone eBio1D3 (1D3)  Conjugate Super Bright™ 702  Excitation/Emission Max  Excitation/Emission	ol (eBR2a), Super Bright™ 702, eBioscience™
Type Antibody  Clone eBio1D3 (1D3)  Conjugate Super Bright™ 702  Excitation/Emission Max  413/702 nm	
Clone eBio1D3 (1D3)  Conjugate Super Bright™ 702  Excitation/Emission Max  413/702 nm	
ConjugateSuper Bright™ 702Excitation/Emission Max413/702 nm	
Excitation/Emission Max 413/702 nm	
Max 413/702 nm	
Form Liquid	
Elquid	
Concentration 0.2 mg/mL	
Purification Affinity chromatography	
Storage buffer PBS, pH 7.2, with BSA	
Contains 0.09% sodium azide	
Storage conditions 4° C, store in dark, DO NOT FRE	EEZE!
RRID AB_2717134	

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	1 Publication
Flow Cytometry (Flow)	0.5 μg/test	24 Publications

#### **Product Specific Information**

Description: The eBio1D3 monoclonal antibody reacts with mouse CD19, a 95 kDa transmembrane glycoprotein. CD19 is expressed by B cells during all stages of development excluding the terminally differentiated plasma cells. Follicular dendritic cells also express CD19. Together CD21, CD81, MHC class II, and CD19 form a multimolecular complex that associates with the BCR. Signaling through CD19 induces tyrosine phosphorylation, calcium flux and proliferation of B cells.\t

Applications Reported: This eBio1D3 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This eBio1D3 antibody has been tested by flow cytometric analysis of mouse splenocytes. This may be used at less than or equal to 0.5  $\mu$ g per test. A test is defined as the amount ( $\mu$ g) of antibody that will stain a cell sample in a final volume of 100  $\mu$ L. Cell number should be determined empirically but can range from 10^5 to 10^8 cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Super Bright 702 is a tandem dye that can be excited with the violet laser line (405 nm) and emits at 702 nm. We recommend using a 710/50 bandpass filter. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer for more information.


Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.

Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222) (100  $\mu$ L of cell sample + 100  $\mu$ L of IC Fixation Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 405 nm; Emission: 702 nm; Laser: Violet Laser

Super Bright Polymer Dyes are sold under license from Becton, Dickinson and Company.

## Product Images For CD19 Monoclonal Antibody (eBio1D3 (1D3)), Super Bright™ 702, eBioscience™

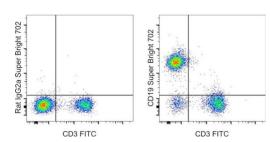

### CD19 Antibody (67-0193-80) in Flow

BALB/c mouse splenocytes were stained with CD3e Monoclonal Antibody, FITC (Product # 11-0031-82) and of Rat IgG2a kappa Isotype Control, Super Bright 702 (Product # 67-4321-82) (left) or of CD19 Monoclonal Antibody, Super Bright 702 (right). Cells in the lymphocyte gate were used for analysis.

View more figures on thermofisher.com

#### □ 25 References

# Immunohistochemistry (1)

#### MedComm

Single-cell transcriptomics reveals distinct cell response between acute and chronic pulmonary infection of *Pseudomonas aeruginosa*.

**Year** 2022

"Published figure using CD19 monoclonal antibody (Product # 67-0193-82) in Immunohistochemistry"

Authors: Hu X,Wu M,Ma T,Zhang Y,Zou C,Wang R,Zhang Y,Ren Y,Li Q,Liu H,Li H,Wang T,Sun X,Yang Y,Tang M,Li X,Li J,Gao X,Li T,Zhou X

#### Flow Cytometry (24)

Cells

Cissus quadrangularis (Hadjod) Inhibits RANKL-Induced
Osteoclastogenesis and Augments Bone Health in an Estrogen-Deficient
Preclinical Model of Osteoporosis Via Modulating the Host
Osteoimmune System.

**Year** 2023

"Published figure using CD19 monoclonal antibody (Product # 67-0193-82) in Flow Cytometry" Authors: Azam Z,Sapra L,Baghel K,Sinha N,Gupta RK,Soni V,Saini C,Mishra PK,Srivastava RK

Investigative ophthalmology & visual science

# Role of FGF10/FGFR2b Signaling in Homeostasis and Regeneration of Adult Lacrimal Gland and Corneal Epithelium Proliferation.

"Published figure using CD19 monoclonal antibody (Product # 67-0193-82) in Flow Cytometry"

Authors: Finburgh EN, Mauduit O, Noguchi T, Bu JJ, Abbas AA, Hakim DF, Bellusci S, Meech R, Makarenkova HP, Afshari

**Year** 2023

View more Flow references on thermofisher.com

#### More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production obcumentation, specifications and/or accompanying package inserts ("Documentation"), No claim of suitability for use in applications regulated by EDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, the warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, the warranty provided herein is valid only when used by properly trained individuals. Unless otherwise sealing and in the product is used in the Product is used in the provided and provided in the provided provided in the provided is used to sample in the provided in the provided provided in the provided provided in the provided pro