

MERTK Monoclonal Antibody (DS5MMER), Super Bright™ 702, eBioscience™

Product Details	
Size	100 µg
Species Reactivity	Mouse
Host/Isotype	Rat / IgG2a, kappa
Recommended Isotype Control	Rat IgG2a kappa Isotype Control (eBR2a), Super Bright™ 702, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	DS5MMER
Conjugate	Super Bright™ 702
Excitation/Emission Max	413/702 nm
Immunogen	FC-tagged full protein expressed in insect cells
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2, with BSA
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2717175

Applications	Tested Dilution	Publications
Immunohistochemistry (IHC)	-	1 Publication
Flow Cytometry (Flow)	1 µg/test	2 Publications

Product Specific Information

Description: This DS5MMER monoclonal antibody recognizes mouse MerTK, a 170-210 kDa member of the TAM family of tyrosine kinase receptors that also includes Axl and Tyro3. MerTK is expressed on tissue macrophages and is involved in the removal of apoptotic cells. This process relies on two soluble ligands of MerTK, Protein S and Gas6 that bind to phosphatidylserine found on the outer leaflet of the plasma membrane of cells undergoing apoptosis. Upon binding these ligands, MerTK undergoes autophosphorylation at multiple tyrosine residues that activate the PI3K and Akt pathways. This results in the phagocytosis of apoptotic cells and also results in the direct inhibition of TLR-induced production of pro-inflammatory cytokines. In addition, MerTK may function as a putative entry receptor for filoviruses. Deficiency of MerTK causes general autoimmunity, inflammation and accumulation of apoptotic bodies. MerTK is constitutively released from the cell surface by metalloproteinases and thus is present in the serum and culture medium. This process can be enhanced by stimulation with LPS. MerTK can be used to help discriminate macrophages from dendritic cells. MerTK is often expressed on malignant cells and may be implicated in immune evasion.

The DS5MMER antibody will work on paraformaldehyde-fixed cells.

Applications Reported: This DS5MMER antibody has been reported for use in flow cytometric analysis.

Applications Tested: This DS5MMER antibody has been tested by flow cytometric analysis of resident peritoneal

macrophages. This may be used at less than or equal to 1.0 µg per test. A test is defined as the amount (µg) of antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.

Super Bright 702 is a tandem dye that can be excited with the violet laser line (405 nm) and emits at 702 nm. We recommend using a 710/50 bandpass filter. Please make sure that your instrument is capable of detecting this fluorochrome.

When using two or more Super Bright dye-conjugated antibodies in a staining panel, it is recommended to use Super Bright Complete Staining Buffer (Product # SB-4401) to minimize any non-specific polymer interactions. Please refer to the datasheet for Super Bright Staining Buffer for more information.


Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.

Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222) (100 µL of cell sample + 100 µL of IC Fixation Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 405 nm; Emission: 702 nm; Laser: Violet Laser

Super Bright Polymer Dyes are sold under license from Becton, Dickinson and Company.

Product Images For MERTK Monoclonal Antibody (DS5MMER), Super Bright™ 702, eBioscience™

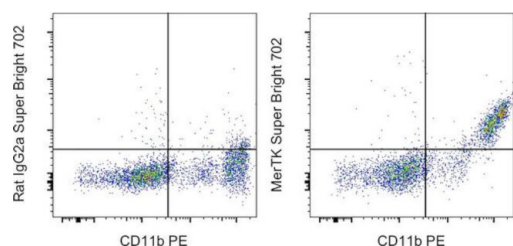

MERTK Antibody (67-5751-82) in Flow

Staining of mouse resident peritoneal exudate cells with CD11b Monoclonal Antibody, PE (Product # 12-0112-82) and 1.0 µg of Rat IgG2a K Isotype Control, Super Bright 702 (Product # 67-4321-82) (left) or 1.0 µg of MerTK Monoclonal Antibody, Super Bright 702 (right). Total viable cells were used for analysis.

[View more figures on thermofisher.com](https://www.thermofisher.com)

3 References

Immunohistochemistry (1)

MedComm

Single-cell transcriptomics reveals distinct cell response between acute and chronic pulmonary infection of *Pseudomonas aeruginosa*.

"Published figure using MERTK monoclonal antibody (Product # 67-5751-82) in Immunohistochemistry"

Authors: Hu X,Wu M,Ma T,Zhang Y,Zou C,Wang R,Zhang Y,Ren Y,Li Q,Liu H,Li H,Wang T,Sun X,Yang Y,Tang M,Li X, Li J,Gao X,Li T,Zhou X

Year
2022

Flow Cytometry (2)

Cellular and molecular gastroenterology and hepatology

AXL Expression on Homeostatic Resident Liver Macrophages Is Reduced in Cirrhosis Following GAS6 Production by Hepatic Stellate Cells.

"Published figure using MERTK monoclonal antibody (Product # 67-5751-82) in Flow Cytometry"

Authors: Pop OT,Geng A,Flint E,Singanayagam A,Ercan C,Possamai L,Patel VC,Kuenzler P,Meier MA,Soysal S,Hruz P,Kollmar O,Tatham KC,Ward JK,Müllhaupt B,Weber A,Wendon J,Niess JH,Heim M,Semela D,Weston C,Antoniades CG,Terracciano LM,Triantafyllou E,Brenig RG,Bernsmeier C

Year
2023

Cell death & disease

Ligand-dependent kinase activity of MERTK drives efferocytosis in human iPSC-derived macrophages.

"Published figure using MERTK monoclonal antibody (Product # 67-5751-82) in Flow Cytometry"

Authors: Wanke F,Gutbier S,Rümmelin A,Steinberg M,Hughes LD,Koenen M,Komuczki J,Regan-Komito D,Wagage S,Hesselmann J,Thoma R,Brugger D,Christopeit T,Wang H,Point F,Hallet R,Ghosh S,Rothlin CV,Patsch C,Geering B

Year
2021

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON-INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.